


# MEDRAD

MCT and MCT *Plus*  
Injection Systems

Service Manual

KMP 860P

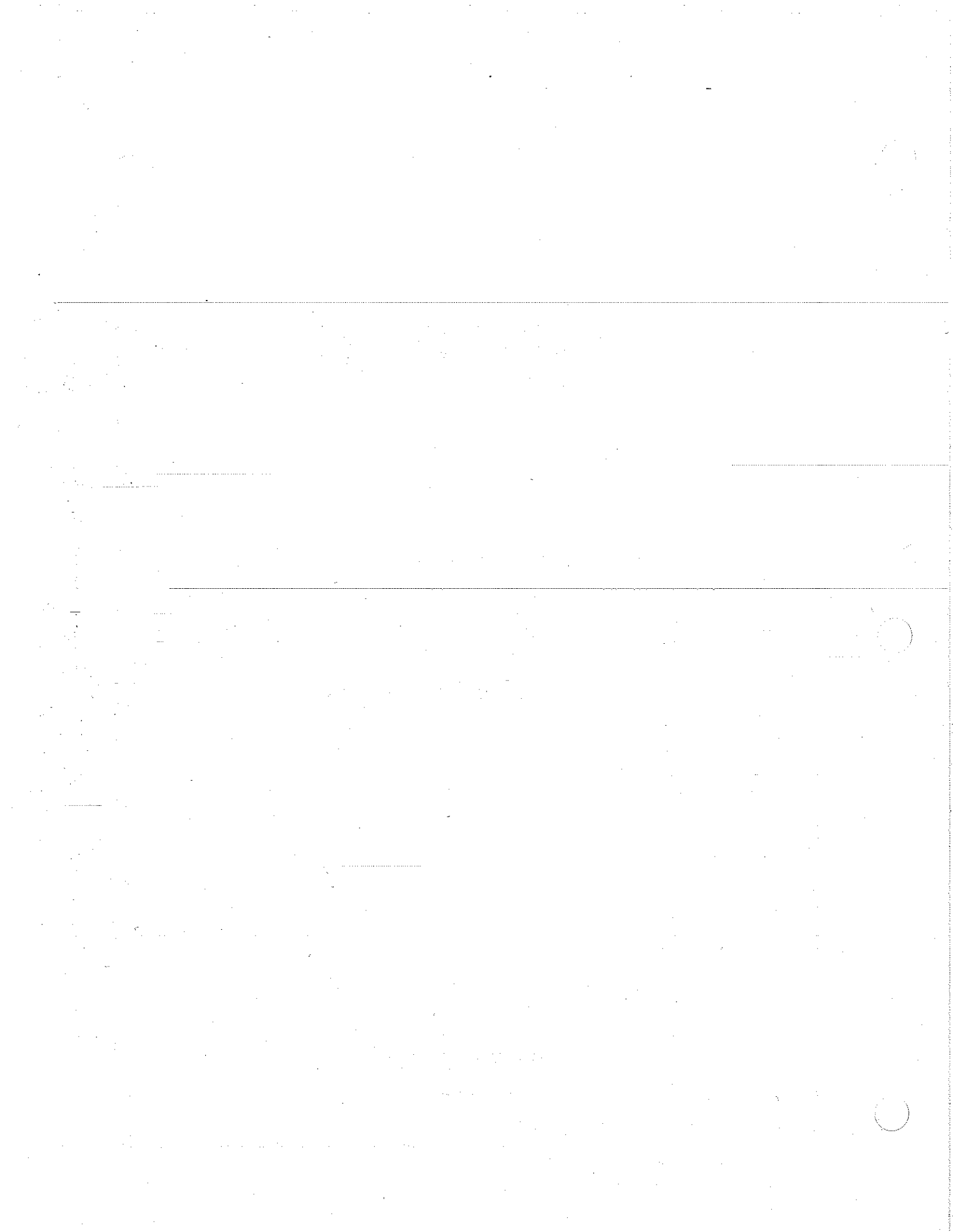Revision C



# MCT and MCT *Plus* SERVICE MANUAL

<b>Table of Contents</b>	1-1
<b>Table of Illustrations (Line-Art)</b>	1-5
<b>Table of Illustrations (Drawings)</b>	1-7
<b>Preface</b>	
Copyright Notice	2-1
Disclaimer	2-1
Warnings/Cautions	2-2
<b>Condensed Operator Guide</b>	
Control Panel	3-1
OP Message Summary	3-3
<b>MCT / MCT Plus Maintenance Procedures</b>	
Recommended Preventive Maintenance	4-1
Recommended Preventive Maintenance Schedule	4-1
Inspection Procedures	4-2
Mounting Systems Inspection	4-3
Cleaning Procedures	4-6
Electrical Leakage Check	4-6
<b>System Specifications</b>	
Specifications	5-1
Dimensions	5-3
<b>MCT / MCT Plus Operational Checkout</b>	
Entering Injections	6-1
Multi-Level Programming	6-3
Storing an Injection Program	6-4
Recalling an Injection Program	6-5
Editing Injection Programs	6-5
The HOLD Feature	6-6
<b>CPU Card</b>	
CPU Card Type B - Summary of Operation	7-2
Block Diagram: CPU Card Type B	7-3
CPU Card Type C - Summary of Operation	7-4
Block Diagram: CPU Card Type C	7-5

<b>Input/Output Card</b>		
	I/O Card - Summary of Operation	8 - 1
	Block Diagram: I/O Card	8 - 2
<b>Servo Control Card</b>		
	Servo Control Card - Summary of Operation	9 - 1
	Block Diagrams: Servo Control Card (2)	9 - 4
<b>Pre-Programmed Injection Card</b>		
	PPI Card - Type B (number 25664 and below) Summary of Operation	10 - 2
	Block Diagram: PPI Card - Type B	10 - 3
	PPI Card - Type C (number 25665 and above) Summary of Operation	10 - 4
	Block Diagram: PPI Card - Type C	10 - 5
<b>Control Panel Circuit</b>		
	Control Panel - Summary of Operation	11 - 1
	Block Diagram: Control Panel	11 - 2
<b>Power Drive Circuit</b>		
	Description of Power Drive Circuit	12 - 1
	PDC Circuit Description	12 - 1
	Block Diagrams: PDC and Relay Cards (2)	12 - 5
	PDCP Circuit Description	12 - 7
	Block Diagrams: PDCP and Relay Cards (2)	12 - 12
<b>Mechanical Stop Drive</b>		
	Mechanical Stop Drive Card - Summary of Operation	13 - 1
	Block Diagram: Mechanical Stop Drive Card	13 - 2
<b>Power Supply Circuit</b>		
	Power Supply Card - Summary of Operation	14 - 1
	Block Diagram: Power Supply Card	14 - 2
<b>Injector Head</b>		
	Description of Injector Head	15 - 1
	Description of Motor and Drive System	15 - 1
	Description of Feedback Pots	15 - 2
	Description of Syringe Heater	15 - 2
	Circuit Descriptions:	
	MCT Injector Head	15 - 2
	MCT <i>Plus</i> Injector Head Card A, MPH	15 - 4
	Block Diagram: Head Card A, MPH	15 - 7
	MCT <i>Plus</i> Injector Head Card B, H5D	15 - 8
	Block Diagram: Head Card B, H5M	15 - 11

## **Imaging System Interface - ISI 100**

Description of Imaging System Interface ISI 100	16 - 1
Cautions	16 - 1
Analog Output Signals	16 - 3
Block Diagram: Analog Output Signals	16 - 4
Imaging System Interface Circuits	16 - 5
Block Diagram: Imaging System Interface	16 - 8

## **Remote Panel**

Description of Remote Panel	17 - 1
Remote Panel Controller Card	17 - 2
Digital Logic	17 - 2
Block Diagram: Remote Panel Digital Logic	17 - 4
DC to DC Converter	17 - 5
Block Diagram: Remote Panel DC to DC Converter	17 - 6

## **Head Mounting Systems**

Head Mounting Systems Overview	18 - 1
OCS 115	18 - 2
OCS 125	18 - 4
KMA 320-T	18 - 6

## **Troubleshooting Guide**

Tools and Test Equipment Required	19 - 1
Replacement Parts:	
Suggested Spare Parts	19 - 3
Critical Parts	19 - 4
How to Order Parts	19 - 7
Service Access Card (SAC)	19 - 8
Troubleshooting Guide	
Introduction	19 - 11
General Troubleshooting Guidelines	19 - 12
Suggested Troubleshooting Sequence	19 - 13
Message Faults:	19 - 14
List of Error Message Codes	19 - 14
Non-Message Faults:	19 - 24
Dead Unit	19 - 25
No Drive	19 - 26
Will Not Inject	19 - 27
Garbled Message	19 - 27

## **Control Unit Disassembly and Reassembly**

Control Unit	20 - 1
Control Unit Cards Removal and Replacement	20 - 1
Control Panel Overlay Removal and Replacement	20 - 2
Injector Head Cover Removal and Replacement	20 - 3
Syringe Heater Removal and Replacement	20 - 4
Plunger Position Pot Removal and Replacement	20 - 5
Mechanical Stop Pot Removal and Replacement	20 - 8
Mechanical Stop Belt Removal and Replacement	20 - 10
Plunger Motor Drive Belt Removal and Replacement	20 - 10
Plunger Drive Motor Removal and Replacement	20 - 11
Mechanical Stop Motor Removal and Replacement	20 - 13
Head Control Panel Removal and Replacement	20 - 19
Turret Limit Switch Removal and Replacement	20 - 21
Mechanical Stop Limit Switch Removal and Replacement	20 - 22
Injector Head Cable Removal and Replacement	20 - 23
Power Switch Removal and Replacement	20 - 24
Oil Seal Removal and Replacement	20 - 25
How to Tighten the Pivot Knuckle	20 - 27
Adapting to Different Line Voltages	20 - 28
Adapting to Different Line Frequencies	20 - 30
Line Frequency Checkout	20 - 31

<b>Service Parts List</b>	21 - 1
---------------------------	--------

## **Interchangeability Information**

Mark V 1000 Series	22 - 1
Mark V 20000 Series	22 - 1
Mark V <i>Plus</i>	22 - 1
MCT	22 - 2
MCT <i>Plus</i>	22 - 2
Interchangeability Notes	22 - 2
Acronym Definitions	22 - 3
MCT / MCT <i>Plus</i> / Mark V / Mark V <i>Plus</i> Hardware	
Interchangeability Matrix of PC Cards	22 - 8
Interchangeability Matrix of PC Cards (Multinational)	22 - 9

<b>Drawings</b>	23 - 1
-----------------	--------

## Table of Illustrations - *Line Art*

<i>Figure Number</i>	<i>Figure Name</i>	<i>Page</i>	
Figure	4.1	Head Support Arm Inspection Points	4 - 3
	4.2	OCS Inspection Points	4 - 4
	4.3	OCS Inspection Points (cont.)	4 - 5
Figure	5.1	Control Console Dimensions	5 - 3
	5.2	Remote Panel Dimensions	5 - 3
	5.3	Front Load Injector Head Dimensions	5 - 4
	5.4	Turret Style Injector Head Dimensions	5 - 4
	5.5	Injector Head Pedestal KMA 320-T	5 - 5
Figure	7.1	Block Diagram: CPU Card Type B	7 - 3
	7.2	Block Diagram: CPU Card Type C	7 - 5
Figure	8.1	Block Diagram: I/O Card	8 - 2
Figure	9.1	Block Diagram: Servo Control Card (1 of 2)	9 - 3
	9.2	Block Diagram: Servo Control Card (2 of 2)	9 - 4
Figure	10.1	Block Diagram: PPI Card Type B	10 - 3
	10.2	Block Diagram: PPI Card Type C	10 - 5
Figure	11.1	Block Diagram: Control Panel	11 - 2
Figure	12.1	Block Diagram: PDC and Relay Cards (1 of 2)	12 - 5
	12.2	Block Diagram: PDC and Relay Cards (2 of 2)	12 - 6
	12.3	Block Diagram: PDCP and Relay Cards (1 of 2)	12 - 12
	12.4	Block Diagram: PDCP and Relay Cards (2 of 2)	12 - 13
Figure	13.1	Block Diagram: Mechanical Stop Drive Card	13 - 2
Figure	14.1	Block Diagram: Power Supply Card	14 - 2
Figure	15.1	Block Diagram: Head Card A, MPH	15 - 7
	15.2	Block Diagram: Head Card B, HSD	15 - 11

<i>Figure Number</i>	<i>Figure Name</i>	<i>Page</i>	
Figure	16.1	Rear Connector Panel	16 - 2
	16.2	Block Diagram: Analog Output Signals	16 - 4
	16.3	Block Diagram: Imaging System Interface	16 - 8
Figure	17.1	Block Diagram: Remote Panel Digital Logic	17 - 4
	17.2	Block Diagram: DC to DC Converter	17 - 6
Figure	18.1	OCS 115	18 - 2
	18.2	OCS 125	18 - 4
	18.3	Injector Head Pedestal KMA 320-T	18 - 6
Figure	19.1	Test Point on the Service Access Card	19 - 8
	19.2	Suggested Troubleshooting Sequence	19 - 13
Figure	20.1	Syringe Heater Cable and J29 Connector	20 - 4
	20.2	Plunger Pot Terminal Identification	20 - 6
	20.3	Head Covers and Cable Brackets	20 - 7
	20.4	Mechanical Stop Pot Terminal Locations	20 - 9
	20.5	Plunger Motor Drive Belt Tension Adjustment	20 - 10
	20.6	Plunger Motor Terminal Identification	20 - 11
	20.7a	Mechanical Stop Motor Mount Configuration (Types 1 and 2)	20 - 13
	20.7b	Mechanical Stop Motor Mount Configuration (Type 3)	20 - 15
	20.7c	Mechanical Stop Motor Mount Configuration (Type 4)	20 - 17
	20.8	Turret Limit Switch Modification	20 - 21
	20.9	Power Switch Terminal Orientation	20 - 24
	20.10	Front Plate Types	20 - 25
	20.11	Seal / O-Ring Installation	20 - 26



## Table of Illustrations - Drawings

---

Drawing Name	Drawing Number	Page No.
<b>TAB 7: CPU Card</b>		
CPU BOARD ASSEMBLY	78101-10-AP-14-1	7 - 7
CPU BOARD ASSEMBLY	78101-10-AP-01-1	7 - 9
<b>TAB 8: I/O Card</b>		
I/O BOARD ASSEMBLY	78101-10-AP-02	8 - 3
<b>TAB 9: SCC Card</b>		
SERVO CONTROL CARD ASSEMBLY	78101-10-AP-03-1	9 - 5
<b>TAB 10: PPI Card</b>		
PPI CARD ASSEMBLY	78101-10-AP-06-1	10 - 7
PPI CARD REVISION HISTORY	78101-10-RH-01	10 - 9
<b>TAB 11: Control Panel</b>		
DISPLAY BOARD CARD ASSEMBLY	87101-01-AP-01	11 - 3
<b>TAB 12: Power Drive Circuit</b>		
POWER DRIVE CARD ASSEMBLY	78101-06-AP-09	12 - 15
POWER DRIVE CARD ASSEMBLY	78101-06-AP-06	12 - 18
POWER DRIVE RELAY CARD ASSEMBLY	85106-04-AP-02	12 - 19
<b>TAB 13: MSD Card</b>		
MECHANICAL STOP POWER CARD ASSEMBLY	78101-04-AP-05	13 - 3
MECHANICAL STOP POWER CARD ASSEMBLY	78101-04-AP-05	13 - 5
<b>TAB 14: Power Supply Circuit</b>		
POWER SUPPLY CARD ASSEMBLY	78101-07-AP-03	14 - 3
POWER SUPPLY CARD ASSEMBLY	78101-07-AP-01	14 - 5

Drawing Name	Drawing Number	Page No.
<b>TAB 15: Injector Head</b>		
MCT PLUS FL INJ. HEAD MAJOR PARTS DETAIL	87101-15-LD-04	15 - 13
MK V/MCT PLUS INJECTOR HEAD CARD ASSY.	78101-15-AP-01	15 - 16
MCT PLUS INJECTOR HEAD CARD ASSEMBLY	87101-15-AP-05	15 - 18
MCT PLUS INJECTOR HEAD MJR. PARTS DTL.	87101-15-LD-03	15 - 20
MCT FL INJECTOR HEAD MJR PARTS DETAIL	87101-03-LD-01	15 - 23
MCT INJECTOR HEAD MAJOR PARTS DETAIL		15 - 26
INJECTOR HEAD INTERCONNECT CARD ASSY.	78101-04-AP-07	15 - 27
INJECTOR HEAD PISTON/PLUNGER CARD ASSY.	78101-04-AP-09	15 - 28
INJECTOR HEAD PLUNGER POS. CARD ASSY.	78101-04-AP-02	15 - 30
INJECTOR HEAD MECH. STOP INDICATOR CARD	78101-04-AP-03	15 - 31
INJECTOR HEAD SCALE INDICATOR CARD ASSY.	78101-04-AP-04	15 - 32
MCT / MCT PLUS HEAD CABLE WIRING		15 - 33
<b>TAB 16: ISI 100</b>		
ISI 100 CARD ASSEMBLY	78101-10-AP-09	16 - 9
ISI 100 CARD ASSEMBLY	78101-10-AP-09	16 - 11
<b>TAB 17: Remote Panel</b>		
REMOTE PANEL CONTROLLER CARD ASSY.	85102-01-AP-02-1	17 - 7
REMOTE PANEL CONTROLLER CARD ASSY.	85102-01-AP-01-1	17 - 10
REM. PAN. CONTROLLER/KEYBOARD ASSY.	87101-15-AM-06	17 - 12
REMOTE PANEL ASSEMBLY	87101-01-AM-03	17 - 14
REMOTE PANEL STARTSWITCH CABLE ASSY.	85102-02-AC-01	17 - 15
REM. PAN. ARMED INDICATING HANDSWITCH	85102-02-AC-02	17 - 16
<b>TAB 23: Drawings</b>		
CABLING DIAGRAMS		23 - 1
MK V ANALOG INTERCONNECTION DIAGRAM	78101-00-WI-02	23 - 3
MOTHER BOARD ASSEMBLY (MN)	85106-07-AP-01	23 - 4
MN MCT/CRC MOTHERBOARD ASSY.	85106-07-AP-02-1	23 - 6
MOTHER BOARD ASSEMBLY (CRC)	85106-03-AP-04	23 - 8
ISI ISOLATION CARD ASSEMBLY	85106-03-AP-03	23 - 9
MCT PLUS CONNECTOR LOCATION DETAIL	87101-15-LD-01	23 - 10
MCT PLUS LABEL LOCATION DETAIL	87101-15-LD-02	23 - 11
MCT LABEL LOCATION DETAIL	87101-00-LD-02	23 - 13
MCT PLUS WORKSTATION 1 ASSEMBLY	87101-15-AM-01	23 - 15
MCT PLUS WORKSTATION 2 ASSEMBLY	87101-15-AM-02	23 - 20
MCT PLUS WORKSTATION 3 ASSEMBLY	87101-15-AM-03	23 - 25
MCT WORKSTATION 1 ASSEMBLY	85106-00-AM-01	23 - 30
MCT WORKSTATION 2 ASSEMBLY	85106-00-AM-02	23 - 35
MCT WORKSTATION 3 ASSEMBLY	85106-00-AM-03	23 - 40
MCT PLUS CHASSIS WIRING DIAGRAM	87101-15-WI-01	23 - 45
CRC CHASSIS WIRING DIAGRAM	85106-02-WI-01	23 - 46

## SECTION 2 - PREFACE

### Copyright Notice

Copyright © 1993 by Medrad, Inc. All rights reserved. No part of this manual may be reproduced in any form without prior written permission of Medrad.

FluiDot is a registered trademark of Medrad, Inc.

**Restricted sale:** Federal (USA) law restricts the sale of this device to, on, or by the order of a physician.

### Disclaimers

Medrad makes no warranties on the contents of this manual, and specifically disclaims any implied warranties of merchantability or fitness for any purpose.

Medrad reserves the right to change specifications and the contents of this manual without obligation.

**External wiring modification disclaimer:** Medrad cannot be liable for damage done by unauthorized interfaces or unauthorized modifications to the MCT/MCT Plus. Such unauthorized action could jeopardize injector operation, safety, or reliability. Medrad will review and, if appropriate, authorize modifications, repairs and special interfaces. To obtain this service, call or write to Medrad's Factory Service Department (address, phone and telex listed below).

All drawings in this manual are for reference purposes only, and may not reflect the construction of units produced prior to the publication of this manual. Reproduction quality of these drawings may have been effected by the level of reduction required. Call Medrad Factory Service if assistance in drawing interpretation is required.

### Effectivity

This manual applies to the MCT and MCT Plus Injection Systems in the following configurations: Control Room Console MCT 305/310, and MCT 305P/310P, Remote Panel Systems RPS 350 and 350P.

*Injection Systems equipped to operate at voltage requirements of 95-125 VAC, will be designated throughout this manual as -1 units. Injection Systems equipped to operate at voltage requirements of 190-250 VAC, will be designated throughout this manual as -2 units.*

If you have any comments on this manual, please call or send the details to:

Technical Support Department  
MEDRAD, INC.  
P.O. Box 780  
Indianola, PA 15051-0780

(412) 967-9700  
1-800-MEDRAD-1  
(1-800-633-7231)  
TLX: 199114  
FAX: (412) 963-0859 USA

If you encounter technical problems with the MCT or MCT Plus Unit or want to order parts, contact your local dealer, or:

#### Domestic:

Medrad Factory Service  
MEDRAD, INC.  
240 Alpha Drive  
Pittsburgh, PA 15238-2870  
(412) 967-9700  
1-800-MEDRAD-S  
1-800-633-7237  
TELEX: 199114  
FAX: (412) 963-1964

#### International:

MEDRAD INTERNATIONAL, B.V.  
Postbus 3084  
6202 NA Maastricht  
The Netherlands  
Phone: (31) (0) 43 640808  
FAX: (31) (0) 43 650020

## Preface

### ABOUT THIS MANUAL

#### Purpose

This is the Service Manual for MCT and MCT *Plus* CT Injection Systems. This manual is intended to provide instructions for servicing the injectors as safely and thoroughly as possible.

#### Important Safety Notice

The information in this manual is intended for people with adequate backgrounds and experience in electronics and electromechanical devices. Any attempt to repair a sophisticated medical device such as the injector may result in personal injury, property damage, or patient injury. Medrad and other distributors of the injector cannot be responsible for the interpretation of this information, nor can they assume any liability in connection with its use.

#### Intended purpose of the MCT/MCT *Plus* Injector

This device is sold for the specific purpose of injecting contrast medium into humans and animals for diagnostic studies. Do not attempt to use the injector for any other purpose such as chemotherapy or drug infusion.

#### Audience

This manual is written for those qualified to service the injector, whether they be Medrad Service Personnel, Laboratory Service Technicians, or Medrad International Dealers.

#### Software

The information in this operation manual applies only to MCT and MCT *Plus* CT Injectors equipped with software version CU2.1 or below, and Remote Panel software version EP2.0 or below.

**MCT *Plus* units designed as Class 1 Type BF equipment and are identified throughout this manual as -I units.**

MCT *Plus* rack mount units manufactured after July 30, 1991, equipped to operate at 220-250 VAC, are designed to, and in compliance with, RF emission standard for Class 1 Type BF equipment, and meet IEC Standard 601-1 for Class 1 Type BF equipment. Class 1 Type BF equipment is identified by the following symbol:



## WARNINGS and CAUTIONS

PLEASE READ CAREFULLY

#### WARNINGS:

**Ensure that FluidDot Labels are present and clearly visible on all pressure jackets on non-Front Load Systems. FluidDot labels are intended to help avoid air embolization, and should be clearly visible on the pressure jacket.**

**Hazardous voltages exist within the MCT/MCT *Plus*. To avoid shock, the injector should only be opened by qualified service personnel. Always disconnect the injector from line power before cleaning or attempting to perform any maintenance.**

**Possible explosion hazard. Do not use the MCT/MCT *Plus* in the presence of flammable anesthetics.**

## Preface

**Do not use the injector if a system monitor message is displayed.** System monitor messages are displayed on the Control Panel, and are preceded by "SC". This type of message indicates serious system malfunctions. Do not use the injector until the cause of these messages can be determined. Contact a Medrad Factory Service Representative for assistance.

**Do not use the injector if any worn cords or cables are detected. Personal Injury or equipment damage/malfunction may result.** Contact the Medrad Factory Service Department for replacement information.

**Check for proper voltage and frequency before plugging the injector into an electrical outlet.** Check the voltage and frequency marked on the serial number tag on the back of the injector. Typical markings are 110 VAC 60 Hz or 220 VAC 50 Hz. Be sure that the outlet providing power to the injector supplies this voltage and frequency. Personal injury or equipment damage/malfunction may result.

**Do not apply voltage to external start lines.** If the injector is being started by an external start connection, do not apply voltage to the start circuit. Provide only a switch closure. If the external circuit contains excessive line frequency noise or transients, personal injury or equipment damage/malfunction may result.

**Use remote disarm switch with external starts.** If the injector is wired to be started externally, the operator does not have full control of the injector with the standard handswitch. The handswitch may be used to start the injector as normal, but once the external start contacts are closed, releasing the handswitch may not abort the injection. Therefore, we suggest providing an external disarm switch if an external start circuit is used. However, any injection may be aborted by pressing a Control Panel key.

**Possible injury to the patient and technician can occur if the injector head stand is moved by pulling on the injector head.** This practice will cause the pivot pin the wear and may cause the injector head to fall.

**Do not position the system by pulling on the injector head or cabling if using an Overhead Counterpoise System.** If the "Head Mount Knob" is not tight and the user pushes up on the head to raise the system to its highest position, the potential exists for the injector head to lift out of its socket and fall. This could cause user or patient injury. Use the Overhead Counterpoise "Vertical Support Arm Assembly" to move and position the head. Also, ensure that the "Head Mount Knob" which secures the injector head to the Vertical Support Arm is tight at all times. Refer to the Overhead Counterpoise System Operation, Installation and Service Manuals for additional Warning and Caution information.

### **CAUTIONS:**

**Electrostatic Discharge (ESD).** Static protection practices must be followed when servicing any component of this system. Before touching any of the circuit cards in the injector, discharge yourself to grounded metal. If the memory chips are to be shipped, place the components in conductive carriers (such as that used by Medrad).

**Do not short PPI Card batteries.** When handling the PPI Card, or when measuring or replacing the batteries on the PPI Card, be careful not to short the batteries, which can be damaged by even a momentary short. Do not set the card on any metal surface; this will short out the connections on the bottom of the card, and discharge the batteries.

**Do not remove or replace circuit boards with power applied to the injector.** Sensitive circuits on the cards can be damaged by abrupt interruption or application of supplies. Remove power before servicing.

**Total length of the head extension cable should not exceed 100 feet (30.5 m).** If cables are connected in series or their total length exceeds 100 feet (30.5 m), injector performance and specifications will be adversely affected.

## Preface

**CAUTION:** Do not manually advance the piston plunger if the injector is armed. If the standard procedure is to manually move the piston plunger to ensure that no air is in the connector tube perform this step before arming. Movement of the plunger up to 2 ml in the forward or reverse direction after the injector is armed will result in an inaccurate volume delivery. The injection bolus volume will be program volume minus the preliminary plunger movement. (i.e., If program volume is 10 ml and the plunger is advanced 2 ml while the injector is armed, the actual volume delivered to the patient during the injection will be 8 ml. This volume displacement will not be displayed when Status values are displayed.)

**Allow injector temperature to stabilize before use.** When the injector is exposed to an extreme temperatures change (heat or cold), allow it to stabilize to room temperature before servicing.

**Use only accessories and options designed specifically for the MCT and MCT Plus.** For proper operation, use only accessories and options provided by Medrad which are designed specifically for the MCT and MCT Plus system. This ensures compatibility with the injector. Do not use an accessory or option designed for another system on the MCT and MCT Plus.

**Perform regular preventive maintenance.** To ensure that your MCT/MCT Plus remains properly calibrated, and that all primary and backup circuits are functioning properly, yearly preventive maintenance is suggested. Annual preventive maintenance (offered in the US and Canada) is not covered by the new machine warranty. Contact your local Medrad Factory Service Representative for details.

### VOLUME REMAINING Keys and Display

If Programmed Injection Volume exceeds the volume remaining in the syringe, the Insufficient Volume Light illuminates. Press the Insufficient Volume key to override, then proceed with arming sequence. The Insufficient Volume Light flashes, and the injector can be armed.

- Volume Remaining - ml
- Insufficient Volume Key
- Insufficient Volume Light

### FIRST PHASE Keys and Displays

The First Phase Block displays the programmed Flow Rate and Volume for the first phase of an injection. To change either of these values, use the respective arrow keys located below each display.

- First Phase Flow Rate Display - ml/sec
- First Phase Volume Display - ml
- First Phase Indicator

### ARMED Keys and Displays

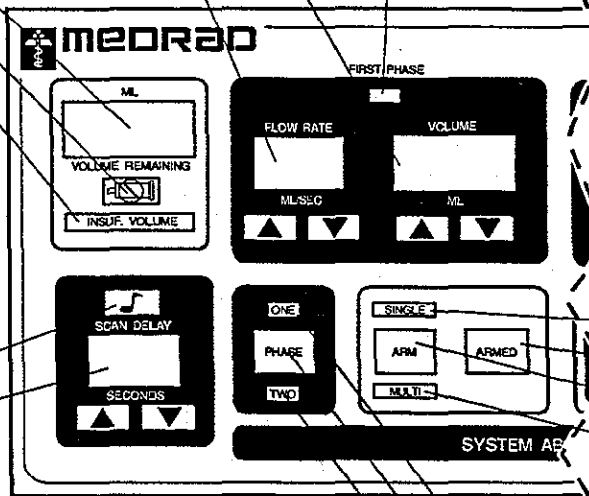
**ARM:**  
Press the Arm key. The OP 100 message is displayed - **Confirm no air in fluid path.** If no air is detected, press the Arm key again. The Armed Indicator Light flashes, and the injector is armed for a Single injection.

**MULTIPLE ARM:**  
To arm the MCT or MCT Plus for a multiple injection, press the Armed key while the injector is armed. The Multiple Injection indicator will light and the injector will remain armed after each injection as long as sufficient volume remains in the syringe.

### SCAN DELAY Keys and Display

External scanning equipment can be delayed by entering a time value in the Scan Delay Window. Use the Scan Delay arrow keys to select a delay time - up to 99 seconds. The Scanner Tone Alert key is a toggle switch which activates a tone alarm whenever scanning begins.

- Scanner Tone Alert Key
- Scan Delay - Sec.



- Single Injection Light
- Armed Indicator Light
- Arming Key
- Multiple Injection Light

### ARROW Keys

The arrow keys below the Scan Delay, Flow Rate, and Volume Displays all perform in a similar manner. To increase a value press the UP arrow key; to decrease a value press the DOWN arrow key. If any arrow key is held in for longer than 2 seconds, the value displayed will increase or decrease at 5 counts/second.

### PHASE SELECTION Keys and Display

If a single phase injection is currently displayed it can be changed to a bi-phasic injection by pressing the Phase key. The Phase Indicator Lights change from ONE to TWO and values in the Second Phase block are displayed.

- Single Phase Injection Light
- Injection Phase Key
- Bi-phasic Injection Light

### PROGRAMMING Keys

**Program Selection**

- Press the Program Select key until the desired Program Number is lit
- Program Number 6 is used for Multi-Level programming

**Program Change**

- When any program is selected, change any programmed value with arrow keys
- Program Number light will flash left to right

**Reprogramming**

- Press Reprogram Key
- All values flash
- Change values with arrow keys

**Program Store**

- Press Program Lock key

**Total Volume Programmed**

- Press and hold Program Lock/Total Volume Programmed key when not reprogramming
- Total Volume Programmed will be displayed in Volume Delivered Window

Program Numbers (1-6)

Program Select Key

Reprogram Key

Program Lock Key  
or Total Volume Programmed

### SECOND PHASE Keys and Displays

The Second Phase Block displays the programmed Flow Rate and Volume for the second phase of an injection. Use the arrow keys to change the values displayed.

Second Phase Flow Rate Display - ml/sec

Second Phase Volume Display - ml

Second Phase Indicator

### VOLUME DELIVERED and TOTAL VOLUME Window

The Volume Delivered Window displays the Total Volume Delivered in ml. This window also displays the Total Volume Selected. Press and hold the Program Lock key to check Total Volume Selected.

Volume Delivered - ml  
or Total Volume Selected

### INJECTION DURATION and PROGRAM LEVEL Window

This window displays Injection Duration. When a Multi-Level program is displayed, the Injection Duration window displays the current Program Level. While injecting, this window displays Injection Duration. While injecting in the multi-level mode, the PL number is briefly displayed.

Injection Duration  
or Program Level

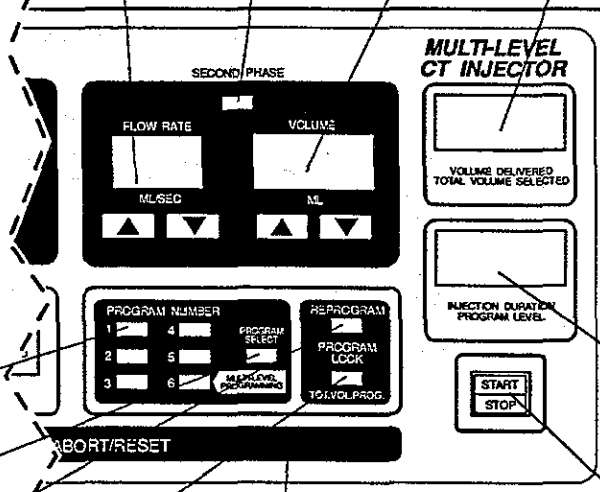
### START/STOP Key

Functions the same as the Hand Start Switch. When the MCT/MCT Plus is armed for a Single injection, use the Start/Stop key to begin the injection, or HOLD an injection during a procedure. If the injector is armed for Multiple injections, HOLD feature is disabled.

### SYSTEM ABORT /RESET

The System Abort/Reset key is used to abort an injection, or reset the MCT after an OP message. When power is first applied, press and hold this key immediately after the countdown to display the version number of the software in your unit.

**TO ABORT: PRESS ANY KEY.**





**OP 100 CHECK SYRINGE FOR AIR, THEN ARM**

This message is displayed after arming the unit for an injection. It is intended to remind the operator about the dangers of injecting an air embolism. Check the syringe and connector tubing for air, then press the ARM key again to proceed.

**OP 101 UNIT DISARMED (applies to integral or remote panel keyboard)**

An error condition has occurred after the unit was armed. Three conditions that may have caused this message are: forward or reverse load button on Injector Head was activated, a key on MCT / MCT Plus panel, (other than START/STOP), was pressed, or syringe turret was moved out of the injection position (non-Front Load systems only). Press the SYSTEM ABORT/RESET bar to continue.

**OP 102 CHECK TURRET, THEN RESET (Non-Front Load injector heads only)**

Turret is not in position for an injection. Re-position the turret, press the SYSTEM ABORT/RESET bar, then re-arm the unit.

**OP 102 MECHANICAL STOP PLATE IMPACTED, MANUALLY REVERSE PLUNGER**

**(Front Load injector heads only)** If the Mechanical Stop microswitch is activated before the injector is armed, this message code will appear. Press the SYSTEM ABORT/RESET bar, then use the manual knob to reverse the plunger to continue.

**OP 103 FORWARD MOTOR MOVEMENT**

After the unit was armed, the manual load knob on the Injector head was moved in the forward direction. Press SYSTEM ABORT/RESET bar to continue.

**OP 104 REVERSE MOTOR MOVEMENT**

After the unit was armed, the manual load knob on the Injector head was moved in the reverse direction. Press SYSTEM ABORT/RESET bar to continue.

**OP 105 REMOTE CONTROL PANEL COILED START SWITCH SHORTED**

This message alerts the operator that the Remote Control Panel Coil Start Switch was depressed or held down prior to or during the arming sequence. Press the SYSTEM ABORT/RESET bar then re-arm.

**OP 106 MAIN UNIT START SWITCH SHORTED**

This message indicates a Start Switch or Coiled Start Switch at the main control unit is depressed during the arming sequence. OP 106 does NOT apply to the Control Panel START/STOP button. Un-short the switch, press the SYSTEM ABORT/RESET bar, then re-arm.

**OP 107 CHECK INJECTOR HEAD CABLE (PLUNGER)**

The feedback signal from the plunger positioning potentiometer is not within allowable limits. Check cable connection J21 on back of MCT / MCT Plus main console for a secure and true connection. Press the SYSTEM ABORT/RESET bar and re-arm the MCT / MCT Plus.

**OP 108 STALL TIME EXCEEDED**

Some obstruction has stopped the plunger in its course of movement. If this condition persists longer than 3 seconds, the OP 108 message code is displayed. Check connector tubing for blockages, sharp bends, or kinks. If stopcocks were used ensure they are opened. Make sure the plunger path is free from obstacles. Press the SYSTEM ABORT/RESET bar to proceed.

**OP 109 EXTERNAL EQUIPMENT DISARMING THE INJECTOR**

The Imaging System Interface has signaled a problem, thereby rendering the MCT / MCT Plus Injector incapable of being armed. Determine which piece of external equipment is causing the problem, correct it, and press the SYSTEM ABORT/RESET bar to proceed.

**OP 110 INSUFFICIENT VOLUME REMAINING IN THE SYRINGE**

Indicates to operator that the programmed volume of the current injection exceeds the actual volume remaining in the syringe. Press the SYSTEM ABORT/RESET bar to erase OP 110 message. Notice the lighted Insufficient Volume indicator in the Volume Remaining area of the MCT / MCT Plus panel. To override the arming lockout and avoid the OP 110 message, press the Syringe key just above the Insufficient Volume indicator. This causes the Insufficient Volume light to flash. The unit can now be armed.

**OP 111 FORWARD/REVERSE LOAD SWITCH DEPRESSED**

Indicates that the forward or reverse load button on the Injector Head was being pressed while the system was being armed. Release the forward or reverse load button, press the SYSTEM ABORT/RESET bar, and repeat the MCT / MCT Plus arming sequence.

**OP 112 CHECK INJECTOR HEAD CABLE (MECHANICAL STOP)**

The feedback signal from the mechanical stop positioning potentiometer is not within allowable limits. Make sure the cable connecting the Injector Head to the main console is connected securely. Press the SYSTEM ABORT/RESET bar to continue.

**OP 113 (MESSAGE NOT USED)**

**OP 114 IMAGING SYSTEM INTERFACE, START COMMAND ACTIVATED**

This message appears when an external device attempts to activate the injector prior to completion of the arming sequence. The external device is connected to J20 or J40 on the back of MCT / MCT Plus Main Console. Check for a secure cable connection. Press the SYSTEM ABORT/RESET bar to clear message.

**OP 115 HOLD TIME EXCEEDED**

The MCT / MCT Plus Injector has remained in the "Hold" mode for longer than 5 minutes. To avoid this problem, the operator must resume a held injection before 5 minutes elapse. Press the SYSTEM ABORT/RESET bar to proceed.

**OP 116 PREMATURE START SWITCH RELEASE**

Whenever the START/STOP switch or external hand start switch is held in longer than 5 seconds to start the injection, it must be held in until the injection is complete. The OP 116 message is displayed if these conditions are not met. Press the SYSTEM ABORT/RESET bar to clear the message, then repeat normal MCT / MCT Plus arming sequence.

**OP 117 MOVE PLUNGER 2ML BEFORE ARMING**

This message was added in software version 1.4 and above as part of additional self-diagnostic tests the MCT / MCT Plus performs when power is first applied. In order to check the circuitry monitoring plunger movement - the plunger must be moved. Normally, this requirement is fulfilled when the syringe is loaded and this message is not displayed. However, if power is interrupted between procedures and an attempt is made to re-arm the unit before the plunger has been moved, the OP 117 message is displayed. Move the plunger a minimum 2 ml FWD or REV before proceeding.

## SECTION 4 - PREVENTIVE MAINTENANCE

This section contains recommended procedures for preventive maintenance of MCT and MCT *Plus* injection systems in the following configurations: Control Room Console MCT 305/310, and MCT 305P/310P, Remote Panel Systems RPS 350 and 350P. Following these simple procedures can:

- Ensure the best performance from your injector.
- Help avoid major problems through regular inspections and checkout.

### Recommended Preventive Maintenance

Preventive maintenance of the injection system should consist of four procedures: Inspection, cleaning, performance checks, and electrical leakage checks. This section contains guidelines, recommended methods, and expected results for each of these procedures:

1. *Inspection:* This first step should encompass inspection of the entire system, looking for obvious signs of damage, such as; cracks in the housing, frayed or worn cords, and contrast spills that may have leaked into the injector.
2. *Cleaning:* This cleaning procedure involves thorough cleaning of the console and head to remove any deposits of contrast medium. If any substances have leaked into any part of the unit, the subassembly should be disassembled and thoroughly cleaned.
3. *Performance Check:* A complete functional performance checkout of the injection system.
4. *Electrical Leakage Check:* To insure the safety of the patient and hospital personnel in injector operations, an electrical leakage check is part of this preventive maintenance procedure.

### Recommended Preventive Maintenance Schedule

Although Medrad recommends the following schedule, your individual maintenance schedule depends upon how your injector is used; the type of procedures performed, and frequency of use. The following guidelines represent a suggested maintenance schedule:

#### Monthly:

Once a month, the system should be inspected and cleaned thoroughly.

#### Semi-Annually:

Every six months, both the Electrical Leakage and Performance Checks should be performed.

**NOTE:** Local regulations or hospital protocol may require electrical leakage checks at *more frequent* intervals.

If this applies, local regulations for leakage tests must be followed.

#### Annually:

Medrad recommends that a complete system calibration and performance checkout, by a qualified Medrad Service Representative, be performed once a year. Contact your local Medrad office for details.

In the United States and Canada, the Medrad Service Department offers Preventive Maintenance Programs. These annual programs greatly assist in the maintaining of the accuracy and reliability of the unit, and can extend the life of the injection system. Contact Medrad for details. For information in Europe, contact Medrad International B.V.. Refer to Section 2 for address and telephone numbers

# Preventive Maintenance

## Inspection Procedures

The following procedures are recommended for the Injector Head, Control Unit, Remote Panel (if the system is so equipped), and mounting systems.

If defects are detected, either repair the system, or call Medrad for service. Do not use the unit until the problem has been corrected.

### Injector Head Inspection

1. Inspect the injector head case for any cracks that could allow fluid to leak inside, or weaken the case.
2. Inspect the injector head cable for cuts, cracks, or worn areas.
3. Inspect the injector head connector for cracks, loose pins, or a loose strain relief.
4. If applicable, inspect the Pressure Jackets for small cracks, crazing, or discoloration. As a matter of preventive maintenance, we recommend replacement of the pressure jackets annually.
5. Inspect the syringe heater and cable; look for cuts, cracks, or any worn areas.
6. Inspect the Control Panel for cuts or cracks that could allow fluid to leak inside.

### Control Unit Inspection

1. Inspect the Control Unit case for cracks that could allow fluid to leak inside, or weaken the case.
2. Inspect the power cord for cuts, cracks, or worn spots. Inspect the plug for cracks, loose prongs, loose wires, or a loose strain relief.
3. Inspect the handswitch and cord: Look for cuts, cracks, or worn spots in the cable; look for cracks and loose parts in the switch and housing.
4. Inspect the film changer connector and cord: Look for cuts, cracks, or worn spots in the cable; look for loose pins, or a loose strain relief on the connector.
5. Inspect any other cables connected to the control unit: Look for cuts, cracks, or worn spots in the cables; look for cracks, loose pins, or loose strain reliefs on the connectors.
6. Inspect the Control Panel for cuts or cracks that could allow fluid to leak inside.

### Remote Panel Inspection

1. Inspect the remote unit case for cracks that could allow fluid to leak inside, or weaken the case.
2. Inspect the interface cable for cuts, cracks, or worn spots. Inspect the connector for cracks, loose prongs, or a loose strain relief.
3. Inspect the handswitch, cord, and connector: Look for cuts, cracks, or worn spots in the cable; look for cracks and loose parts in the switch, housing, and connector.
4. Inspect the Control Panel for cuts or cracks that could allow fluid to leak inside.

# Preventive Maintenance

## Mounting Systems Inspection

### Pedestal Mount - KMA 320/KMA 320-T (See Figure 4.1)

- Inspect the stand, base, and support arm for cracks and other defects that could weaken the structure.
- Ensure that all mounting bolts and screws are secure.
- Ensure that the casters roll smoothly, with no binding or scraping.
- Ensure that the locking mechanism on all locking casters is functional.

### Injector Support Arm

The injector head support arm is typically associated with the free standing injector head stand (KMA-320 or KMA-320T).

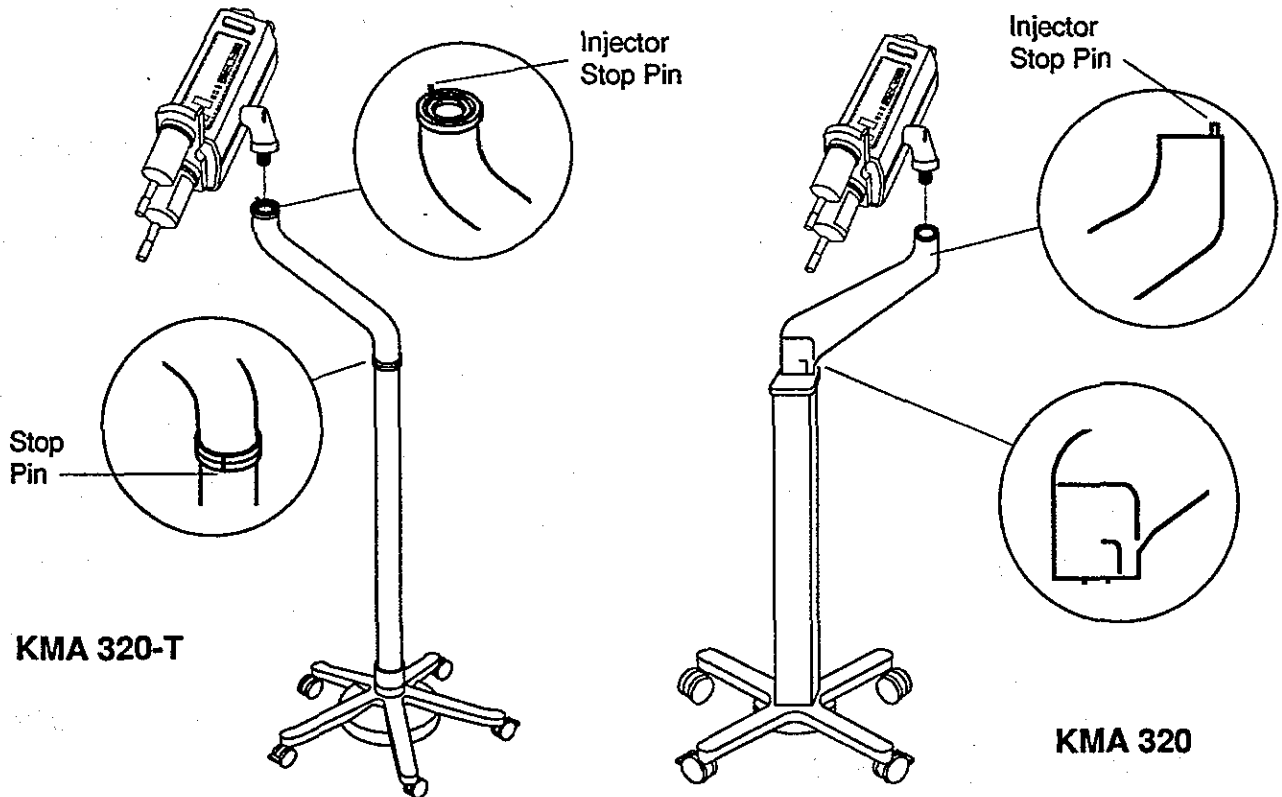


Figure 4.1: Head Support Arm Inspection Points

# Preventive Maintenance

The injector head pivot knuckle should sit flat and rotate smoothly on the support arm. The injector head should not rotate more than 360 degrees. If the head is able to rotate more than 360 degrees, carefully lift the head from the support arm and check for the presence of the injector head stop pin. If the pin is not in place, the head cable can wrap around the arm causing damage to the head cable and possibly lifting the injector head off of the support arm.

Additionally, the support arm should be inspected for cracks where the head rests on the arm and at the base where it is attached to the injector head stand. If a crack is present, the crack can lengthen, compromising the integrity of the arm. If any cracks are found, Medrad Service should be contacted.

## Overhead Counterpoise Systems, Wall and Ceiling Mount (See Figures 4.2 and 4.3)

Medrad provides both Ceiling and Wall mounted overhead counterpoise systems. Each system consists of an articulating vertical and horizontal arm that require regular maintenance inspections.

On a daily basis the black head mount knob should be checked to ensure that it is fully tightened. This knob is intended to secure the injector head on the vertical support arm and minimize the chance of the injector head being accidentally lifted from the support. Rotation of the injector head should be accomplished by turning the entire vertical support arm and not by rotating the injector head. Additionally, the following items should be regularly inspected:

**Long Pin and Spring Pin - (Figure 4.2)** As a safety back up to the retaining clip, there is a long pin which passes through the vertical and horizontal support arm where they join. The through-hole at the bottom of the long pin should have a spring pin through it. The spring pin must be in place for the safety backup to be effective. On early OCS systems the long pin and spring pin will not be present. These early systems do not have a through hole for the long pin.

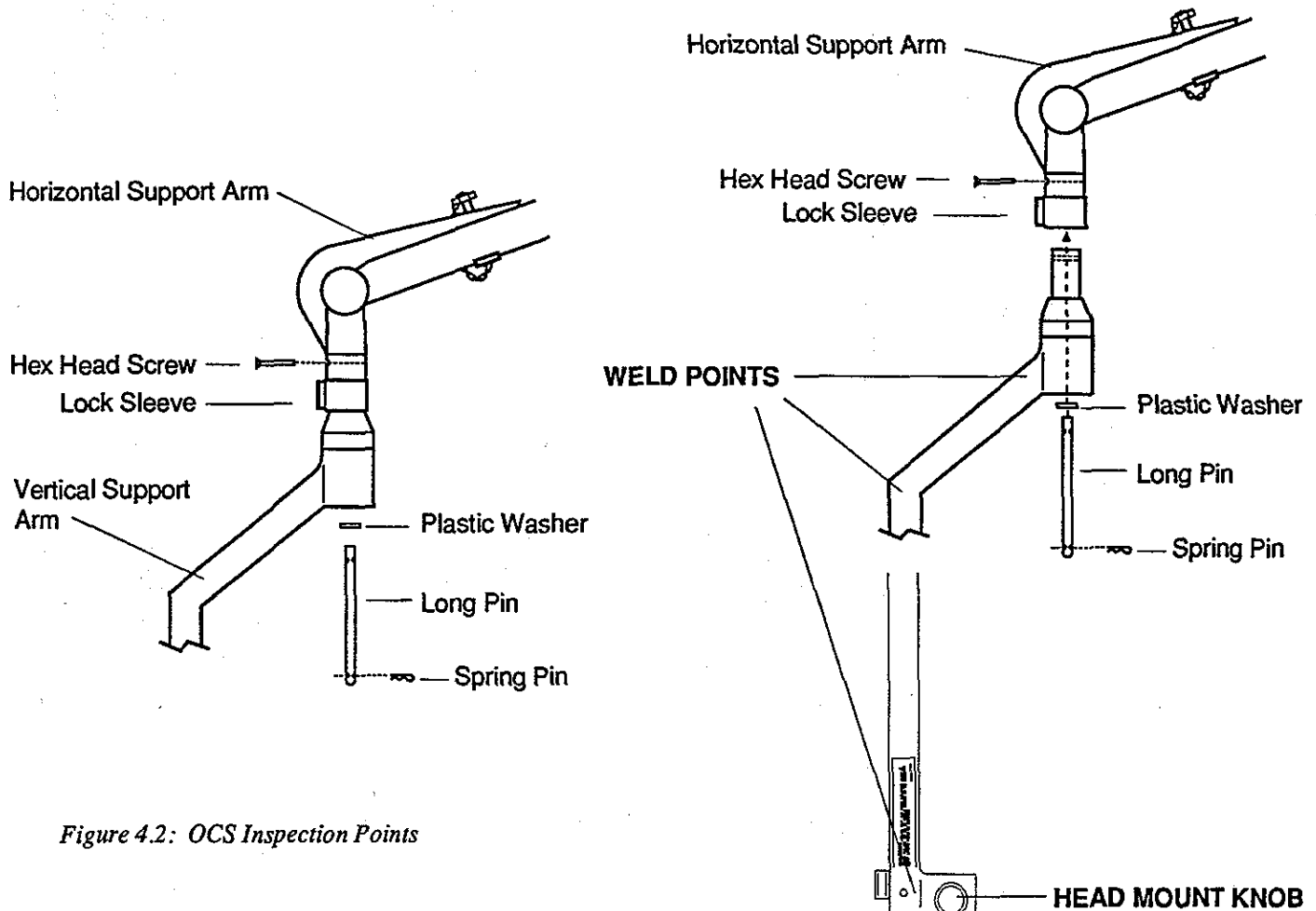


Figure 4.2: OCS Inspection Points

## Preventive Maintenance

**Vertical Support Arm - (Figure 4.2)** The vertical support arm weld joints should be inspected for cracks. Cracks can lead to failure of the support arm.

**Retaining Clip and Plastic Locking Sleeve - (Figure 4.3)** The plastic lock sleeve should be inspected regularly to ensure it is not cracked or showing other signs of damage. It ensures the retainer clip remains in place. If the lock sleeve is missing or not in place, the retainer clip may come out resulting in the vertical arm releasing from the horizontal support. A damaged or missing plastic lock sleeve should be replaced immediately.

**Horizontal Arm - (Figure 4.3)** The horizontal arm should be inspected overall for cracks and any other loose or damaged parts. Any missing or damaged parts should be replaced.

**Ceiling and Wall Mount Inspection - (Figure 4.3)** The ceiling or wall mount should be inspected for the following:

- All bolts and set screws must be fully tightened and secure.
- Ensure the horizontal arm is level. An unlevel horizontal arm may indicate the mounting pin is bent or fractured or the mount is not properly secured to the wall or ceiling.
- The horizontal and vertical arm assembly should remain in the position in which it was placed. Any drifting, drooping or swaying of the arm may indicate damage to the mounting pin. Examine the mounting pin and if loose or damaged, replace immediately.
- The horizontal arm should be inspected overall for cracks and any other loose or damaged parts. Any missing or damaged parts should be replaced.

Refer to the OCS Installation and Service Manual for additional information or contact Medrad Service.

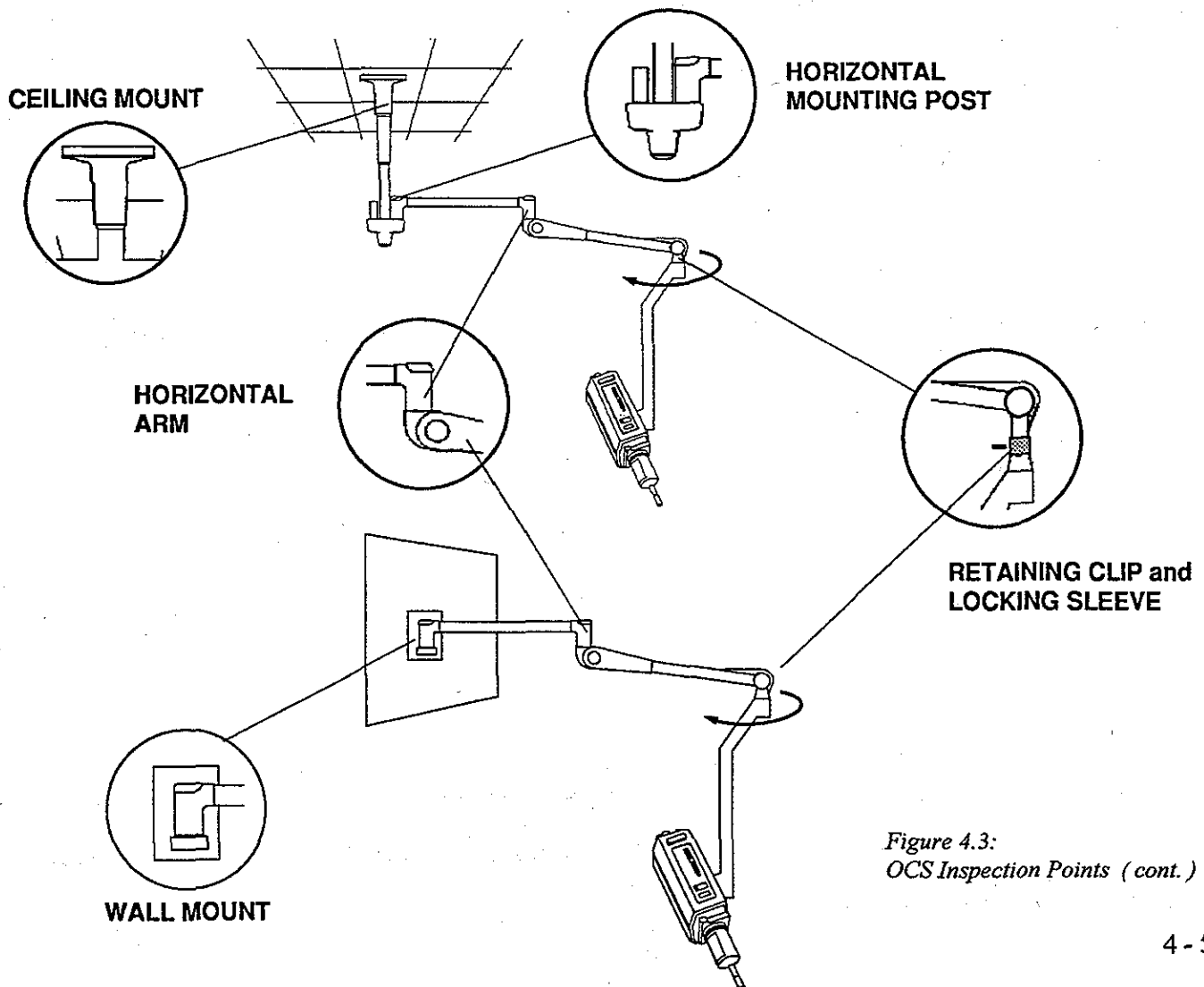


Figure 4.3:  
OCS Inspection Points (cont.)

# Preventive Maintenance

## Cleaning Procedure

This procedure is recommended for the cleaning of a MCT/MCT Plus system. If not removed, contrast medium can interfere with proper operation.

**WARNING:** Remove the power cord from the power source before cleaning any part of the system. Failure to do this could result in the exposure of lethal voltages, causing injury or death.

- Using warm water and soft cloths, clean the injector console and head thoroughly to remove contrast medium or other deposits.
- *Non-Front Load head* Remove the turret, then clean the front casting where the piston moves in and out; Clean the turret and its pivot point; Remove any contrast from the entire length of the piston rod. Replace the turret.
- *Front Load head* Remove the Adaptor Plate and clean. With a warm cloth, clean the front plate and the entire length of the piston rod. Check the O-ring on the Front Plate for cuts, and to ensure that it is seated properly in the groove. Replace the Adaptor Plate.

**CAUTION:** Do not soak or immerse any parts of the injector in water. While cleaning the outside of the unit, avoid letting any water seep inside injector components.

- Allow time for the injector to dry thoroughly before using it. Do not dry by applying heat from any source.
- If contrast medium has leaked inside the injector head or console, the effected subassembly should be disassembled and cleaned. This cleaning procedure can be done in the field by trained Medrad Service personnel, or returned to Medrad Factory Service. If the cleaning will be performed in the field, do not disturb any internal wiring or components. Ensure that the system is completely dry before applying power.

## Electrical Leakage Check

To insure safe operation of the injector, a chassis electrical leakage check must be part of preventive maintenance.

Use a commercial leakage tester such as one of the following:

Manufacturer	Model
Bio-Tek Instruments, Inc. Electrical Safety Analyzer	Model 501
Dynatech Nevada, Inc. Electrical Safety Analyzer	PEI Model 2000A
Ohmic Instruments Co.	BET-300 Series or HSM-200 Series

### Procedure:

1. With the AC ground open, power applied, and the line at normal, leakage should be less than 100 micro amps (500 micro amps with -I units).
2. With the AC ground open, power applied and the line reversed, leakage should be less than 100 micro amps (500 micro amps with -I units).
3. Disconnect the leakage test box.



## SECTION 5 - SYSTEM SPECIFICATIONS

**Note:** Medrad reserves the right to change product designs and specifications at any time in the continuing effort to improve their products. The following specifications refer to new units operating at room temperature (approximately 21 degrees C or 70 degrees F).

### Control Console:

Dimensions:	Weight:	MCT	MCT Plus
19" W (48 cm)	<i>Console</i>	42 lbs (19.0 kg)	<i>Console</i> 44 lbs. (20.0 kg)
10 1/2" H (27 cm)	<i>Head</i>	18 lbs (8.0 kg)	<i>Head</i> 18 lbs. (8.0 kg)
13" D (32 cm)	<i>Total System</i>	60 lbs (27.0 kg)	<i>Total System</i> 62 lbs. (28.0 kg)

### Remote Panel:

Dimensions:	Weight:
8" H (19 cm)	5 lbs (2.3 kg)
17" W (42 cm)	without cable
4" D (10 cm)	

### Cord Lengths:

Power Cord	10 ft (3.0m) [MCT Plus -2 7ft (2.1m)]
Handswitch	25 ft (7.6m)
Injector Head	5 ft (1.5m)
Remote Panel	50 ft (15.2m) Std. 100 ft (30.5m) Max.

### Electrical Requirements:

Standby:	less than 1 Amp
Available:	95-105 VAC, 50 or 60 Hz, 4 A maximum
	105-120 VAC, 50 or 60 Hz, 4 A maximum
	120-125 VAC, 50 or 60 Hz, 4 A maximum
	190-210 VAC, 50 or 60 Hz, 2 A maximum
	210-240 VAC, 50 or 60 Hz, 2 A maximum
	240-250 VAC, 50 or 60 Hz, 2 A maximum

### Electrical Leakage:

Chassis	<100 microamps	[MCT Plus -2 < 500 microamps]
---------	----------------	-------------------------------

### Heat Dissipation:

Standby over 99% of the time):	300 BTU
Intermittently (<1% of the time):	1650 BTU
Average:	< 320 BTU

### Load Rate:

MCT Plus	Forward/reverse: 8 to 30 seconds to move the entire length of the plunger.
MCT	Forward/reverse: 15 to 30 seconds to move the entire length of the plunger.

### Syringe Heater:

37C (98F) nominal

### Syringe:

Disposable: 200ml

# System Specifications

## **Automatic Safety Stop:**

Sets and locks itself automatically before each injection.  
Sets and locks itself over the selected volume by 4.0 ml.

## **System Monitor:**

Redundant circuitry terminates injection if any of these conditions occur:

Overrate	10ml/s for 100 ms
Overvolume	5ml or 50% whichever is less
Overpressure	100 psi

## **Flow Rate:**

0.1 ml/s to 9.9 ml/s in 0.1 ml/s increments.

## **Total Volume Selection:**

1 ml to 200 ml in 1 ml increments

## **Pressure Limit:**

Factory Set at 300 psi. The pressure limiting for the MCT and MCT *Plus* is 300 psi, 21kg/cm<sup>2</sup>, or 2069 kPa.

## **Status Readouts:**

<b>Volume Remaining</b>	<b>Volume Delivered</b>	<b>Injection Duration</b>
Volume of Contrast available in the syringe prior to injecting.	Amount of Contrast delivered during the last injection cycle	Length of last injection in 1 sec. increments.
Read in 1 ml increments.	Read in 1 ml increments.	

## **Scan Delay:**

1 - 99 seconds

## **Program Memory:**

A total of 6 memory positions of: 6 single or bi-phasic injections  
*or*  
5 single or bi-phasic, and 1 multi-level injection

## **Level Control - MCT w/software CU 1.4 and below:**

Up to 2 phases (Volume and Flow Rate changes) per injection

## **Level Control - MCT/MCT *Plus* w/software CU2.0A and 2.1:**

Up to 2 phases (Volume and Flow Rate changes) per injection with exception of memory program 6. Program 6 can have up to 10 phases (or 5 program levels).

## **Analog Output Signals (optional; included with ISI 100):**

Injection Duration	Standby = < 0.1 VDC
	Injection = +5.0 VDC
Injection Volume	200 ml syringe: 0.6 VDC/10 ml
Injection Flow Rate	200 ml syringe: 2.0 VDC/10 ml /sec

# System Specifications

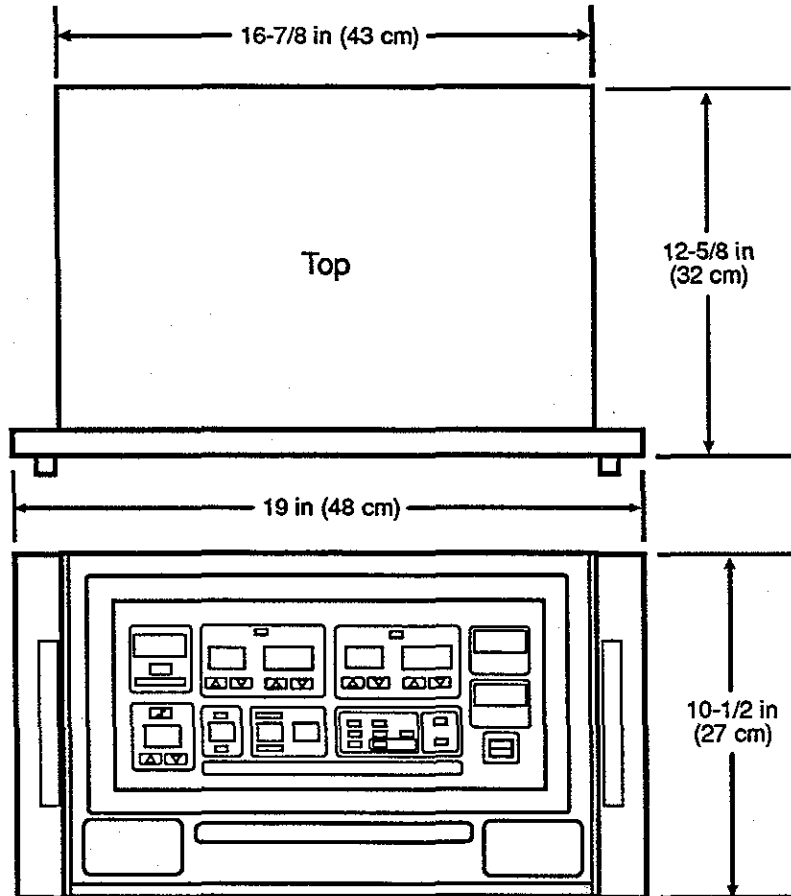


Figure 5.1:  
CRC Dimensions

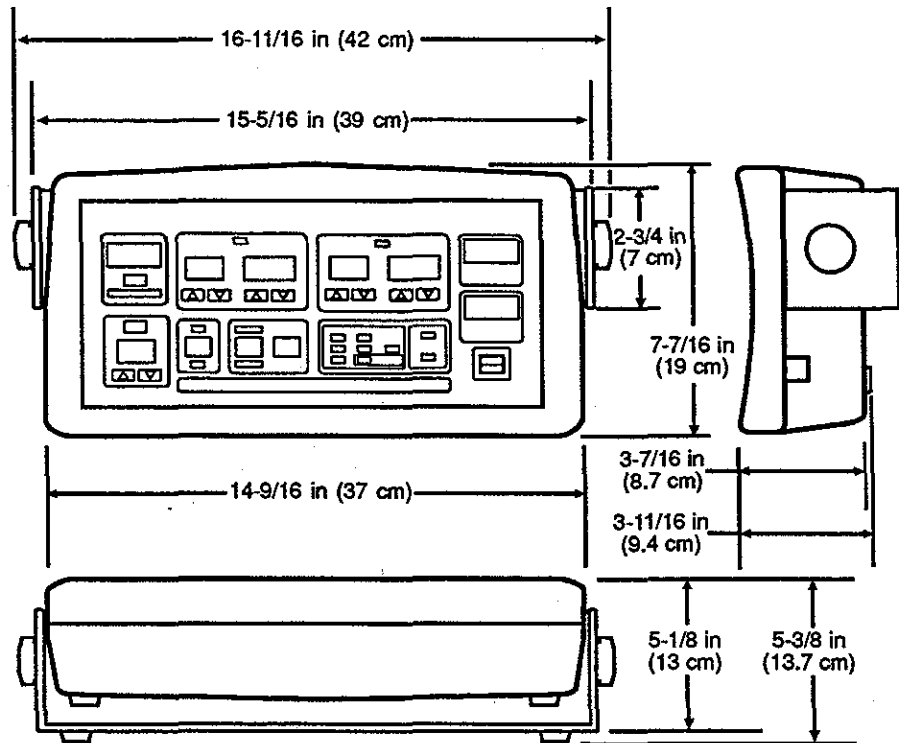


Figure 5.2:  
Remote Panel Dimensions

# System Specifications

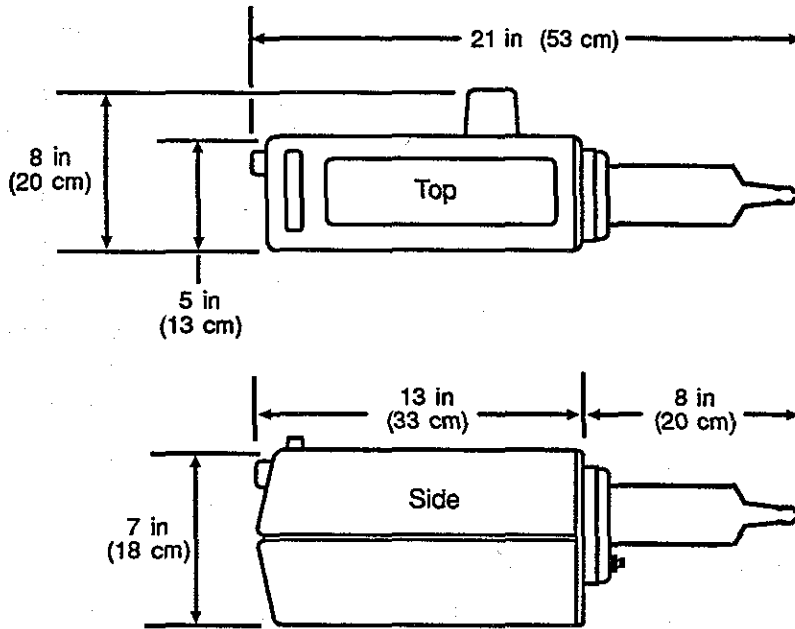


Figure 5.3:  
Front Load Injector Head Dimensions

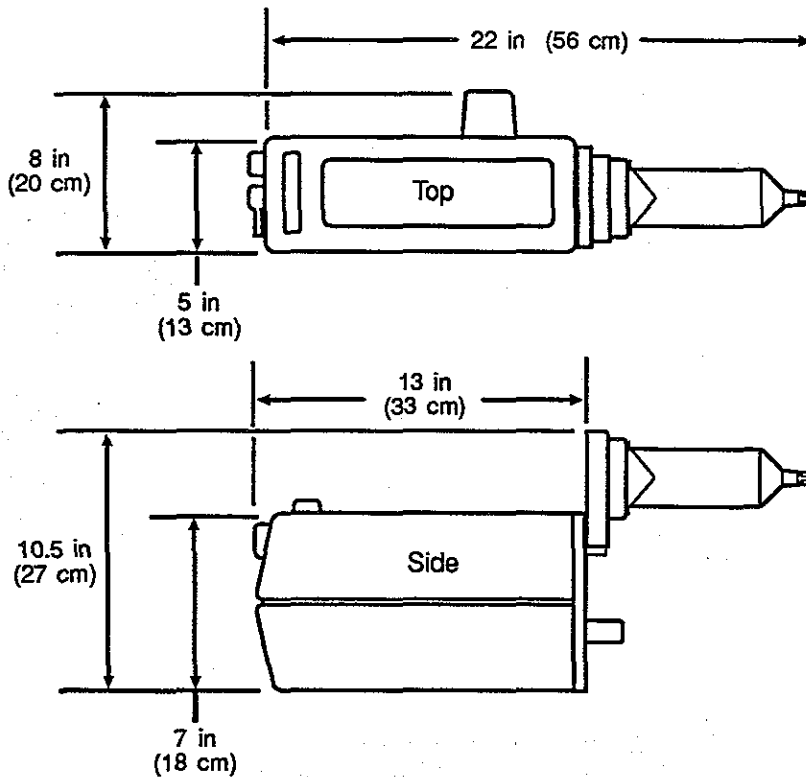


Figure 5.4:  
Turret Style Injector Head Dimensions

# System Specifications

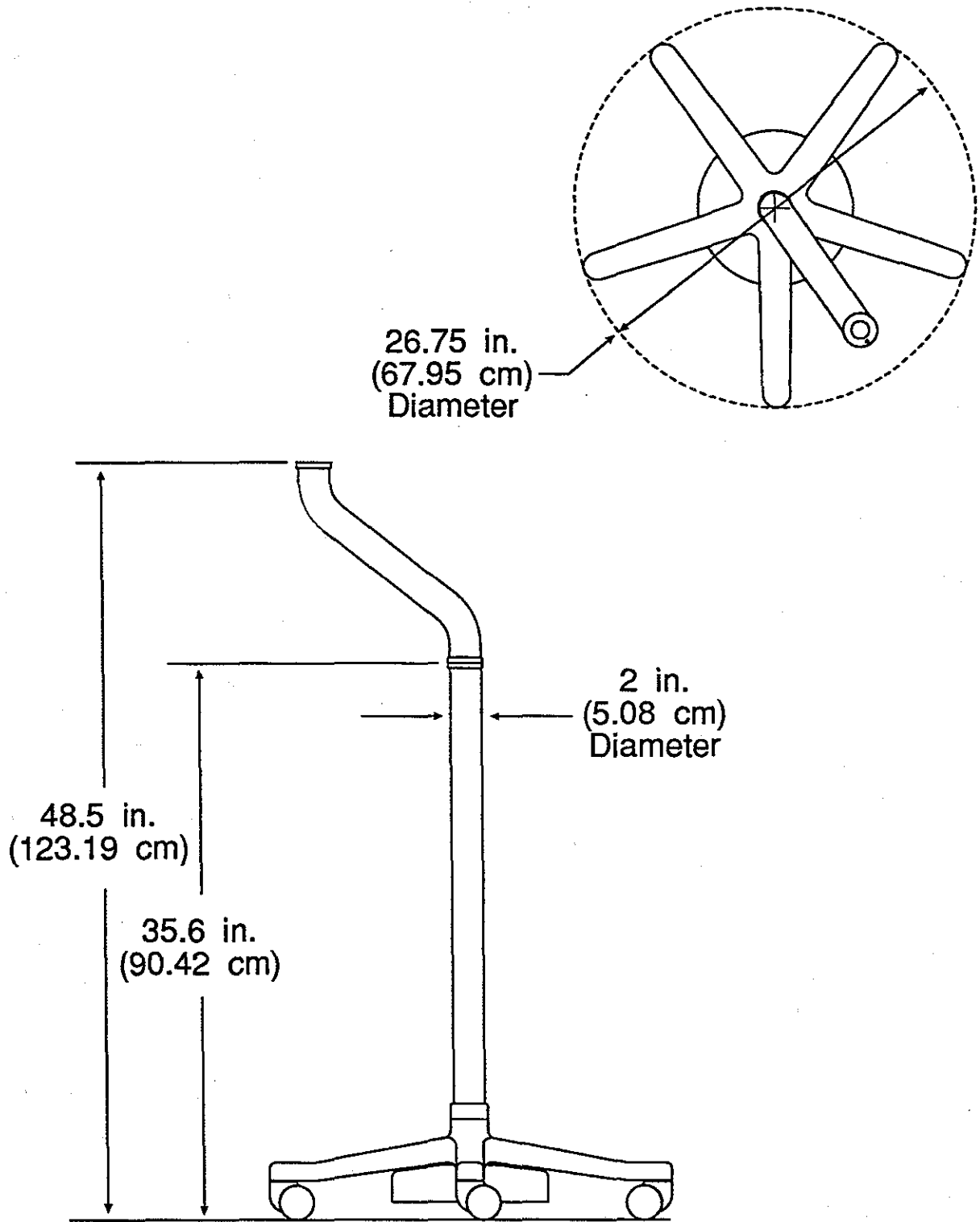


Figure 5.5:  
KMA-320T Dimensions

# System Specifications

Notes:

## SECTION 6 - OPERATIONAL CHECKOUT

This section is an operational checkout procedure to be used in the verification of proper operation of the following units: Control Room Console MCT 305/310, and MCT 305P/310P, Remote Panel Systems RPS 350 and 350P. The following basic procedures are included:

- Entering information on the Control Panel.
- Saving a new program.
- Editing a program previously stored in memory.
- Recalling a program stored in memory.
- Changing a displayed value without affecting the stored program.

### Entering an Injection Program

- 1.) Prior to programming the injector, apply power to the unit and reverse the head plunger completely by activating the **REVERSE** load on the head.
- 2.) Press the **PROGRAM SELECT** button until the desired **PROGRAM NUMBER** lights.
  - Each time **PROGRAM SELECT** is pressed, the numbers advance in order from 1 to 6, then repeat. Program number 6 is the exclusive storage location for multi-level injections. Single and Bi-phasic injections can be stored under this Program Number, but the Program Number 6 memory address is the only one with the capacity to store multi-level injections. As many as 5 levels can be programmed per multi-level injection. Each level contains 2 phases, providing a total of ten phases available to the user.
  - When programming the injector for the first time, we suggest starting with program number 1, then progressing in sequential order.
- 3.) If a multi-level injection is required, advance to storage location 6 using the **PROGRAM SELECT** key. Additional programming instructions for multi-level injections can be found on page 6 - 3.
- 4.) Press the **REPROGRAM** button. The **FIRST PHASE**, **SCAN DELAY**, and **SECOND PHASE** (if a two phase injection is selected) windows will flash.
- 5.) Select either a one or two phase injection by pressing the **PHASE** button, illuminating the proper indicator (**ONE** or **TWO**).
  - **ONE PHASE** Injection: One selected flow rate and one selected volume.
  - **TWO PHASE** Injection: Two selected flow rates and two selected volumes, one after the other, with no delay between phases.

## Operational Checkout

6.) Enter the **FIRST PHASE** settings.



- Press either button under the **FLOW RATE ML/SEC** window until the desired flow rate between 0.1 and 9.9 ml/sec. is displayed in the window.




- Press either button under **VOLUME ML** until the desired volume between 1 ml and 200 ml is displayed in the window.

7.) Use the guidelines provided above in step 6 to program **SECOND PHASE** flow rate and volume..

8.) Enter a scan delay as follows:



- Press either button under **SCAN DELAY SECONDS** until the desired scan delay of up to 99 seconds is displayed in the window.

Normally five beeps will sound when the scan delay reaches zero during the countdown, alerting the operator to trigger the scanner. To eliminate the sound of the beep, press the **SCANNER TONE ALERT** [  ] key once.

To restore the beep, press the **SCANNER TONE ALERT** button. The unit will beep to indicate that the tone has been restored.

9.) To store the injection values currently displayed on the Control Panel, press the **PROGRAM LOCK** key. These values are now stored under the current **PROGRAM NUMBER**.



# Operational Checkout

## Multi-Level Programming

Software version CU 2.0 and EP 2.0 and later added the capability of multiple level injections to MCT injectors. The following procedures do not apply to software releases prior to CU 2.0.

### Introduction

- 1.) Press the PROGRAM SELECT key until program number 6 is illuminated. This is the only program number with the capability to store Multi-Level programs.
- 2.) If a Multi-Level program has been previously stored, the INJECTION DURATION window will display PL 1.
- 3.) Press the REPROGRAM key. All programmable values on the Control Panel will flash. Multi-Level programming must start at PL 1.

### Entering a Multi-Level Program

- 4.) Use the arrow keys to change Flow Rate, Volume, and Scan Delay. Each level must contain two phases in order to advance to the next higher level.

**Note:** If Scan Delay is used, it must be programmed in the first level, and can not be set or changed in any program level above PL 1.

- 5.) Press the PROGRAM SELECT key to advance to the next level. Values in this level will also be flashing. Enter the desired values using the arrow keys. When both phases of this level are programmed, advance to the next level by pressing the PROGRAM SELECT key. A maximum of five two-phase levels can be programmed.
- 6.) When all values have been entered, press the PROGRAM LOCK key.

### Previewing a Multi-Level Program

- 7.) To preview a Multi-Level program, advance to program number 6. The INJECTION DURATION (program level) window displays PL 1. Pressing the PROGRAM SELECT key will advance the program to Level 2 (PL 2). Once the PROGRAM SELECT key has advanced to the highest level in the Multi-Level program, pressing the PROGRAM SELECT key again reverts to program number 1.

### Arming Sequence

- 8.) When arming the injector for Multi-Level injections, ensure that the highest PL level of the program is displayed in the INJECTION DURATION window. If this is not done, the message "SEE PL x" will be displayed in the Phase windows. Press the PROGRAM SELECT key to advance to the highest injection level, then arm.

### Editing a Multi-Level Program

- 9.) Multi-Level programs are edited in the same way as single level programs. Page 6 - 5 of this section explains how to edit stored injection programs.

## Operational Checkout

### Erasing a Multi-Level Program

- 10.) To erase one level of a Multi-Level program, advance to the level to be deleted, then press the REPROGRAM key. While the displayed values are flashing, press the SYSTEM ABORT/RESET key.

**CAUTION:** *The current level will be eliminated ALONG WITH ALL LEVELS ABOVE THE DELETED LEVEL. A single level can be deleted ONLY WHEN THE LEVEL TO BE DELETED IS THE HIGHEST LEVEL OF THE MULTI-LEVEL INJECTION.*

Press the PROGRAM LOCK key to complete this sequence.

### To UN-DO an Erase

- 11.) If it is decided to retain the original Multi-Level program before the PROGRAM LOCK button is pressed, press the SYSTEM ABORT/RESET bar to scroll downward until PL 1 is displayed in the Injection Duration window. Pressing the SYSTEM ABORT/RESET bar at this level will reset the Control Panel so that it stops flashing and the original Multi-Level program has been reinstated. The original Multi-Level program can be reviewed by pressing the PROGRAM SELECT key.

### Total Volume Indication

The VOLUME DELIVERED window will display Total Volume Selected if the PROGRAM LOCK key is pressed when the unit is not in the reprogramming mode.

While reprogramming a Multi-Level injection, the VOLUME DELIVERED window will display the total volume programmed in all levels.

### Storing an Injection Program

- 1.) After entering desired values in the reprogram mode, press the PROGRAM LOCK button, saving the program as displayed.
  - The program displayed on the Control Panel is now stored in the PPI memory.
  - Any values that were flashing will now stop.
  - The injection program is now stored in the program location illuminated in the PROGRAM NUMBER display.
  - To confirm that the injection has been stored, press the Program Select button until the desired program number is selected.
- 2.) To enter and save another program, press the PROGRAM SELECT button to advance to the next unused PROGRAM NUMBER (or any other storage location), then repeat the steps as described under *Entering and Saving programs*

## Operational Checkout

### Recalling a Stored Injection Program

To recall a stored injection program, press the PROGRAM SELECT button until the program location of the desired program is illuminated on the Control Panel. The Control Panel will then display the values stored for the selected injection sequence.

### Editing a Stored Injection Program

Once an injection sequence is stored in the MCT/MCT Plus, the injection parameters can be edited to either temporarily or permanently modify the stored program.

To permanently modify the stored program:

- 1.) Press the PROGRAM SELECT button until the desired PROGRAM NUMBER is illuminated. The Control Panel windows will display the injection values that have been stored.

Select the phase that is desired, ONE or TWO.

- 2.) Press the REPROGRAM button: All injection parameters will flash. Enter the new injection parameters and delay settings to be stored. The VOLUME DELIVERED window will display the TOTAL VOLUME SELECTED for the entire injection, be it single or multi-level.

Any value can now be selected for the FIRST PHASE, SECOND PHASE or SCAN DELAY.

- 3.) To store the injection program as displayed, press the PROGRAM LOCK button.

Any windows that were flashing will stop, and the program is now permanently stored in the PPI memory.

### Editing the Display Without Affecting the Stored Program

When a program is recalled, any value on the Control Panel can be changed to alter the injection for a particular procedure, without affecting the permanently stored program.

- 1.) Press the PROGRAM SELECT button until the desired PROGRAM NUMBER is illuminated.

- The Control Panel will display the injection values that are stored.

**DO NOT PRESS THE REPROGRAM BUTTON**

## Operational Checkout

2.) Enter the desired value changes.

- Enter any desired value changes for **FIRST PHASE**, **SECOND PHASE**, or **SCAN DELAY**. When a value change occurs, the **PROGRAM NUMBER** indicator light will flash from left to right. This indicates that the currently modified program no longer displays the stored injection values.
- In order to create a two-phase injection from an existing one-phase injection, press the **PHASE** button. The **TWO** phase indicator will illuminate, and the second phase windows will display default values. The desired second phase values can then be programmed.
- Likewise, in order to create a one-phase injection from an existing two-phase injection, press the **PHASE** button. The **ONE** phase indicator will illuminate and second phase values will no longer be displayed on the Control Panel.

3.) When the Control Panel is programmed as desired, the user can proceed with the other details of the injection. The stored program is not affected by the changes made.

4.) To restore the original recalled program, press the **PROGRAM SELECT** button once. The program number light will then stop flashing and remain illuminated.

### Hold Feature (Single Arm Only)

During a contrast enhanced CT procedure, patient discomfort occasionally dictates that the injection/scan be stopped. The operator can then resolve the patient's discomfort and resume the procedure. The "Hold" feature allows the operator to stop the injection for up to 5 minutes. Leaving the **MCT/MCT Plus** in "Hold" for greater than five minutes will disarm the injector, and trigger the "SEE OP 115" Message to appear on the Control Panel.

**Note:** When the injector is set for **MULTIPLE** arm, the Hold feature is disabled. If the **START** key is depressed during a **MULTIPLE** arm injection, the current injection will be aborted and the injector will re-arm if sufficient volume remains in the syringe.

If the start switch is depressed after a single arm injection has begun, the Hold feature is initialized and the following sequence will occur:

- 1.) The injection stops or "freezes."
- 2.) The **VOLUME DELIVERED** and **INJECTION DURATION** values "freeze."
- 3.) The **SCAN DELAY** countdown stops.
- 4.) The Phase light above the currently active **FLOW-VOLUME** window will blink.
- 5.) The first phase **FLOW-VOLUME** window will display and flash the word "Hold."
- 6.) All second phase information (if used) will blank out.

If the **START** switch is then pressed again, the **MCT/MCT Plus** will beep five times if the scan delay has previously reached zero, alerting the operator to trigger the scanner. If the scan delay has not reached zero, the procedure will resume from where the injection left off. Each time the **START** switch is depressed, the sequence in steps 1 through 6 will be repeated until the entire programmed CT injection is completed.

When in the "Hold" mode, if it is decided to terminate the procedure, pressing **RESET** or any other key will disarm the injector, triggering the "SEE OP 101" Message code to appear. Pressing **RESET** again will return the **MCT/MCT Plus** to the standby mode.

## SECTION 7 - CPU CARD SUMMARY OF OPERATION

### CPU Card Summary of Operation

There are two versions of the CPU Card, identified by the card revision level, part number, and serial number. Use the table below to determine which card type the injector has, then refer to the appropriate description on the following pages.

Assembly No.	Revision	Serial Number	CPU Card Type
78101-10-PC-01	Rev. E and above	20000 to 27467	Type B
78101-10-PC-14	Rev. A and above	Above 27467	Type C

# CPU Card

## CPU Card Type B (Part Number 78101-10-AP-01, serial numbers 20000 - 27467)

The CPU Card contains the main injector operating program. Refer to the figure on page 7 - 3 for a block diagram of the CPU Card type B. The CPU Card contains the following circuits:

### CPU

The Z-80 CPU U10 is supported by 32 Kbytes of EPROM, and 8 Kbytes of RAM. This processor provides a 16-bit address bus, an 8-bit data bus, and several control lines for Read/Write, Interrupts, etc.

### Clock

Crystal Y2 (8 MHz) and clock generator chip U2, produce the 4 MHz clock system for the CPU, and provide a Reset pulse when power is first applied. This pulse resets the CPU, and ensures that the program starts at the beginning.

### Memory

The CPU memory is comprised of four EPROM chips (U11 through U14), and one RAM chip (U16). These chips connect the address and data buses, and are controlled by the enable lines from the memory / I/O decoder, and the Read/Write lines from the CPU. If these chips need to be replaced, contact Medrad Service for current revision EPROMS.

### Memory / Input/Output Decoder

The Memory / Input/Output Decoder, comprised of U22 and U23, provides:

- Control of the memory chips during read/write cycles
- Two major decode lines (MEMIO1 and MEMIO2: Pins 16 and 17) that are buffered and sent to the Input/Output Card
- Direction control for the data buffer

### Bus Buffers

The control lines from the CPU are buffered by U17, then sent to the motherboard bus. The address bus is buffered by U18 and U19, then sent to the motherboard bus. The flow of data is controlled by bidirectional buffer U20, which is controlled by the Memory / I/O Decoder. Control and addressing lines flow away from the CPU, while the buffered data bus, can flow in either direction.

### Watchdog 1

The CPU watchdog (U3) must receive pulses from the I/O Card on the WATCHDOG 1 line, or CPU operations will be interrupted. The pulses are generated on the I/O Card, by decoding a CPU address. A loss of these pulses is most likely the result of a problem with the CPU or the EPROMs.

### Test Section

The circuit surrounding U7 is used for in-house testing, and the enabling of Utility programs.

# CPU Card

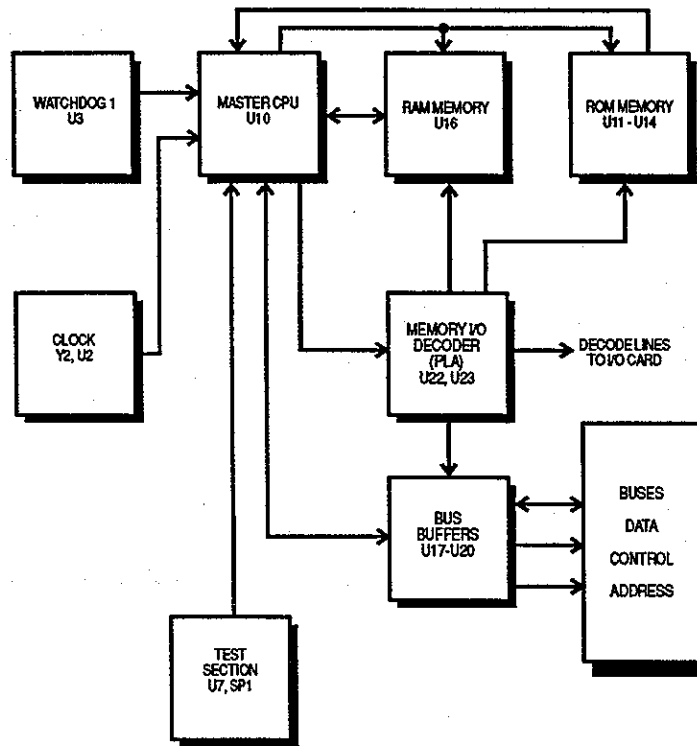


Figure 7.1: Block Diagram: CPU Card Type B

# CPU Card

## CPU Card Type C (Part Number 78101-10-AP-14, serial numbers above 27467)

The CPU Card contains the main injector operating program. Refer to the figure on page 7 - 5 for a block diagram of the CPU Card type C. The CPU Card contains the following circuits:

### CPU

The Z-80 CPU U10 is supported by 64 Kbytes of EPROM, and 8 Kbytes of RAM. This processor provides a 16-bit address bus, an 8-bit data bus, and several control lines for Read/Write, Interrupts, etc.

### Clock

Crystal Y2 (8 MHz) and clock generator chip U2, produce the 4 MHz clock system for the CPU, and provide a Reset pulse when power is first applied. This pulse resets the CPU, and ensures that the program starts at the beginning.

### Memory

The CPU memory is comprised of one EPROM chip (U24), and one RAM chip (U16). These chips connect the address and data buses, and are controlled by the enable lines from the memory / I/O decoder, and the Read/Write lines from the CPU. EPROM chip U24 is only accessed when switch SP2-2 is set to the "U24" position. If this chip must be replaced, contact Medrad Service for a current version EPROM. In addition to standard memory, a separate EPROM chip (U25) may be installed to access special utility functions. This may be accessed by setting switch SP2-2 to the "U25" position. Contact Medrad Factory Service for more information.

### Memory / Input/Output Decoder

The Memory / Input/Output Decoder, comprised of U23 and U26, provides:

- Control of the memory chips during read/write cycles
- Two major decode lines (MEMIO1 and MEMIO2: Pins 16 and 17) that are buffered and sent to the Input/Output Card
- Direction control for the data buffer

### Bus Buffers

The control lines from the CPU are buffered by U17, then sent to the motherboard bus. The address bus is buffered by U18 and U19, then sent to the motherboard bus. The flow of data is controlled by bidirectional buffer U20, which is controlled by the Memory / I/O Decoder. Control and addressing lines flow away from the CPU, while the buffered data bus, can flow in either direction.

### Watchdog 1

The CPU watchdog (U3) must receive pulses from the I/O Card on the WATCHDOG 1 line, or CPU operations will be interrupted. The pulses are generated on the I/O Card, by decoding a CPU address. A loss of these pulses is most likely the result of a problem with the CPU or the EPROMs.

### Test Section

The circuit surrounding U7 is used for in-house testing, and the enabling of Utility programs.



# CPU Card


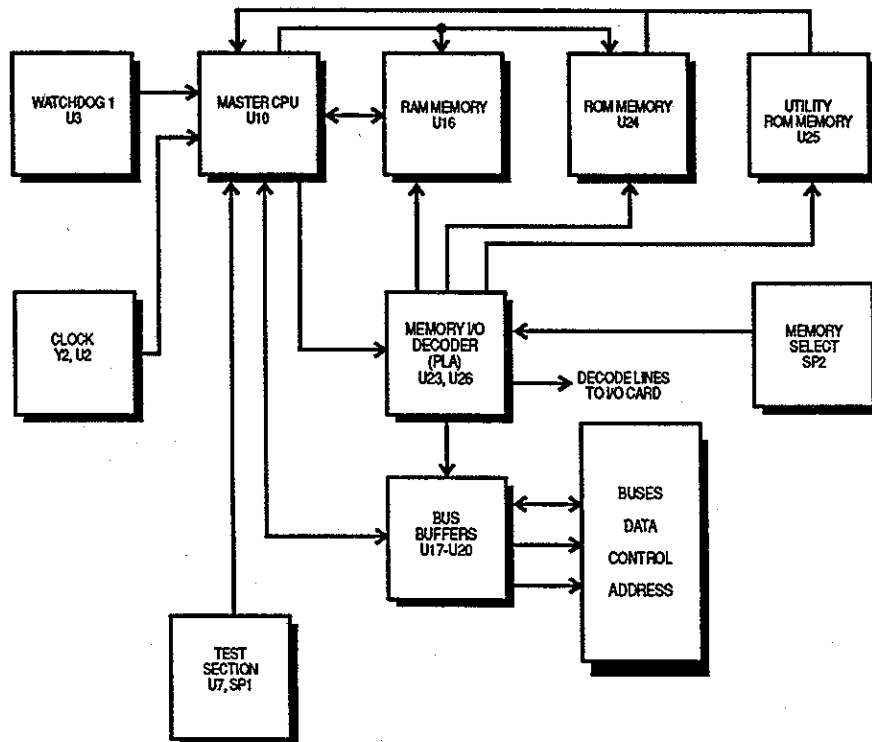


Figure 7.2: Block Diagram: CPU Card Type C

# CPU Card
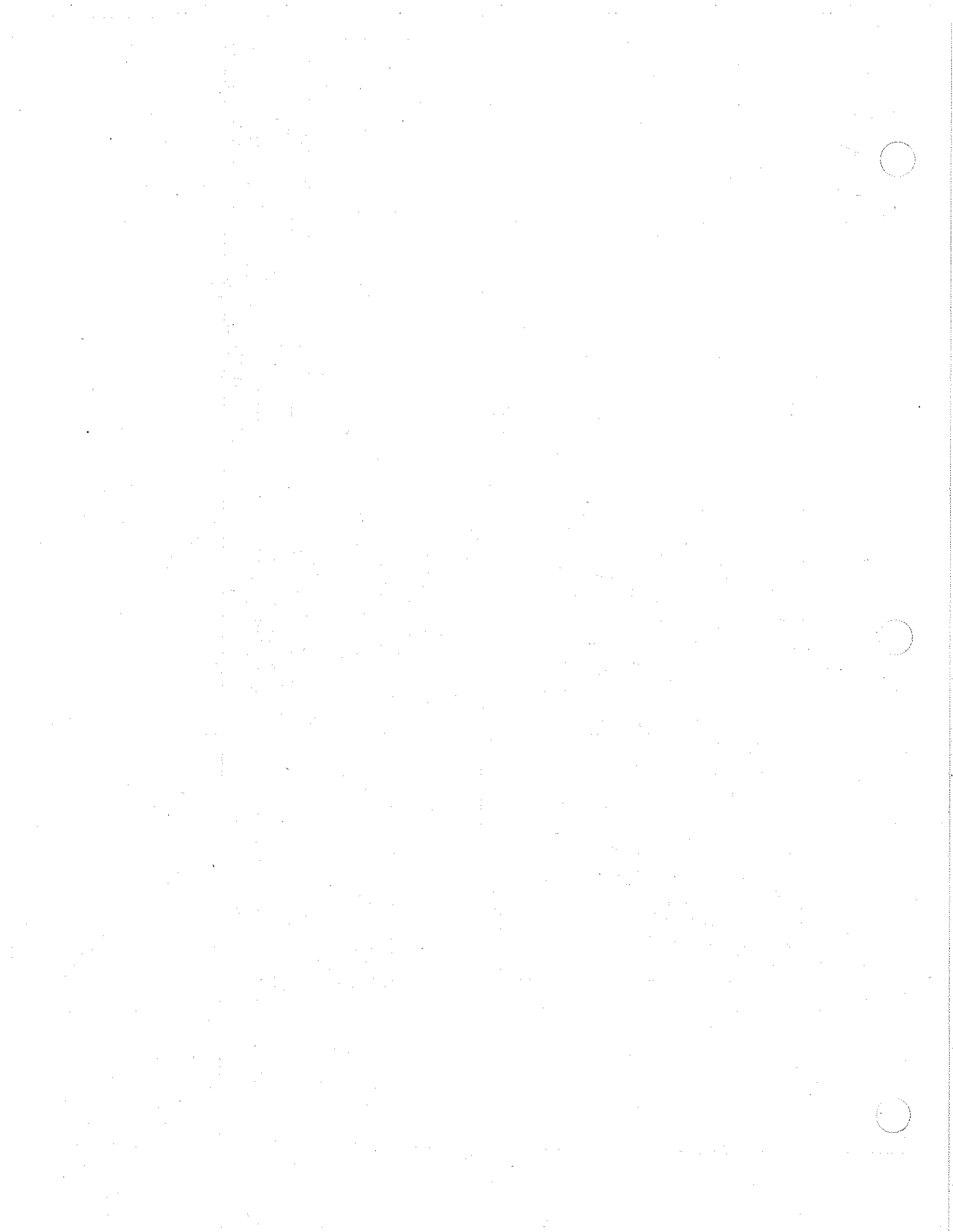
<i>Serial Number</i>	<i>Card Part Number</i>	<i>Chip Designation</i>	<i>Chip Part Number</i>
CPU 2XXXX to 27467	78101-10-AP-01 (-4)	U11	CU 2.1 U11
		U12	CU 2.1 U12
		U13	CU 2.1 U13
		U14	CU 2.1 U14
		U22	M5JA22
		U23	M5JA23
		CPU 27468 and Above	78101-10-AP-14 (-4)
U23	M5JA23		
U26	M5JA26		





# BILL OF MATERIAL

ITEM	DESCRIPTION	PART NUMBER	QTY.	REF. DESIG.	PROCESS CODE	PREP. CODE	ITEM	DESCRIPTION	PART NUMBER	QTY.	REF. DESIG.	PROCESS CODE
1	PC CARD REV. C. AND UP	78101-10-PC-14	1		LG	37						
2	CARD EJECTOR, WHITE	510-7000-000	2		AW	38						
3	CARD EJECTOR PIN	510-7001-000	2		AW	39						
4	TUBULAR TERMINAL	642-1238-000	3		BW	40	IC. CLOCK CONTROLLER, Z8581	070-8581-000	1		U2	DIP
5	LABEL, BLANK SERIAL NO	550-0093-652	1	CG. SYS CLK. NMT	BW	41	IC. MONOSTABLE, 74LS122	074-74LS-122	1		U3	DIP
6	IC SOCKET, 20 PIN	545-8201-020	2	U26, U23	LG	42	IC. HEX. INVERTER, 74PC04	080-74PC-004	1		U5	DIP
7	IC SOCKET, 28 PIN	545-8281-028	2	U24, U25	LG	43	IC. BUFFER, 74LS541	074-74LS-541	4		U7, U17-19	DIP
8	IC SOCKET, 40 PIN	545-8401-040	1	U10	LG	44	IC. SRAM, 8K X. 8	071-8464-000	1		U16	LG
9	IC SOCKET, 40 PIN	545-8401-040	1		LG	45	IC. BUS TRANS., 74LS645	074-74LS-645	1		U20	DIP
10						46	IC. 7800 CPU	070-8400-000	1		U10	AW
11												
12												
13												
14												
15												
16												
17												
18												
19												
20	CAP, 10 UF, 35V	204-0106-135	2	C32, 33	LG	1-LS						
21	CAP, 2.2 UF, 35V	206-0225-350	1	C2	LG	1-LS						
22	CAP, .01 UF, 50V	201-0103-050	14	C3-5, 10, 11, 13, 18, 24-30	LG	1-LS						
23												
24	CAP, 33PF, 100V	201-0330-100	2	C14, 15	LG	1-LS						
25	SWITCH, DIP, 16 PIN SPST	450-4358-021	1	SP1	BW							
26	SWITCH, KEY, SPST MOMENTARY	450-1152-000	1	S1	LG							
27	SWITCH, DIP, 2 POS	450-7802-002	1	SP2	BW							
28	CRYSTAL, 8.0 MHZ	180-0080-020	1	Y2	LG							
29												
30	RESISTOR, 200K OHM 1/4W 5%	303-0204-000	1	R1	AXI	5BEND						
31	RESISTOR, 100 OHM 1/4W 5%	303-0101-000	3	R16, 17, 18	AXI	5BEND						
32	RESISTOR, 4.7K OHM 1/4W 5%	303-0472-000	1	R7	AXI	5BEND						
33	RESISTOR, 1K OHM 1/4W 5%	303-0102-000	2	R14, R15	AXI	5BEND						
34												
35	RESISTOR ARRAY, 10K OHM SIP	340-4310-103	6	RN1-6	LG	5BEND						
36												

**MEDRAD INC.**  
 271 KAPPA DRIVE  
 PITTSBURGH, PA 15238

**MKV/MCT INJECTOR  
 CPU BOARD ASSEMBLY**

DWG. NUMBER: **78101-10-AP-14-1**  
 REV. **F**

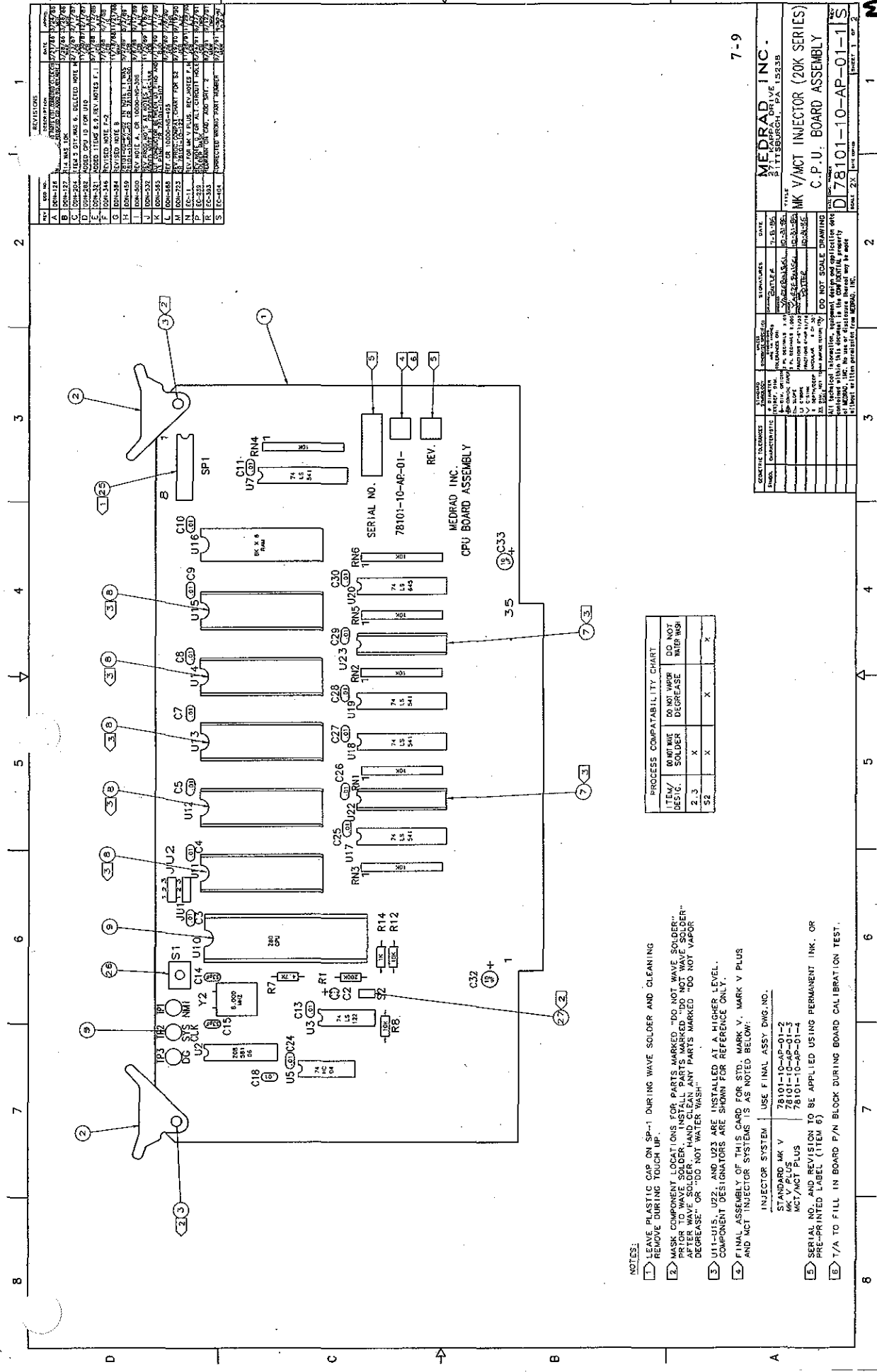
DATE COPIED: \_\_\_\_\_ SHEET 2 OF 2

All technical information equipment design and application data contained within this document is the CONFIDENTIAL Property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.

SIGNATURES	DATE	SEE
DRAWN: E. J. PAVLOCK	1/30/92	SHT. 1
CHECKED: B. J. O'HELLA	1/30/92	FOR
ENGR: R. T. BOLLINGER	2/10/92	REV.
PROJ. MGR: G. F. HINTERMEYER	2/10/92	LEVEL



REV.	DATE	DESCRIPTION
A	11/27/78	INITIAL DESIGN
B	12/11/78	REVISED TO ADD COMPONENTS
C	01/11/79	REVISED TO ADD COMPONENTS
D	02/11/79	REVISED TO ADD COMPONENTS
E	03/11/79	REVISED TO ADD COMPONENTS
F	04/11/79	REVISED TO ADD COMPONENTS
G	05/11/79	REVISED TO ADD COMPONENTS
H	06/11/79	REVISED TO ADD COMPONENTS
I	07/11/79	REVISED TO ADD COMPONENTS
J	08/11/79	REVISED TO ADD COMPONENTS
K	09/11/79	REVISED TO ADD COMPONENTS
L	10/11/79	REVISED TO ADD COMPONENTS
M	11/11/79	REVISED TO ADD COMPONENTS
N	12/11/79	REVISED TO ADD COMPONENTS
O	01/11/80	REVISED TO ADD COMPONENTS
P	02/11/80	REVISED TO ADD COMPONENTS
Q	03/11/80	REVISED TO ADD COMPONENTS
R	04/11/80	REVISED TO ADD COMPONENTS
S	05/11/80	REVISED TO ADD COMPONENTS




PROCESS COMPATIBILITY CHART

ITEM/DESIGN	DO NOT SOLDER	DO NOT DECREASE	DO NOT INCREASE	DO NOT INCREASE (RIB)
2, 3	X		X	X
5, 6		X		

- NOTES:
- LEAVE PLASTIC CAP ON SP-1 DURING WAVE SOLDER AND CLEANING REMOVE DURING TOUCH UP.
  - MASK COMPONENT LOCATIONS FOR PARTS MARKED "DO NOT WAVE SOLDER" PRIOR TO WAVE SOLDER. INSTALL PARTS MARKED "DO NOT WAVE SOLDER" AFTER WAVE SOLDER. WAVE SOLDER PARTS MARKED "DO NOT WAVE SOLDER" OR "DO NOT WATER WASH" OR "DO NOT WATER WASH".
  - U1, U15, U18, AND U23 ARE INSTALLED AT A MICRER LEVEL. COMPONENT DESIGNATORS ARE SHOWN FOR REFERENCE ONLY.
  - FINAL ASSEMBLY OF THIS CARD FOR STD. MARK V, MARK V PLUS AND MCT INJECTOR SYSTEMS IS AS NOTED BELOW:
- | INJECTOR SYSTEM | USE FINAL ASSY DWG. NO. |
|-----------------|-------------------------|
| STANDARD MK V   | 78101-10-AP-01-2        |
| MCT/ACT PLUS    | 78101-10-AP-01-5        |
| MCT/ACT PLUS    | 78101-10-AP-01-4        |
- SERIAL NO. AND REVISION TO BE APPLIED USING PERMANENT INK. OR PRE-PRINTED LABEL (ITEM 6)
  - T/A TO FILL IN BOARD P/N BLOCK DURING BOARD CALIBRATION TEST.

SYMBOLS	DATE	DESCRIPTION
1-6	11/27/78	INITIAL DESIGN
7-9	12/11/78	REVISED TO ADD COMPONENTS
10-12	01/11/79	REVISED TO ADD COMPONENTS
13-15	02/11/79	REVISED TO ADD COMPONENTS
16-18	03/11/79	REVISED TO ADD COMPONENTS
19-21	04/11/79	REVISED TO ADD COMPONENTS
22-24	05/11/79	REVISED TO ADD COMPONENTS
25-27	06/11/79	REVISED TO ADD COMPONENTS
28-30	07/11/79	REVISED TO ADD COMPONENTS
31-33	08/11/79	REVISED TO ADD COMPONENTS
34-36	09/11/79	REVISED TO ADD COMPONENTS
37-39	10/11/79	REVISED TO ADD COMPONENTS
40-42	11/11/79	REVISED TO ADD COMPONENTS
43-45	12/11/79	REVISED TO ADD COMPONENTS
46-48	01/11/80	REVISED TO ADD COMPONENTS
49-51	02/11/80	REVISED TO ADD COMPONENTS
52-54	03/11/80	REVISED TO ADD COMPONENTS
55-57	04/11/80	REVISED TO ADD COMPONENTS
58-60	05/11/80	REVISED TO ADD COMPONENTS
61-63	06/11/80	REVISED TO ADD COMPONENTS
64-66	07/11/80	REVISED TO ADD COMPONENTS
67-69	08/11/80	REVISED TO ADD COMPONENTS
70-72	09/11/80	REVISED TO ADD COMPONENTS
73-75	10/11/80	REVISED TO ADD COMPONENTS
76-78	11/11/80	REVISED TO ADD COMPONENTS
79-81	12/11/80	REVISED TO ADD COMPONENTS
82-84	01/11/81	REVISED TO ADD COMPONENTS
85-87	02/11/81	REVISED TO ADD COMPONENTS
88-90	03/11/81	REVISED TO ADD COMPONENTS
91-93	04/11/81	REVISED TO ADD COMPONENTS
94-96	05/11/81	REVISED TO ADD COMPONENTS
97-99	06/11/81	REVISED TO ADD COMPONENTS
100-102	07/11/81	REVISED TO ADD COMPONENTS
103-105	08/11/81	REVISED TO ADD COMPONENTS
106-108	09/11/81	REVISED TO ADD COMPONENTS
109-111	10/11/81	REVISED TO ADD COMPONENTS
112-114	11/11/81	REVISED TO ADD COMPONENTS
115-117	12/11/81	REVISED TO ADD COMPONENTS
118-120	01/11/82	REVISED TO ADD COMPONENTS
121-123	02/11/82	REVISED TO ADD COMPONENTS
124-126	03/11/82	REVISED TO ADD COMPONENTS
127-129	04/11/82	REVISED TO ADD COMPONENTS
130-132	05/11/82	REVISED TO ADD COMPONENTS
133-135	06/11/82	REVISED TO ADD COMPONENTS
136-138	07/11/82	REVISED TO ADD COMPONENTS
139-141	08/11/82	REVISED TO ADD COMPONENTS
142-144	09/11/82	REVISED TO ADD COMPONENTS
145-147	10/11/82	REVISED TO ADD COMPONENTS
148-150	11/11/82	REVISED TO ADD COMPONENTS
151-153	12/11/82	REVISED TO ADD COMPONENTS
154-156	01/11/83	REVISED TO ADD COMPONENTS
157-159	02/11/83	REVISED TO ADD COMPONENTS
160-162	03/11/83	REVISED TO ADD COMPONENTS
163-165	04/11/83	REVISED TO ADD COMPONENTS
166-168	05/11/83	REVISED TO ADD COMPONENTS
169-171	06/11/83	REVISED TO ADD COMPONENTS
172-174	07/11/83	REVISED TO ADD COMPONENTS
175-177	08/11/83	REVISED TO ADD COMPONENTS
178-180	09/11/83	REVISED TO ADD COMPONENTS
181-183	10/11/83	REVISED TO ADD COMPONENTS
184-186	11/11/83	REVISED TO ADD COMPONENTS
187-189	12/11/83	REVISED TO ADD COMPONENTS
190-192	01/11/84	REVISED TO ADD COMPONENTS
193-195	02/11/84	REVISED TO ADD COMPONENTS
196-198	03/11/84	REVISED TO ADD COMPONENTS
199-201	04/11/84	REVISED TO ADD COMPONENTS
202-204	05/11/84	REVISED TO ADD COMPONENTS
205-207	06/11/84	REVISED TO ADD COMPONENTS
208-210	07/11/84	REVISED TO ADD COMPONENTS
211-213	08/11/84	REVISED TO ADD COMPONENTS
214-216	09/11/84	REVISED TO ADD COMPONENTS
217-219	10/11/84	REVISED TO ADD COMPONENTS
220-222	11/11/84	REVISED TO ADD COMPONENTS
223-225	12/11/84	REVISED TO ADD COMPONENTS
226-228	01/11/85	REVISED TO ADD COMPONENTS
229-231	02/11/85	REVISED TO ADD COMPONENTS
232-234	03/11/85	REVISED TO ADD COMPONENTS
235-237	04/11/85	REVISED TO ADD COMPONENTS
238-240	05/11/85	REVISED TO ADD COMPONENTS
241-243	06/11/85	REVISED TO ADD COMPONENTS
244-246	07/11/85	REVISED TO ADD COMPONENTS
247-249	08/11/85	REVISED TO ADD COMPONENTS
250-252	09/11/85	REVISED TO ADD COMPONENTS
253-255	10/11/85	REVISED TO ADD COMPONENTS
256-258	11/11/85	REVISED TO ADD COMPONENTS
259-261	12/11/85	REVISED TO ADD COMPONENTS
262-264	01/11/86	REVISED TO ADD COMPONENTS
265-267	02/11/86	REVISED TO ADD COMPONENTS
268-270	03/11/86	REVISED TO ADD COMPONENTS
271-273	04/11/86	REVISED TO ADD COMPONENTS
274-276	05/11/86	REVISED TO ADD COMPONENTS
277-279	06/11/86	REVISED TO ADD COMPONENTS
280-282	07/11/86	REVISED TO ADD COMPONENTS
283-285	08/11/86	REVISED TO ADD COMPONENTS
286-288	09/11/86	REVISED TO ADD COMPONENTS
289-291	10/11/86	REVISED TO ADD COMPONENTS
292-294	11/11/86	REVISED TO ADD COMPONENTS
295-297	12/11/86	REVISED TO ADD COMPONENTS
298-300	01/11/87	REVISED TO ADD COMPONENTS
301-303	02/11/87	REVISED TO ADD COMPONENTS
304-306	03/11/87	REVISED TO ADD COMPONENTS
307-309	04/11/87	REVISED TO ADD COMPONENTS
310-312	05/11/87	REVISED TO ADD COMPONENTS
313-315	06/11/87	REVISED TO ADD COMPONENTS
316-318	07/11/87	REVISED TO ADD COMPONENTS
319-321	08/11/87	REVISED TO ADD COMPONENTS
322-324	09/11/87	REVISED TO ADD COMPONENTS
325-327	10/11/87	REVISED TO ADD COMPONENTS
328-330	11/11/87	REVISED TO ADD COMPONENTS
331-333	12/11/87	REVISED TO ADD COMPONENTS
334-336	01/11/88	REVISED TO ADD COMPONENTS
337-339	02/11/88	REVISED TO ADD COMPONENTS
340-342	03/11/88	REVISED TO ADD COMPONENTS
343-345	04/11/88	REVISED TO ADD COMPONENTS
346-348	05/11/88	REVISED TO ADD COMPONENTS
349-351	06/11/88	REVISED TO ADD COMPONENTS
352-354	07/11/88	REVISED TO ADD COMPONENTS
355-357	08/11/88	REVISED TO ADD COMPONENTS
358-360	09/11/88	REVISED TO ADD COMPONENTS
361-363	10/11/88	REVISED TO ADD COMPONENTS
364-366	11/11/88	REVISED TO ADD COMPONENTS
367-369	12/11/88	REVISED TO ADD COMPONENTS
370-372	01/11/89	REVISED TO ADD COMPONENTS
373-375	02/11/89	REVISED TO ADD COMPONENTS
376-378	03/11/89	REVISED TO ADD COMPONENTS
379-381	04/11/89	REVISED TO ADD COMPONENTS
382-384	05/11/89	REVISED TO ADD COMPONENTS
385-387	06/11/89	REVISED TO ADD COMPONENTS
388-390	07/11/89	REVISED TO ADD COMPONENTS
391-393	08/11/89	REVISED TO ADD COMPONENTS
394-396	09/11/89	REVISED TO ADD COMPONENTS
397-399	10/11/89	REVISED TO ADD COMPONENTS
400-402	11/11/89	REVISED TO ADD COMPONENTS
403-405	12/11/89	REVISED TO ADD COMPONENTS
406-408	01/11/90	REVISED TO ADD COMPONENTS
409-411	02/11/90	REVISED TO ADD COMPONENTS
412-414	03/11/90	REVISED TO ADD COMPONENTS
415-417	04/11/90	REVISED TO ADD COMPONENTS
418-420	05/11/90	REVISED TO ADD COMPONENTS
421-423	06/11/90	REVISED TO ADD COMPONENTS
424-426	07/11/90	REVISED TO ADD COMPONENTS
427-429	08/11/90	REVISED TO ADD COMPONENTS
430-432	09/11/90	REVISED TO ADD COMPONENTS
433-435	10/11/90	REVISED TO ADD COMPONENTS
436-438	11/11/90	REVISED TO ADD COMPONENTS
439-441	12/11/90	REVISED TO ADD COMPONENTS
442-444	01/11/91	REVISED TO ADD COMPONENTS
445-447	02/11/91	REVISED TO ADD COMPONENTS
448-450	03/11/91	REVISED TO ADD COMPONENTS
451-453	04/11/91	REVISED TO ADD COMPONENTS
454-456	05/11/91	REVISED TO ADD COMPONENTS
457-459	06/11/91	REVISED TO ADD COMPONENTS
460-462	07/11/91	REVISED TO ADD COMPONENTS
463-465	08/11/91	REVISED TO ADD COMPONENTS
466-468	09/11/91	REVISED TO ADD COMPONENTS
469-471	10/11/91	REVISED TO ADD COMPONENTS
472-474	11/11/91	REVISED TO ADD COMPONENTS
475-477	12/11/91	REVISED TO ADD COMPONENTS
478-480	01/11/92	REVISED TO ADD COMPONENTS
481-483	02/11/92	REVISED TO ADD COMPONENTS
484-486	03/11/92	REVISED TO ADD COMPONENTS
487-489	04/11/92	REVISED TO ADD COMPONENTS
490-492	05/11/92	REVISED TO ADD COMPONENTS
493-495	06/11/92	REVISED TO ADD COMPONENTS
496-498	07/11/92	REVISED TO ADD COMPONENTS
499-501	08/11/92	REVISED TO ADD COMPONENTS
502-504	09/11/92	REVISED TO ADD COMPONENTS
505-507	10/11/92	REVISED TO ADD COMPONENTS
508-510	11/11/92	REVISED TO ADD COMPONENTS
511-513	12/11/92	REVISED TO ADD COMPONENTS
514-516	01/11/93	REVISED TO ADD COMPONENTS
517-519	02/11/93	REVISED TO ADD COMPONENTS
520-522	03/11/93	REVISED TO ADD COMPONENTS
523-525	04/11/93	REVISED TO ADD COMPONENTS
526-528	05/11/93	REVISED TO ADD COMPONENTS
529-531	06/11/93	REVISED TO ADD COMPONENTS
532-534	07/11/93	REVISED TO ADD COMPONENTS
535-537	08/11/93	REVISED TO ADD COMPONENTS
538-540	09/11/93	REVISED TO ADD COMPONENTS
541-543	10/11/93	REVISED TO ADD COMPONENTS
544-546	11/11/93	REVISED TO ADD COMPONENTS
547-549	12/11/93	REVISED TO ADD COMPONENTS
550-552	01/11/94	REVISED TO ADD COMPONENTS
553-555	02/11/94	REVISED TO ADD COMPONENTS
556-558	03/11/94	REVISED TO ADD COMPONENTS
559-561	04/11/94	REVISED TO ADD COMPONENTS
562-564	05/11/94	REVISED TO ADD COMPONENTS
565-567	06/11/94	REVISED TO ADD COMPONENTS
568-570	07/11/94	REVISED TO ADD COMPONENTS
571-573	08/11/94	REVISED TO ADD COMPONENTS
574-576	09/11/94	REVISED TO ADD COMPONENTS
577-579	10/11/94	REVISED TO ADD COMPONENTS
580-582	11/11/94	REVISED TO ADD COMPONENTS
583-585	12/11/94	REVISED TO ADD COMPONENTS
586-588	01/11/95	REVISED TO ADD COMPONENTS
589-591	02/11/95	REVISED TO ADD COMPONENTS
592-594	03/11/95	REVISED TO ADD COMPONENTS
595-597	04/11/95	REVISED TO ADD COMPONENTS
598-600	05/11/95	REVISED TO ADD COMPONENTS
601-603	06/11/95	REVISED TO ADD COMPONENTS
604-606	07/11/95	REVISED TO ADD COMPONENTS
607-609	08/11/95	REVISED TO ADD COMPONENTS
610-612	09/11/95	REVISED TO ADD COMPONENTS
613-615	10/11/95	REVISED TO ADD COMPONENTS
616-618	11/11/95	REVISED TO ADD COMPONENTS
619-621	12/11/95	REVISED TO ADD COMPONENTS
622-624	01/11/96	REVISED TO ADD COMPONENTS
625-627	02/11/96	REVISED TO ADD COMPONENTS
628-630	03/11/96	REVISED TO ADD COMPONENTS
631-633	04/11/96	REVISED TO ADD COMPONENTS
634-636	05/11/96	REVISED TO ADD COMPONENTS
637-639	06/11/96	REVISED TO ADD COMPONENTS
640-642	07/11/96	REVISED TO ADD COMPONENTS
643-645	08/11/96	REVISED TO ADD COMPONENTS
646-648	09/11/96	REVISED TO ADD COMPONENTS
649-651	10/11/96	REVISED TO ADD COMPONENTS
652-654	11/11/96	REVISED TO ADD COMPONENTS
655-657	12/11/96	REVISED TO ADD COMPONENTS
658-660	01/11/97	REVISED TO ADD COMPONENTS
661-663	02/11/97	REVISED TO ADD COMPONENTS
664-666	03/11/97	REVISED TO ADD COMPONENTS
667-669	04/11/97	REVISED TO ADD COMPONENTS
670-672	05/11/97	REVISED TO ADD COMPONENTS
673-675	06/11/97	REVISED TO ADD COMPONENTS
676-678	07/11/97	REVISED TO ADD COMPONENTS
679-681	08/11/97	REVISED TO ADD COMPONENTS
682-684	09/11/97	REVISED TO ADD COMPONENTS
685-687	10/11/97	REVISED TO ADD COMPONENTS
688-690	11/11/97	REVISED TO ADD COMPONENTS
691-693	12/11/97	REVISED TO ADD COMPONENTS
694-696	01/11/98	REVISED TO ADD COMPONENTS
697-699	02/11/98	REVISED TO ADD COMPONENTS
700-702	03/11/98	REVISED TO ADD COMPONENTS
703-705	04/11/98	REVISED TO ADD COMPONENTS
706-708	05/11/98	REVISED TO ADD COMPONENTS
709-711	06/11/98	

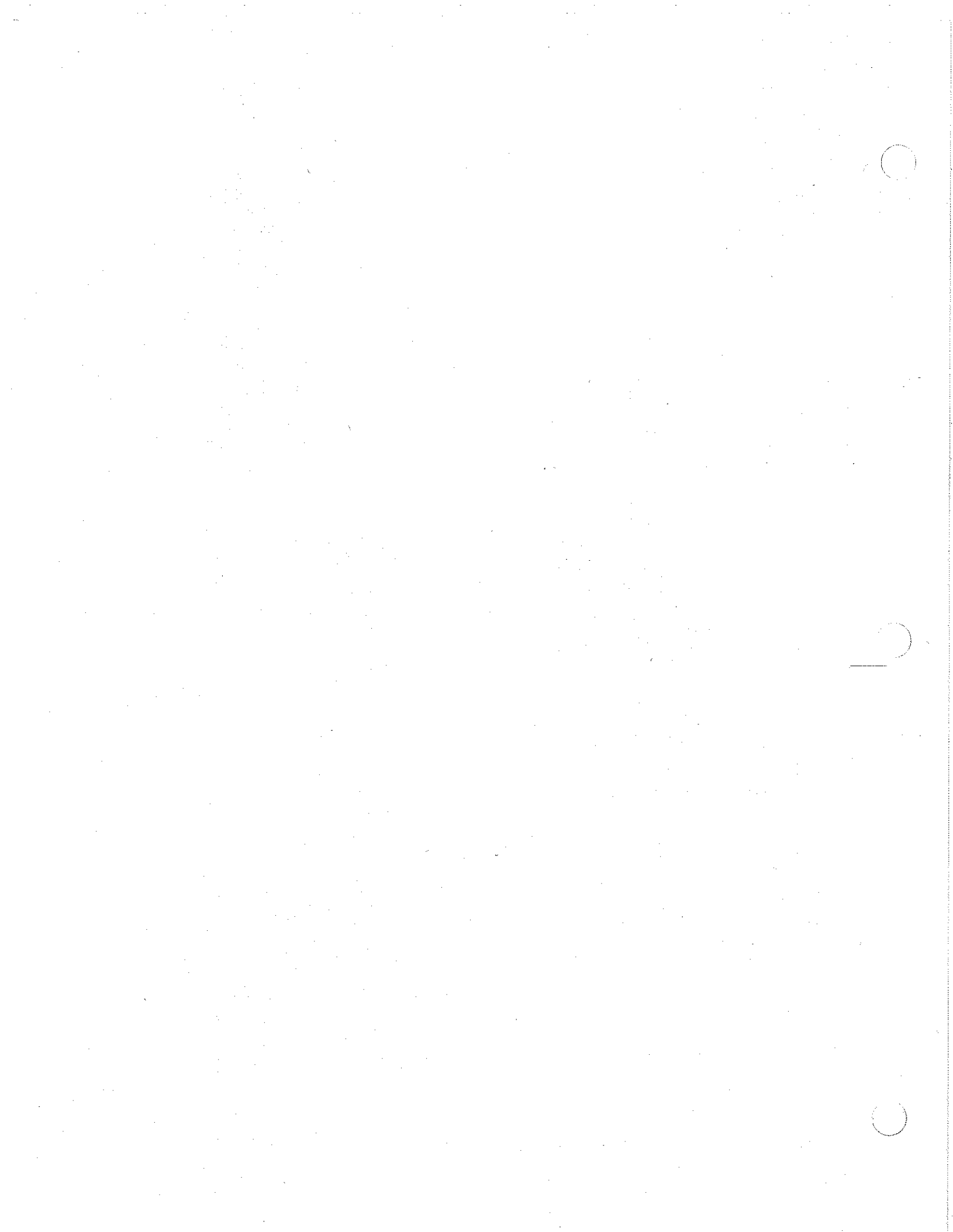



# BILL OF MATERIAL

ITEM	DESCRIPTION	PART NUMBER	QTY.	REF. DESIG.	PROCESS PREP. CODE	ITEM	DESCRIPTION	PART NUMBER	QTY.	REF. DESIG.	PROCESS PREP. CODE
1	PC CARD, REV. H AND UP	78101-10-PC-01	1		LG	37					
2	CARD EJECTOR, WHITE	510-7000-000	2		AW	38					
3	CARD EJECTOR PIN	510-7001-000	2		AW	39					
4	TUBULAR TERMINAL	642-1238-000	3		BW	40	IC CLOCK CONTROLLER, Z8581	070-8581-000	1	U2	G-008
5	ABEL, BLANK SERIAL NO	550-0093-652	1		BW	41	IC MONOSTABLE, 74LS122	074-74LS-122	1	U3	G-013
6	IC SOCKET, 20 PIN	545-8201-020	2	U22, U23	G-016	42	IC HEX INVERTER, 74PC04	080-74PC-004	1	U5	G-010
7	IC SOCKET, 28 PIN	545-8281-028	5	U11, U12, U13, U14, U15	G-002	43	IC BUFFER, 74LS541	074-74LS-541	4	U7, U17-19	G-004
8	IC SOCKET, 40 PIN	545-8401-040	1	U10	G-014	44	IC SRAM, 8K X 8	071-8464-000	1	U16	G-018
9	IC SOCKET, 40 PIN	545-8401-040	1		G-014	45	IC BUS TRANS, 74LS645	074-74LS-645	1	U20	G-006
10						46	IC Z80 CPU	070-8400-000	1	U10	AW
11											
12											
13											
14											
15											
16											
17											
18											
19											
20	CAP, 10 UF, 35V	204-0106-135	2	C32, C33	G-009, I-LS						
21	CAP, 2.2 UF, 35V	206-0225-350	1	C2	G-015, I-LS						
22	CAP, .01 UF, 50V	201-0103-050	17	C3-5, 7-11, 13, 18, 24-30	G-019, I-LS						
23											
24	CAP, 33PF, 100V	201-0330-100	2	C14, C15	G-012, I-LS						
25	SWITCH, DIP, 16 PIN SPST	450-4358-021	1	SP1	G-021						
26	SWITCH, KEY, SPST MOMENTARY	450-1152-000	1	S1	G-020						
27	SWITCH, SPST	450-4356-651	1	S2	BW						
28	CRYSTAL, 8.0 MHZ	180-0080-020	1	Y2	G-011						
29											
30	RESISTOR, 200K OHM 1/4W 5%	303-0204-000	1	R1	G-003, 5BEND						
31	RESISTOR, 10K OHM 1/4W 5%	303-0103-000	2	R8, R12	G-001, 5BEND						
32	RESISTOR, 4.7K OHM 1/4W 5%	303-0472-000	1	R7	G-007, 5BEND						
33	RESISTOR, 1K OHM 1/4W 5%	303-0102-000	1	R14	G-005, 5BEND						
34											
35	RESISTOR ARRAY, 10K OHM SIP	340-4310-103	6	RN1-6	G-017, 5BEND						
36											

**MEDRAD INC.**  
 271 KAPPA DRIVE  
 PITTSBURGH, PA 15238

TITLE  
**MKV/MCT INJECTOR (20K SERIES)  
 CPU BOARD ASSEMBLY**

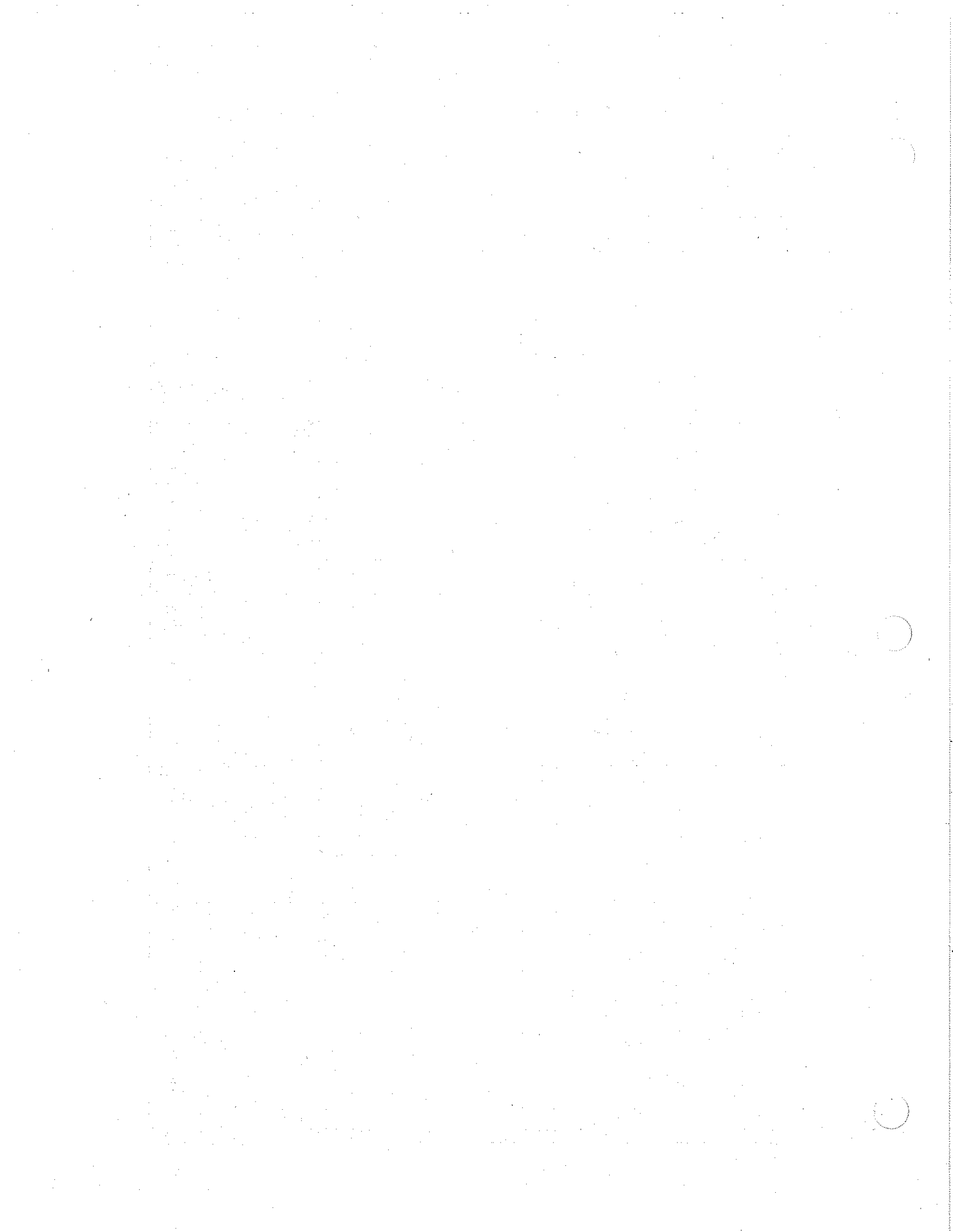
DWC. NUMBER  
**78101-10-AP-01-1**

DATE COPIED  
 SHEET 2 OF 2

All technical information equipment design and application data contained within this document is the CONFIDENTIAL Property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.

SEE SHT. 1 FOR REV. LEVEL

SIGNATURES	DATE
M. W. HALEN	8/2/91
J. P. LANE	7/12/91
A. J. VARZESINSKI	8/12/91
D. B. KODJICKI	7/12/91



## SECTION 8 - INPUT/OUTPUT CARD SUMMARY OF OPERATION

### I/O Card Summary of Operation

Through the I/O Card, the CPU:

- reads the status of the Control Panel, external switches, and system events
- operates the Control Panel displays and external connections
- synchronizes and times the injection

Refer to page 8 - 2 for a block diagram of the I/O Card. The following circuit groups comprise the I/O Card:

#### Bus Buffers

The control bus is buffered by U14, and the address bus is buffered by U15. The bidirectional data bus is buffered and controlled by U22, with direction control from I/O decoder U8.

#### Display Decoder

Signals for the Display Card on the Control Panel are provided by 4:16 decoder U1, while U7 inserts "Wait States" for the display drivers. These signals control the Sentinel and the display drivers for the LED displays.

#### Display Buffers

Two address lines, and eight data lines, are buffered for the Display Card by U16 and U19.

#### I/O Decoder

Decoder U8, controlled by address and control lines, provides an enable to the programmable I/O (PIO) chips U5 and U6, the counter/timer chip (CTC) U4, and, provides data direction control for the data buffer U22.

#### Parallel Input/Output Chips

Parallel I/O chips permit the CPU to access external inputs, and to control outputs.

PIO1 — The inputs to U6 are the hand start switch, remote start switch, and remote disarm input. Through the data buses, the CPU can check the digital status of these inputs.

Outputs from U6 include: The film changer relay enable, two injecting signals, Armed signal, and Head Indicator signals. Through the I/O chip, the CPU can turn these outputs on and off.

PIO2 — This device (U5) is dedicated to decoding of the keypad on the Control Panel. Port A (sense line) rotates a Logic-0 on each line on the keypad matrix. When a key is pressed, one line of Port B will go to Logic-0. The CPU compares Port A to Port B to determine which key was pressed.

Port A lines A0 - A3: for the Main Keypad  
Port A lines A4 - A7: for the ECG/Scope Keypad (not used in any MCT units)

# Input / Output Card

## Counter/Timer Chip (CTC)

The CTC (U4) generates interrupt signals for the CPU. The interrupts are prioritized and asynchronous with the system clock.

- **Trigger 0:** (Not used with MCT units, SYS 500-PI, or RMT 505-PI/510-PI) R-wave detector interrupt request. When an R-wave is detected, an interrupt is generated. Trigger 0 has the highest priority and is used only with the ECG option.
- **Trigger 1:** PPI memory error interrupt request. PPI primary and backup memories are compared each time primary memory is read. If primary and backup data values are not equal, an interrupt request is generated.
- **Trigger 2:** UART (Universal Asynchronous Receiver-Transmitter) interrupt request. This indicates that the Remote Panel is requesting communication with the Main Unit.
- **Trigger 3:** Timed Interrupt. This provides a real-time clock and delay timer function. During an injection, this interrupt program performs Flow Rate and Volume Limit calculations. This interrupt is not externally triggered, but driven from the time divided clock generator.

## Test Decoder

This 4:16 decoder (U2) provides Reset pulses for Watchdog 1, and controls timing of the buzzer.

## Test Section

The circuit surrounding U17 and U18 is for in-house testing.

## ECG Detector (Not used in any MCT units)

This circuitry is used only with the ECG option. The incoming R-wave signal is sent to the multiplying DAC circuit (U20, U21, U23, and U24 or U20, U21, and U26), while sensitivity to this signal is controlled by the Gain Up/Down keys on the ECG Keypad. Gain information is sent through the data bus to U20, which sets a multiplication factor for DAC U21. The output of the DAC is buffered (at U24 or U26A), then sent to the R-wave detector U23 or U26B, of which output is a pulse that indicates an R-wave. This pulse goes through a pulse shaping network (U25 or U27), then to CTC Trigger 0.

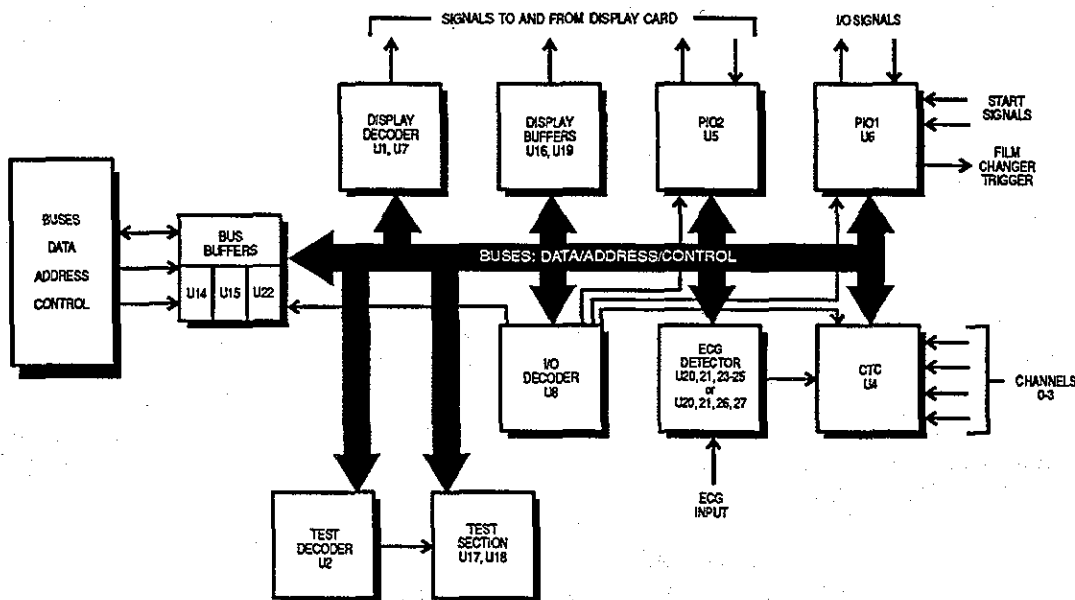
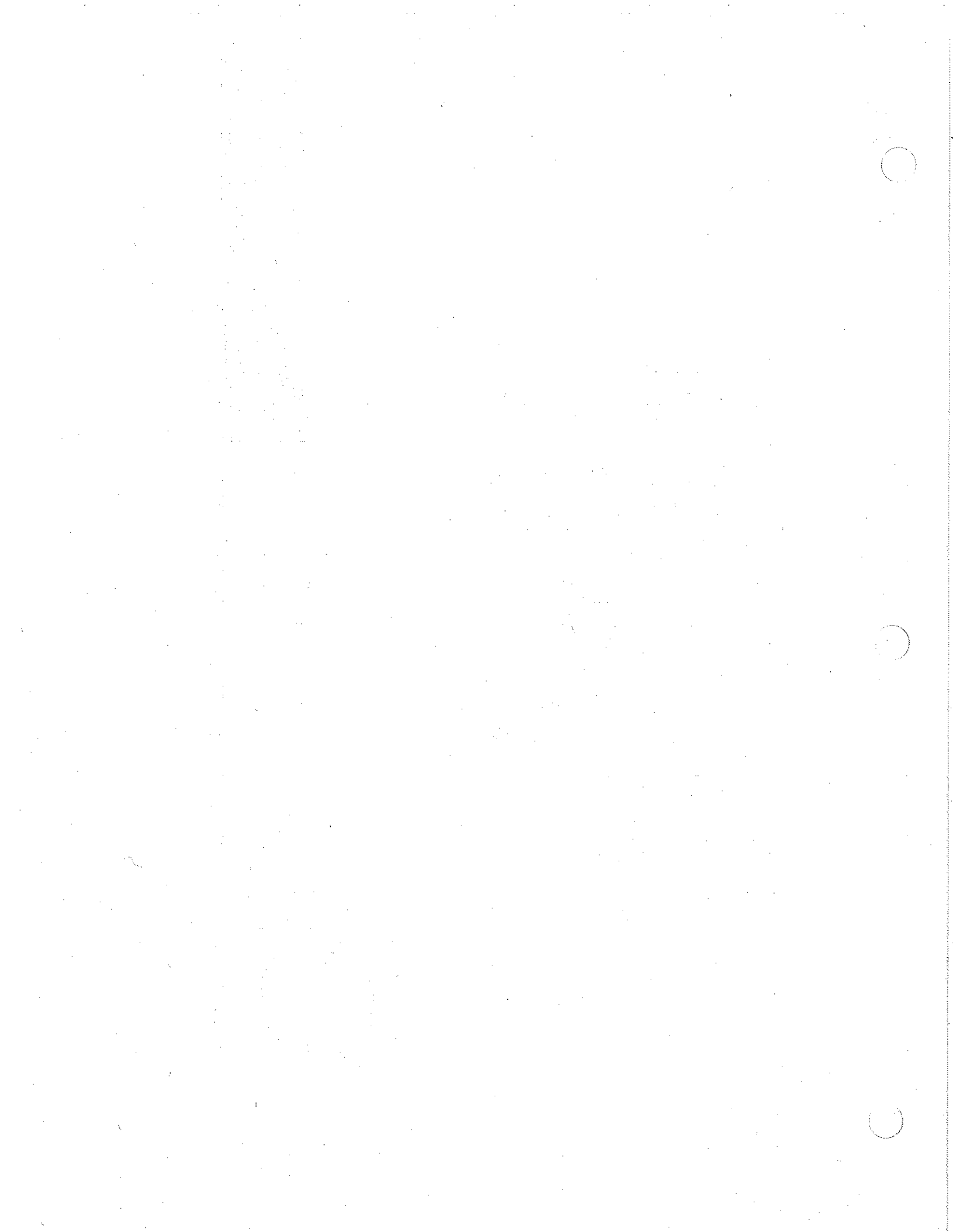



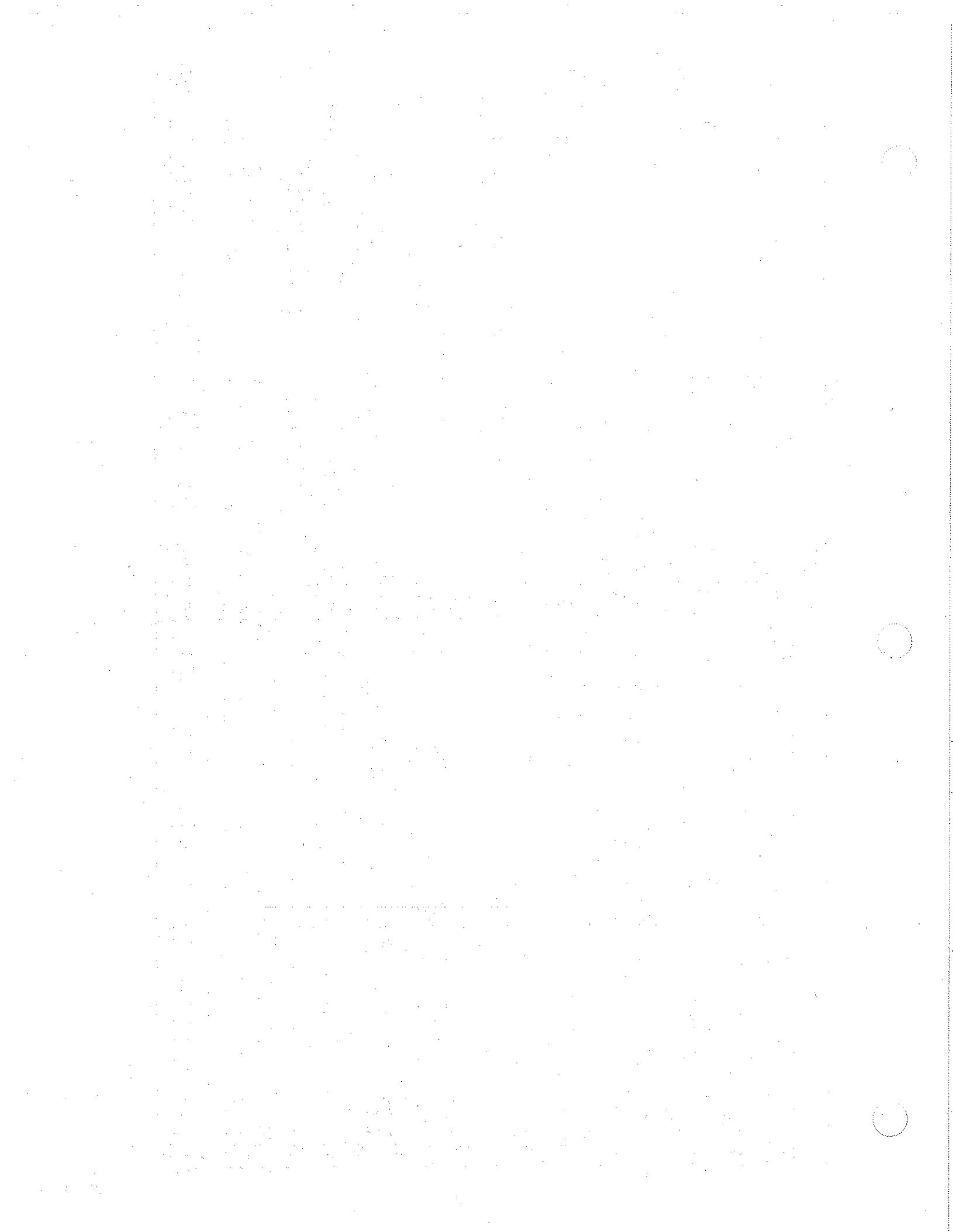
Figure 8.1: Block Diagram: I/O Card






# BILL OF MATERIAL

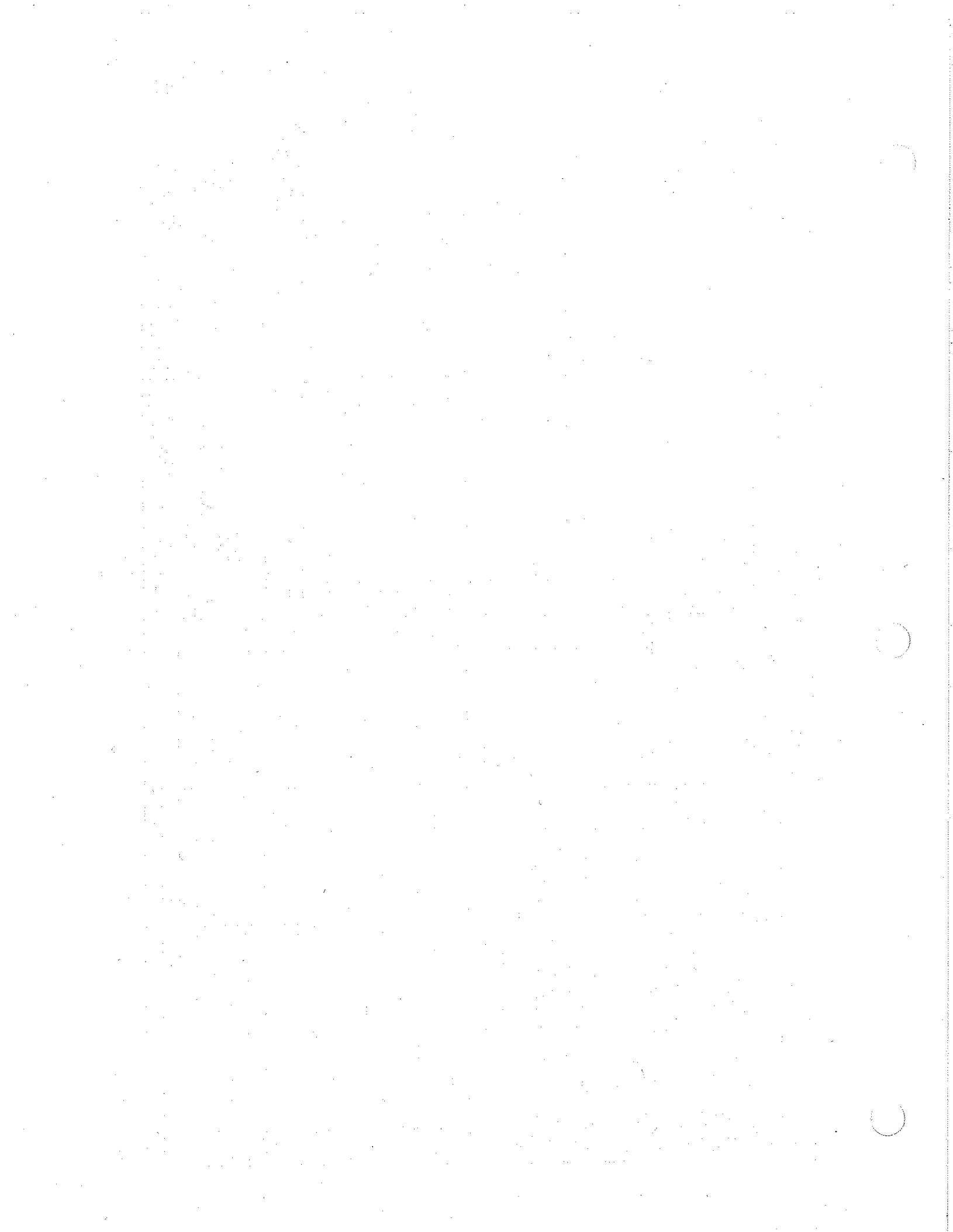ITEM	DESCRIPTION	PART NUMBER	QTY.	REF. DESIG.	PROCESS CODE	PREP CODE	ITEM	DESCRIPTION	PART NUMBER	QTY.	REF. DESIG.	PROCESS CODE	PREP CODE
1	PC CARD, (REV O & UP)	78101-10-PC-02	1		LG		37	DIODE, 1N4148	105-4148-000	7	D4-D10	AMI	5 BEND
2	CARD EJECTOR, WHITE	510-7000-000	1		AW		38						
3	CARD EJECTOR PIN	510-7001-000	2		AW		39						
4	CARD EJECTOR, VIOLET	510-7003-000	1		AW		40	RESISTOR, 2.7 KOHM, 1/4W, 5%	303-0272-000	2	R7,R15	AMI	5 BEND
5	LABEL, BLANK EPROM	550-0017-619	1	U8	AW	NOTE 1	41	RESISTOR, 10 KOHM, 1/4W, 5%	303-0103-000	5	R9-R12,R21	AMI	5 BEND
6							42	RESISTOR, 2.2 MOHM, 1/4W, 5%	303-0225-000	1	R8	AMI	5 BEND
7							43	RESISTOR, 200 OHM, 1/4W, 5%	303-0201-000	2	R23,R24	AMI	5 BEND
8	CONNECTOR, 40 PIN HEADER	518-0208-040	1	J12	BW		44	RESISTOR, 1 KOHM, 1/4W, 5%	303-0102-000	1	R22	AMI	5 BEND
9	CONNECTOR, 34 PIN HEADER	518-0208-034	1	J11	BW		45	RESISTOR, 4.7 KOHM, 1/4W, 5%	303-0472-000	2	R1,R25	AMI	5 BEND
10							46	RESISTOR, 39 KOHM, 1/4W, 5%	303-0393-000	1	R2	AMI	5 BEND
11	LABEL, BLANK SERIAL NUMBER	550-0093-652	1		BW	NOTE 3	47	RESISTOR, 1.5 MOHM, 1/4W, 5%	303-0155-000	1	R26	AMI	5 BEND
12	TUBULAR TERMINAL	642-1238-000	8	46, 80, 800, 100, 201, 202, 115, 100	BW		48						
13	POLARIZING KEY	518-3518-000	4	J11, J12	BW		49						
14	IC SOCKET, 20 PIN	545-8201-020	1	U8	LG		50						
15							51						
16							52	RESISTOR ARRAY, 10 KOHM, SIP	340-4310-103	3	RN3, RN6, RN7	LG	
17							53						
18							54						
19							55						
20	FERRITE BEAD	280-1625-000	1	FB1	LG	5 BEND	56						
21	CAPACITOR, 10µF, 35V	204-0106-135	4	C1, C2, C3, C4	LG	.1HS	57						
22	CAPACITOR, .47µF, 35V	206-0474-350	1	C20	LG	.1HS	58						
23	CAPACITOR, .01µF, 50V	201-0103-050	23	65, 66, 611, 614, 616, 619, 621, 621	LG	.1HS	59	IC, 4-16 DECODER, 74LS154	074-7415-004	2	U1, U2	LG	
24	CAPACITOR, .1µF, 50V	201-0104-050	5	C6, C7, C12, C13, C15	LG	.1HS	60	IC, PIO, Z8420	070-8420-000	2	U5, U6	LG	
25							61	IC, CTC, Z8430	070-8430-000	1	U4	LG	
26							62						
27													
28							All technical information equipment design and application data contained within this document is the CONFIDENTIAL Property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.						
29							TITLE <b>MK V INJECTOR I/O BOARD ASSEMBLY</b>						
30	FUSE, 1/16A, PICO	476-0062-000	2	F1, F2	LG	5 BEND	SIGNATURES DRAWN B. J. O'HELLA CHECKED A. J. YARZEMSKI ENGR. I. H. BRAUNSTEIN		DATE 10-26-90 10-29-90 10-29-90		SEE SHT. 1 FOR REV. LEVEL		
31							N/A		DATE COPIED 78101-10-AP-02		SHEET 2 OF 3		
32							REV. NUMBER 78101-10-AP-02		REV. W		3		
33							120-4401-000		1		Q1		
34							2N4401		1		Q1		
35	TRANSISTOR, 2N4401	120-4401-000	1	Q1	BW		120-4401-000		1		Q1		
36							120-4401-000		1		Q1		









## SECTION 9 - SERVO CONTROL CARD SUMMARY OF OPERATION

### Summary of Operation

The Servo Control Card serves two primary functions:

1. Provide signals for the Flow Rate and Power Drive circuits
2. Read and interpret information from the head, to control injection Flow Rate and Pressure [limit].

Refer to the figure on pages 9 - 3 and 9 - 4 for a block diagram of the Servo Control Card.

### Bus Buffers

The control bus is buffered by U12, and the address bus is buffered by U11. The bi-directional data bus is buffered and controlled by U10, with direction control coming from direction logic.

### Direction Logic

Gates U7-U9 control the direction line of U10, controlling the flow of data to and from the card.

### A/D Converter

Analog-to-digital (A/D) converter U1 inputs a multiplexed analog signal, and converts it into a digital word. This process is controlled by the A/D converter logic (U4, U5) and the DAC decoder (U14). When requested from the CPU, the converted 10-bit value is gated to the data bus by tri-state buffers U2 and U3.

### A/D Converter Logic

Gates U4 and U5 control the A/D converter U1, starting the A/D conversion when the data is ready, and enabling buffers U2 and U3 to transfer the converted word to the data bus.

### DAC Decoder

Decoder U14, with U18 and U6, provide enable signals for the DACs (digital-to-analog converters) and the MUX (multiplexer). These enables allow data bus information to be written to the device. The decoder also provides signals for the A/D converter logic, direction logic, and Reset pulses for WATCHDOG 2.

### WATCHDOG 2

The DAC decoder sends pulses to WATCHDOG 2 (U15) when the mechanical stop is in position, and "Ready to Inject". If the Mechanical Stop takes too long to move into position, the pulses from the DAC decoder will stop, resulting in a disarm condition. The Control Panel will then display the message code SC222 (*Mechanical Stop Position Failure*).

### PIO3

The inputs to PIO3 (U16) are signals from the Injector Head, Power Drive Circuit, and the analog portion of the SCC. These inputs include: Turret switch position, 60 ml indication, and FWD/REV indication from the Injector Head, Aux Monitor from the PDC, and an Overpressure Limit signal from the SCC. The status of PIO3 inputs are directed to the CPU through data busses. Outputs from PIO3 include: Flow Control signals, the Safe Relay Enable, and the System Override signal.

### System Monitor Flip-Flop

The flip-flop (U17) sends alarm signals to disarm the injector, if any of the following input conditions occur: WATCHDOG 2, OPLIM (overpressure signal), or AUX MONITOR (power drive failure).

# Servo Control Card

## Command, Position, and Signal DACs

When enabled and timed by the DAC decoder U14, one of the DACs U21-U24 will read a byte of data bus information. The DAC then converts the byte to an analog voltage between zero and 10 VDC. DAC resolution is 40mV per count.

### Pressure Limit Command DAC

The output of this DAC is the Pressure Limit Command signal (PSILIM SELECTED). This signal is sent to the Multiplexer, the Primary Pressure Limit circuit, and the Overpressure Limit circuit.

### Flow Rate Command DAC

The output of this DAC is the Flow Rate Command signal (FLOW RATE SELECTED). This signal is sent to the Multiplexer, Flow Scale Circuit and Power Drive Card (PDC or PDCP).

### Mechanical Stop Position DAC

The output of this DAC is the Mechanical Stop Position signal (MSPOSCMD). This signal is sent to the Multiplexer, and through Unity Gain Buffer U26C, to the Mechanical Stop Drive Card.

### Velocity Signal DAC

The output of this DAC is a digitized velocity signal, sent to the Imaging System Interface Card. This signal is buffered, becoming the Flow Profile Signal available through connector J15, pin 4.

## Multiplexer (MUX)

The multiplexer circuit is comprised of multiplexer U19, Quad latch U20, and Unity-Gain buffer U26A. This circuit receives eight inputs, and under timing control connects each input to the circuit output, which is sent to the A/D converter.

## TTL/CMOS Converter

U25 converts TTL logic (0/5V) to CMOS logic (0/15V).

## 10V Reference

Chip U28 develops +10 VDC as a reference voltage for the Plunger and Mechanical Stop position feedback pots in the head. This voltage is also applied to the multiplexer to be monitored by the CPU.

## Pot Processor

This circuit, comprised of U27C and associated components, buffers and filters the pot wiper signal. The pot signal is an input to the multiplexer and the error amplifier.

## Primary Pressure Limit Circuit

The Primary pressure limit circuit (U26B), compares the pressure limit command to the actual pressure developed during an injection. The output of this circuit is sent to three Pressure Limiting circuits to perform the following functions:

- Flow Rate Integrator Circuit U29. The output of U26B is summed with the Flow Rate command to effectively reduce the input voltage to U29. This action results in the reduction of the Flow Rate command signal, thereby preventing the selected pressure limit from being exceeded.
- Inverting amplifier, U27B. The output of this amplifier is summed with the output of the Main Flow Integrator, to reduce the voltage sent to the error amplifier U27D, thereby decreasing the error signal.
- Circuit Q4. This circuit is driven to cut-off, sending a pressure limit indication signal (PRESLIMIND) to PIO3 (U16), and eventually to the CPU. This PRESLIMIND signal will also send Q12 into a saturation state, which shunts R81 and reduces the gain of the error amplifier, thereby decreasing the error signal.

# Servo Control Card

## Overpressure Circuit

If the actual pressure exceeds the pressure limit command by greater than 100-150 PSI, the Overpressure Circuit (U26D) forwards an OPLIM signal to stop the injection.

## Flow Scale Circuit

Using the two flow scale signals, circuits Q5, Q6, and Q11 attenuate the flow rate selected signal during ML/MIN and ML/HR flow scale injections. (ML/MIN and ML/HR does not apply to MCT units)

## Standby Reset Circuit

Driven by the INJECT signal, Q8, Q9, and Q10 control the flow rate integrator, U29. During standby, the integrator follows the position pot. The output of the integrator is equal in amplitude, and opposite in polarity to the position pot. When an injection is to occur, this circuit allows U29 to generate the position command.

## Flow Rate Integrator

The master position command is generated by integrator U29. The output of this circuit is a ramp, with a slope proportional to the input signal (flow rate selected).

## Error Amplifier

When the plunger is moving at the desired flow rate, the master position command from U29 and the pot wiper signal from U27C are equal but opposite in polarity. The output of error amplifier U27D is the difference between these two signals, or the position error. This error signal is used by the drive circuits to power the motor. If pressure limiting occurs, a signal from U27B reduces the error signal, and the gain of U27D is reduced by Q12.

## Feed Forward Circuit

This circuit, comprised of Q7 and U27A, provides a boost voltage to the motor to compensate for the internal resistance of the motor. This turning voltage varies with the Flow Rate.

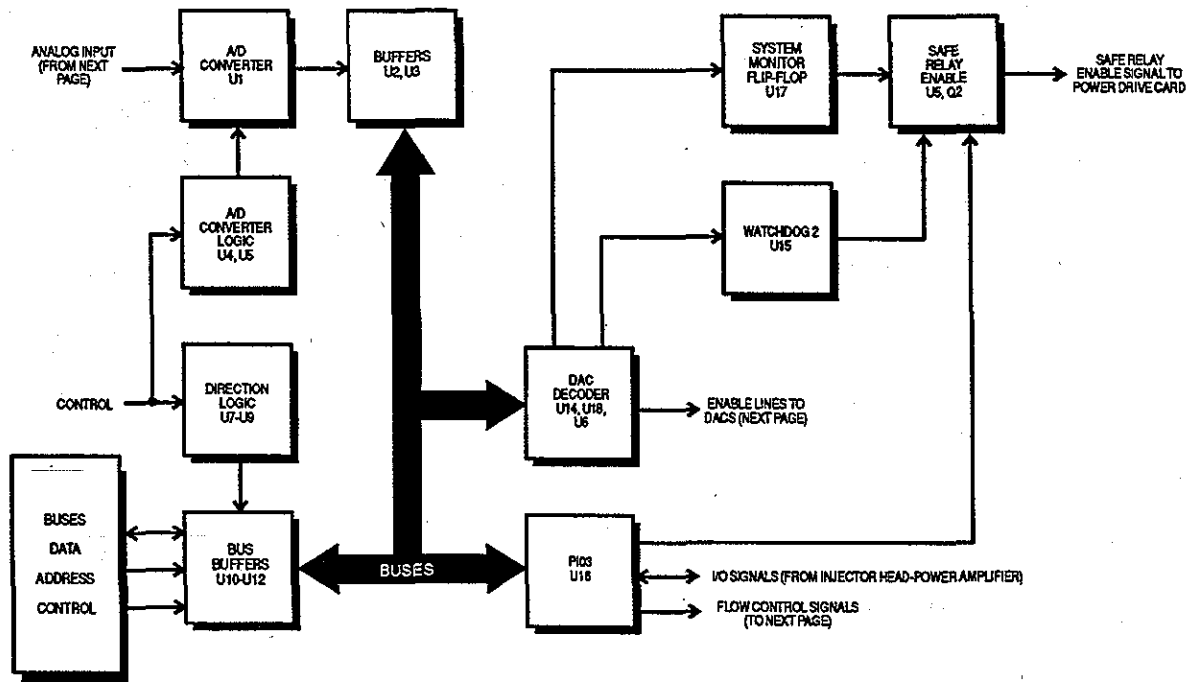


Figure 9.1: Block Diagram: Servo Control Card (sheet 1 of 2)

# Servo Control Card


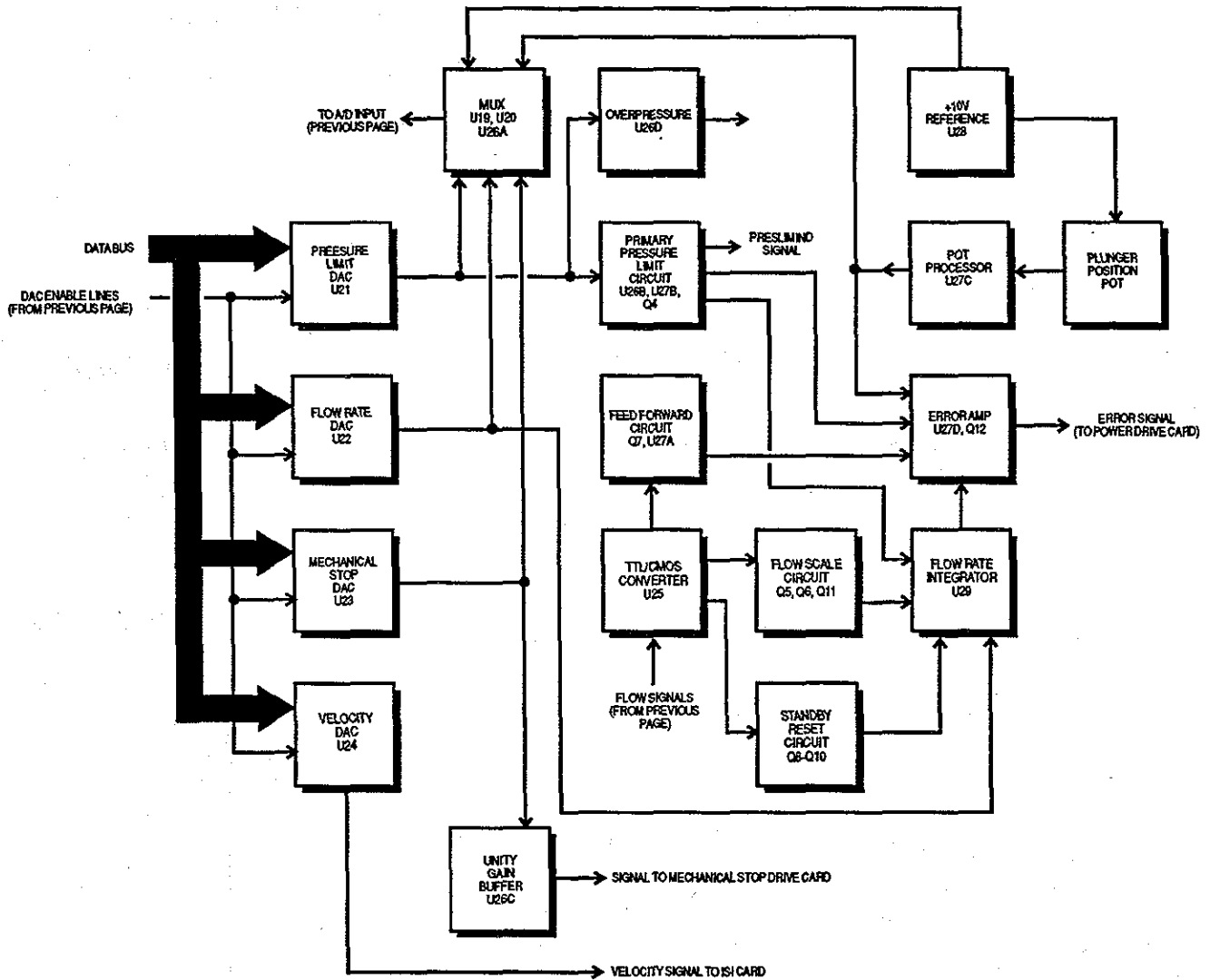


Figure 9.2: Block Diagram: Servo Control Card (sheet 2 of 2)

REV.	DATE	BY	CHKD.
1	11/17/68	W. J. MANN	W. J. MANN
2	1/2/69	J. M. BAK	J. M. BAK
3	1/2/69	J. M. BAK	J. M. BAK
4	1/2/69	J. M. BAK	J. M. BAK
5	1/2/69	J. M. BAK	J. M. BAK
6	1/2/69	J. M. BAK	J. M. BAK
7	1/2/69	J. M. BAK	J. M. BAK
8	1/2/69	J. M. BAK	J. M. BAK

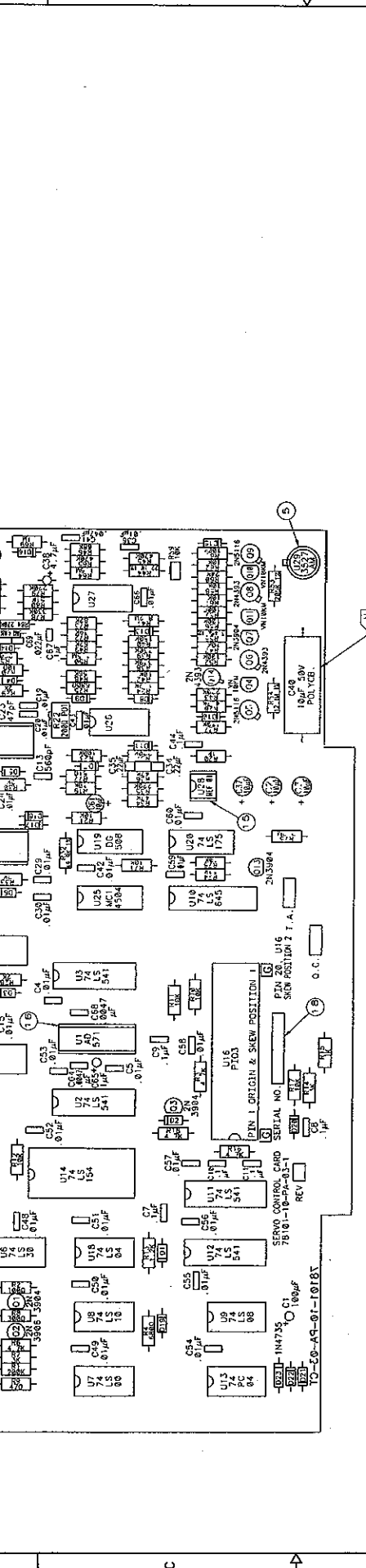REVISIONS  
 1. ORIGINAL DESIGN  
 2. CHANGE TO BOARD FOR TESTING  
 3. CHANGE TO BOARD FOR TESTING  
 4. CHANGE TO BOARD FOR TESTING  
 5. CHANGE TO BOARD FOR TESTING  
 6. CHANGE TO BOARD FOR TESTING  
 7. CHANGE TO BOARD FOR TESTING  
 8. CHANGE TO BOARD FOR TESTING



REV.	DATE	BY	CHKD.
1	11/17/68	W. J. MANN	W. J. MANN
2	1/2/69	J. M. BAK	J. M. BAK
3	1/2/69	J. M. BAK	J. M. BAK
4	1/2/69	J. M. BAK	J. M. BAK
5	1/2/69	J. M. BAK	J. M. BAK
6	1/2/69	J. M. BAK	J. M. BAK
7	1/2/69	J. M. BAK	J. M. BAK
8	1/2/69	J. M. BAK	J. M. BAK

REVISIONS  
 1. ORIGINAL DESIGN  
 2. CHANGE TO BOARD FOR TESTING  
 3. CHANGE TO BOARD FOR TESTING  
 4. CHANGE TO BOARD FOR TESTING  
 5. CHANGE TO BOARD FOR TESTING  
 6. CHANGE TO BOARD FOR TESTING  
 7. CHANGE TO BOARD FOR TESTING  
 8. CHANGE TO BOARD FOR TESTING

REV.	DATE	BY	CHKD.
1	11/17/68	W. J. MANN	W. J. MANN
2	1/2/69	J. M. BAK	J. M. BAK
3	1/2/69	J. M. BAK	J. M. BAK
4	1/2/69	J. M. BAK	J. M. BAK
5	1/2/69	J. M. BAK	J. M. BAK
6	1/2/69	J. M. BAK	J. M. BAK
7	1/2/69	J. M. BAK	J. M. BAK
8	1/2/69	J. M. BAK	J. M. BAK

REVISIONS  
 1. ORIGINAL DESIGN  
 2. CHANGE TO BOARD FOR TESTING  
 3. CHANGE TO BOARD FOR TESTING  
 4. CHANGE TO BOARD FOR TESTING  
 5. CHANGE TO BOARD FOR TESTING  
 6. CHANGE TO BOARD FOR TESTING  
 7. CHANGE TO BOARD FOR TESTING  
 8. CHANGE TO BOARD FOR TESTING

PROCESS COMPATIBILITY CHART

ITEM	DO NOT DESIC.	WAVE SOLDER	DO NOT VAPOR DECREASE	DO NOT WATER WASH
ITEM 2	X			
ITEM 3	X			
ITEM 4	X			
S1.2	X		X	X
U29	X		X	X
C49	X		X	X

- A. INSTALL COMPONENTS MARKED "DO NOT WAVE SOLDER" AFTER WAVE SOLDER. HAND CLEAN COMPONENTS MARKED "DO NOT VAPOR DECREASE OR WATER WASH". MASK ALL COMPONENT LOCATIONS MARKED "DO NOT WAVE SOLDER".
- B. ALL UNMARKED DIODES 1N4148.
- C. INSTALL PROGRAMMING SWITCH (ITEM 1) AS DETILED ABOVE FOR NUMBER "0" DURING ASSEMBLY. A "WALL SET" PROGRAMMING SWITCHES DURING TEST AS REQUIRED TO THE "AP" OR "CT" MODES AS SHOWN BELOW:  
 'AP' 2 0 0 0 0 0 0 1 1 3  
 'CT' 2 0 0 0 0 0 0 1 4
- D. PRODUCTION TO FILL IN REV. BLOCK USING PERMANENT INK WITH THE ASSEMBLY DRAWING REVISION TO WHICH THE BOARD WAS BUILT. SERIAL NO. BLOCK TO BE FILLED USING PERMANENT INK OR LABEL (ITEM 16).
- E. ALL SOCKETED COMPONENTS TO BE INSTALLED AFTER WAVE SOLDER.
- F. LOCATIONS MARKED "G" ARE USED TO INDICATE ORIGIN AND SKEW LOCATIONS.
- G. USE ADHESIVE TO FASTEN C49 TO THE BOARD. USE RTV OR EQUIVALENT REMOVABLE SILICON BASED ADHESIVE.

REV.	DATE	BY	CHKD.
1	11/17/68	W. J. MANN	W. J. MANN
2	1/2/69	J. M. BAK	J. M. BAK
3	1/2/69	J. M. BAK	J. M. BAK
4	1/2/69	J. M. BAK	J. M. BAK
5	1/2/69	J. M. BAK	J. M. BAK
6	1/2/69	J. M. BAK	J. M. BAK
7	1/2/69	J. M. BAK	J. M. BAK
8	1/2/69	J. M. BAK	J. M. BAK

REVISIONS  
 1. ORIGINAL DESIGN  
 2. CHANGE TO BOARD FOR TESTING  
 3. CHANGE TO BOARD FOR TESTING  
 4. CHANGE TO BOARD FOR TESTING  
 5. CHANGE TO BOARD FOR TESTING  
 6. CHANGE TO BOARD FOR TESTING  
 7. CHANGE TO BOARD FOR TESTING  
 8. CHANGE TO BOARD FOR TESTING

REV.	DATE	BY	CHKD.
1	11/17/68	W. J. MANN	W. J. MANN
2	1/2/69	J. M. BAK	J. M. BAK
3	1/2/69	J. M. BAK	J. M. BAK
4	1/2/69	J. M. BAK	J. M. BAK
5	1/2/69	J. M. BAK	J. M. BAK
6	1/2/69	J. M. BAK	J. M. BAK
7	1/2/69	J. M. BAK	J. M. BAK
8	1/2/69	J. M. BAK	J. M. BAK

REVISIONS  
 1. ORIGINAL DESIGN  
 2. CHANGE TO BOARD FOR TESTING  
 3. CHANGE TO BOARD FOR TESTING  
 4. CHANGE TO BOARD FOR TESTING  
 5. CHANGE TO BOARD FOR TESTING  
 6. CHANGE TO BOARD FOR TESTING  
 7. CHANGE TO BOARD FOR TESTING  
 8. CHANGE TO BOARD FOR TESTING

REV.	DATE	BY	CHKD.
1	11/17/68	W. J. MANN	W. J. MANN
2	1/2/69	J. M. BAK	J. M. BAK
3	1/2/69	J. M. BAK	J. M. BAK
4	1/2/69	J. M. BAK	J. M. BAK
5	1/2/69	J. M. BAK	J. M. BAK
6	1/2/69	J. M. BAK	J. M. BAK
7	1/2/69	J. M. BAK	J. M. BAK
8	1/2/69	J. M. BAK	J. M. BAK

REVISIONS  
 1. ORIGINAL DESIGN  
 2. CHANGE TO BOARD FOR TESTING  
 3. CHANGE TO BOARD FOR TESTING  
 4. CHANGE TO BOARD FOR TESTING  
 5. CHANGE TO BOARD FOR TESTING  
 6. CHANGE TO BOARD FOR TESTING  
 7. CHANGE TO BOARD FOR TESTING  
 8. CHANGE TO BOARD FOR TESTING

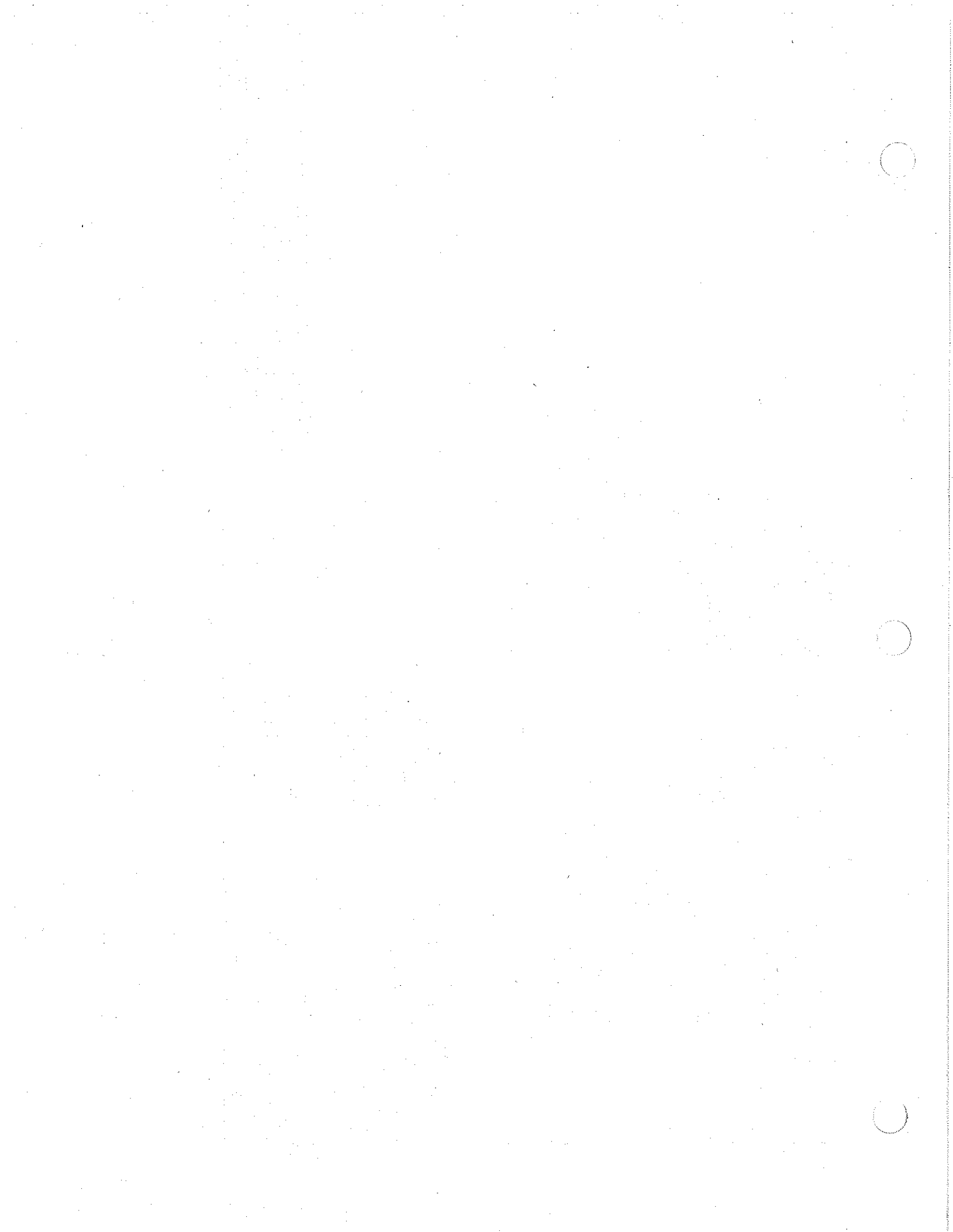
9-5

**MEDRAD, INC.**  
 DRAY 15238  
 PITTSBURGH, PA.

MKX MCT INJECTOR  
 SERVO CONTROL CARD ASSEMBLY

78101-10-AP-03-1

PAGE RTS





# BILL OF MATERIAL

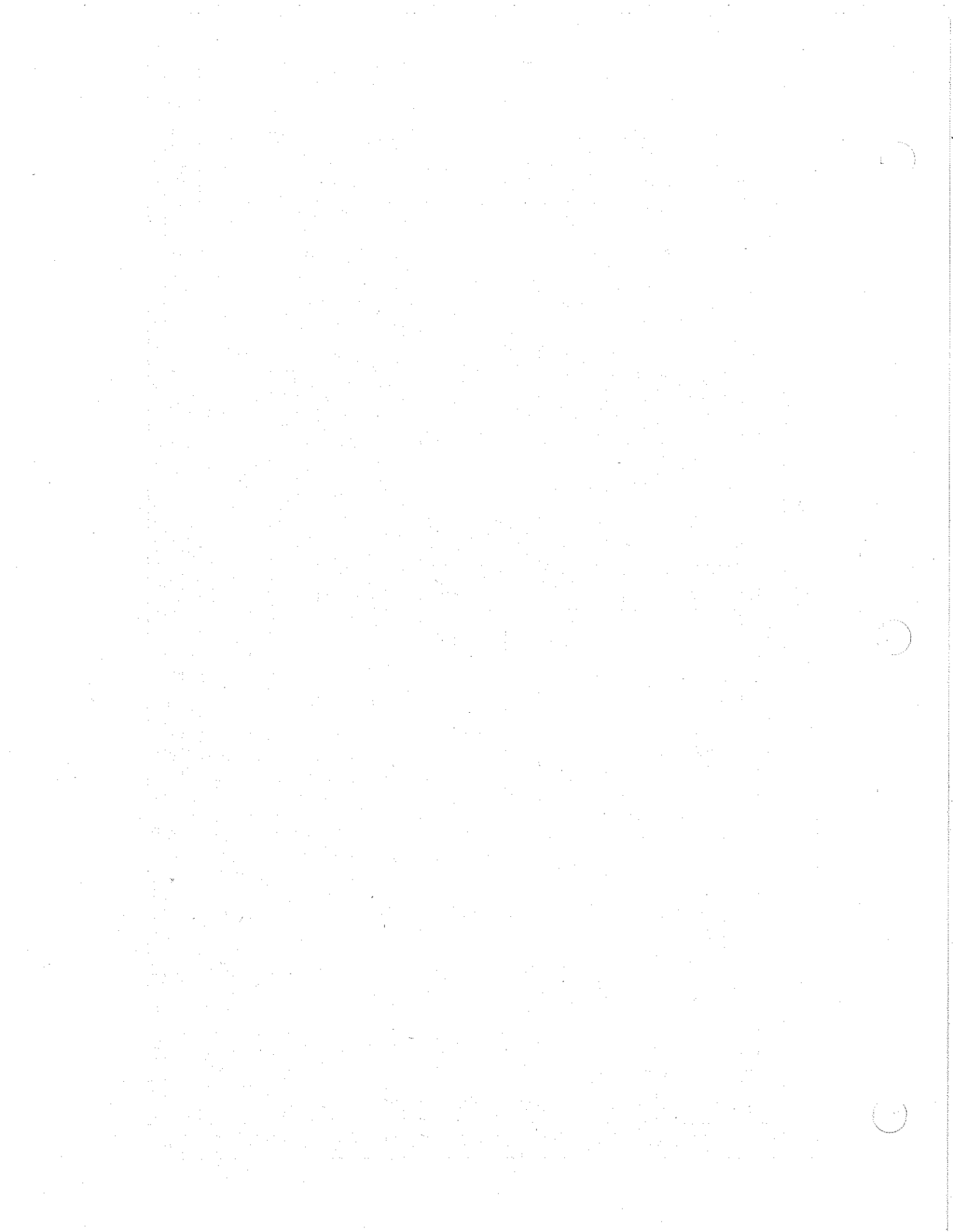
ITEM	DESCRIPTION	PART NUMBER	QTY.	PROCESS CODE	REF. DESIG	ITEM	DESCRIPTION	PART NUMBER	QTY.	PROCESS CODE	REF. DESIG
1	PC CARD (REV L AND UP)	78101-10-PC-03	1	LG		50	RESISTOR, 47 OHM 1/4 W 5%	303-0470-000	1	AXI	5 BEND R9
2	CARD EJECTOR (WHITE)	510-7000-000	1	AW		51	RESISTOR, 200 KOHM 1/4 W 5%	303-0204-000	1	AXI	5 BEND R1
3	CARD EJECTOR PIN	510-7001-000	2	AW		52	RESISTOR, 3 KOHM 1/4 W 5%	303-0302-000	2	AXI	5 BEND R7, 14
4	CARD EJECTOR (GREEN)	510-7004-000	1	AW		53	RESISTOR, 4.7 KOHM 1/4 W 5%	303-0472-000	5	AXI	5 BEND R3, 5, 6, 18, 26
5	SPACER, IC PAD, 8 PIN	600-0615-075	1	AW	U29	54	RESISTOR, 300 KOHM 1/4 W 5%	303-0301-000	1	AXI	5 BEND R8
6	SWITCH, SPST	450-4356-661	2	AW	S1, 2	55	RESISTOR, 100 OHM 1/4 W 5%	303-0101-000	2	AXI	5 BEND R2, 57
7	PROGRAMMING SHUNT	518-5312-202	4	AW	NOTE B	56	RESISTOR, 680 OHM 1/4 W 5%	303-0681-000	2	AXI	5 BEND R4, 25
10	HEADER, 14 PIN	500-1070-700	1	AW	JU1	57	RESISTOR, 1.2 KOHM 1/4 W 5%	303-0122-000	2	AXI	5 BEND R13, 24
11	TUBULAR TERMINAL	642-1238-000	6	BW		58	RESISTOR, 10 KOHM 1/4 W 5%	303-0103-000	13	AXI	5 BEND R11, 12, 14, 22, 47, 58, 61, 62, 71, 76
15	IC SOCKET, 8 PIN	545-0008-000	1	LG	U28	59	RESISTOR, 1 KOHM 1/4 W 5%	303-0102-000	2	AXI	5 BEND R15, 20
16	IC SOCKET, 18 PIN	545-8181-018	1	LG	U1	60	RESISTOR, 2 KOHM 1/4 W 5%	303-0202-000	4	AXI	5 BEND R28, 29, 30, 33
17	IC SOCKET, 22 PIN	545-8221-022	2	LG	U22, 23	61	RESISTOR, 10 OHM 1/4 W 5%	303-0100-000	1	AXI	5 BEND R70
18	LABEL, BLANK SERIAL NUMBER	550-0093-652	1	BW	NOTE E	62	RESISTOR, 220 KOHM 1/4 W 5%	303-0224-000	5	AXI	5 BEND R35, 37, 43, 51, 79
20	CAP, 10uf 50V POLYCARBONATE	209-0106-500	1	AW	14 BEND C40	63	RESISTOR, 100 KOHM 1/4 W 5%	303-0104-000	5	AXI	5 BEND R23, 27, 52, 56, 63
21	CAP, 560pF, 50V	201-0561-050	1	LG	1-LS C13	64	RESISTOR, 51 OHM 1/4 W 5%	303-0510-000	1	AXI	5 BEND R77
22	CAP, .047uf, 50V	201-0473-050	1	LG	1-LS C41	65	RESISTOR, 49.9 KOHM 1%	313-4992-000	1	AXI	5 BEND R32
23	CAP, .0047uf, 50V	201-0472-050	2	LG	1-LS C64, 68	66	RESISTOR, 4.02 KOHM 1%	313-4021-000	1	AXI	5 BEND R82
24	CAP, .022uf, 50V	201-0223-050	1	LG	1-LS C69	67	RESISTOR, 270 OHM 1/4 W 5%	303-0274-000	1	AXI	5 BEND R84
25	CAP, .47PF, 50V	201-0470-050	4	LG	1-LS C18, 23, 28, 33	68	RESISTOR, 1.21 MOHM 1%	313-1214-000	1	AXI	5 BEND R85
26	CAP, .001uf, 50V	201-0102-050	4	LG	1-LS C17, 22, 27, 32	69	RESISTOR, 22.1 KOHM 1%	313-2212-000	2	AXI	5 BEND R44, 83
27	CAP, .22uf, 50V	201-0224-050	3	LG	1-LS C34, 35, 39	70	RESISTOR, 1 MOHM 1/4 W 5%	303-0105-000	2	AXI	5 BEND R69, 81
28	CAP, .1uf, 50V	201-0104-050	8	LG	1-LS C7, 8, 9, 10, 11, 44, 62, 67	71	RESISTOR, 51 KOHM 1/4 W 5%	303-0513-000	2	AXI	5 BEND R74, 80
29	CAP, .01uf, 50V	201-0103-050	35	LG	1-LS C3, 10, 14, 15, 16, 17, 18, 19, 20, 21, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36, 37, 38, 39, 40, 41, 42, 43, 44, 45, 46, 47, 48, 49, 50, 51, 52, 53, 54, 55, 56, 57, 58, 59, 60, 61, 62, 63, 64, 65, 66, 67, 68, 69, 70, 71, 72, 73, 74, 75, 76, 77, 78, 79, 80, 81, 82, 83, 84, 85, 86, 87, 88, 89, 90, 91, 92, 93, 94, 95, 96, 97, 98, 99, 100	72	RESISTOR, 330 KOHM 1/4 W 5%	303-0334-000	3	AXI	5 BEND R38, 39, 78
30	CAP, 2.2uf, 35V	206-0225-350	1	LG	1-LS C6	73	RESISTOR, 150 KOHM 1/4 W 5%	303-0154-000	1	AXI	5 BEND R40
31	CAP, 100uf, 10V	204-0107-010	1	LG	1-LS C1	74	RESISTOR, 8.2 MOHM 1/4 W 5%	303-0825-000	1	AXI	5 BEND R41
32	CAP, 1uf, 35V	206-0105-350	1	LG	1-LS C65	75	RESISTOR, 470 KOHM 1/4 W 5%	303-0474-000	2	AXI	5 BEND R45, 65
33	CAP, 10uf, 35V	206-0106-351	4	LG	1-LS C2, 3, 37, 61		All technical information equipment design and application data contained within this document is the CONFIDENTIAL Property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.				
34	CAP, 4.7uf, 35V	206-0475-350	1	LG	1-LS C38						
40	TRANSISTOR, 2N3904	120-3904-000	4	BW	Q1, 3, 7, 13						
41	TRANSISTOR, 2N3906	121-3906-000	1	BW	Q2						
42	TRANSISTOR, FET, VN10KM	019-0010-000	4	BW	04, 10, 11, 12						
43	TRANSISTOR, 2N5116	175-5116-000	2	BW	Q5, 9						
44	TRANSISTOR, FET, 2N4393	175-4393-000	3	BW	Q6, 8, 14						
45	DIODE, 1N4148	105-4148-000	21	AXI	5 BEND D1-11, 13-22						
46	DIODE, ZENER, 1N4735, 6.2V	102-4735-100	1	BW	5 BEND D23						
47	DIODE, ZENER, 1N4746, 18V	102-4746-100	2	BW	5 BEND D24, D25						

**MEDRAD INC.**  
 271 KAPPA DRIVE  
 PITTSBURGH, PA 15238  
 TITLE

**MARK V INJECTOR SERVO CONTROL CARD ASSEMBLY**

DWO. NUMBER **78101-10-AP-03-1** REV. **JJ**  
 DATE COPIED **SHEET 2 OF 3**

SEE **9-13-90 SHT. 1**  
 FOR **9-25-90**  
 REV. **9-26-90**  
**LEVEL**



# BILL OF MATERIAL

ITEM	DESCRIPTION	PART NUMBER	QTY.	PROCESS CODE	PREP CODE	REF. DESIG	ITEM	DESCRIPTION	PART NUMBER	QTY.	PROCESS CODE	PREP CODE	REF. DESIG
76	RESISTOR, 68 KOHM 1/4 W 5%	303-0683-000	1	AXI	5 BEND	R46	115	IC, D FLIP-FLOP, 74LS175	074-74LS-175	1	DIP		U20
77	RESISTOR, 1.5 MOHM 1/4 W 5%	303-0155-000	2	AXI	5 BEND	R50,64	116	RESISTOR, 4.7K0, 1/4W, 5%	303-0472-000	1	LG		5 BEND R16
78	RESISTOR, 82 KOHM 1/4 W 5%	303-0823-000	1	AXI	5 BEND	R73	117	RESISTOR, 1M0, 1/4W, 5%	303-0105-000	1	LG		5 BEND R66
79	RESISTOR, 39 KOHM 1/4 W 5%	303-0393-000	1	AXI	5 BEND	R72							
80	RESISTOR, 62 KOHM 1/4 W 5%	303-0623-000	1	AXI	5 BEND	R49							
81	RESISTOR, 56 KOHM 1/4 W 5%	303-0563-000	1	AXI	5 BEND	R48							
82	RESISTOR, 47 KOHM 1/4 W 5%	303-0473-000	2	AXI	5 BEND	R34,36							
83	RESISTOR, 24 KOHM 1/4 W 5%	303-0243-000	1	AXI	5 BEND	R60							
84	RESISTOR, 3.48 KOHM 1%	313-3481-000	1	AXI	5 BEND	R55							
85	RESISTOR, 2.2 MOHM 1/4 W 5%	303-0225-000	1	AXI	5 BEND	R75							
86	RESISTOR, 12.1 MOHM 1%	315-1215-000	1	AXI	5 BEND	R54							
87	RESISTOR, 205 KOHM 1%	313-2053-000	1	AXI	5 BEND	R53							
88	POT, 200 OHM 1 TURN	353-0201-010	1	LG		R22							
89	POT, 10 KOHM, 25 TURN	362-0103-010	1	LG		R31							
90	POT, 10 KOHM 1 TURN	353-0103-010	1	LG		R59							
95	IC, A/D CONV, AD571	050-0571-000	1	AW		U1							
96	IC, VOLTAGE REF, REF-01	052-4101-000	1	AW		U28							
97	IC, PRECISION OP AMP, 3527AM	050-3527-070	1	AW		U29							
98	IC, MONOSTABLE, 74LS122	074-74LS-122	1	DIP		U15							
99	IC, QUAD OR GATE, 74LS32	074-74LS-032	1	DIP		U4							
100	IC, D FLIP-FLOP, 74LS74	074-74LS-074	1	DIP		U17							
101	IC, QUAD NAND GATE, 74LS00	074-74LS-000	2	DIP		U5,7							
102	IC, D/A CONV, NE5018	053-5018-000	2	AW		U22,23							
103	IC, D/A CONV, NE5018	053-5018-000	2	DIP		U21,24							
104	IC, 8-INPUT NAND, 74LS30	074-74LS-030	1	DIP		U6							
105	IC, 3-INPUT NAND, 74LS10	074-74LS-010	1	DIP		U8							
106	IC, HEX INVERTER, 74LS04	074-74LS-004	1	DIP		U18							
107	IC, 4-16 DECODER, 74LS154	074-74LS-154	1	LG		U14							
108	IC, BUFFER, 74LS541	074-74LS-541	4	DIP		U2,3,11,12							
109	IC, LEVEL SHIFTER, MC14504B	004-4504-000	1	DIP		U25							
110	IC, ANALOG SWITCH, DG508	051-0508-000	1	DIP		U19							
111	IC, QUAD OP AMP	060-4741-500	2	DIP		U26,27							
112	IC, HEX INVERTER, 74PC04	080-74PC-004	1	DIP		U13							
113	IC, QUAD 2 INPUT AND, 74LS08	074-74LS-008	1	DIP		U9							
114	IC, PIO, Z8420	070-8420-000	1	LG		U16							
115	IC, BUS TRANS, 74LS645	074-74LS-645	1	DIP		U10							

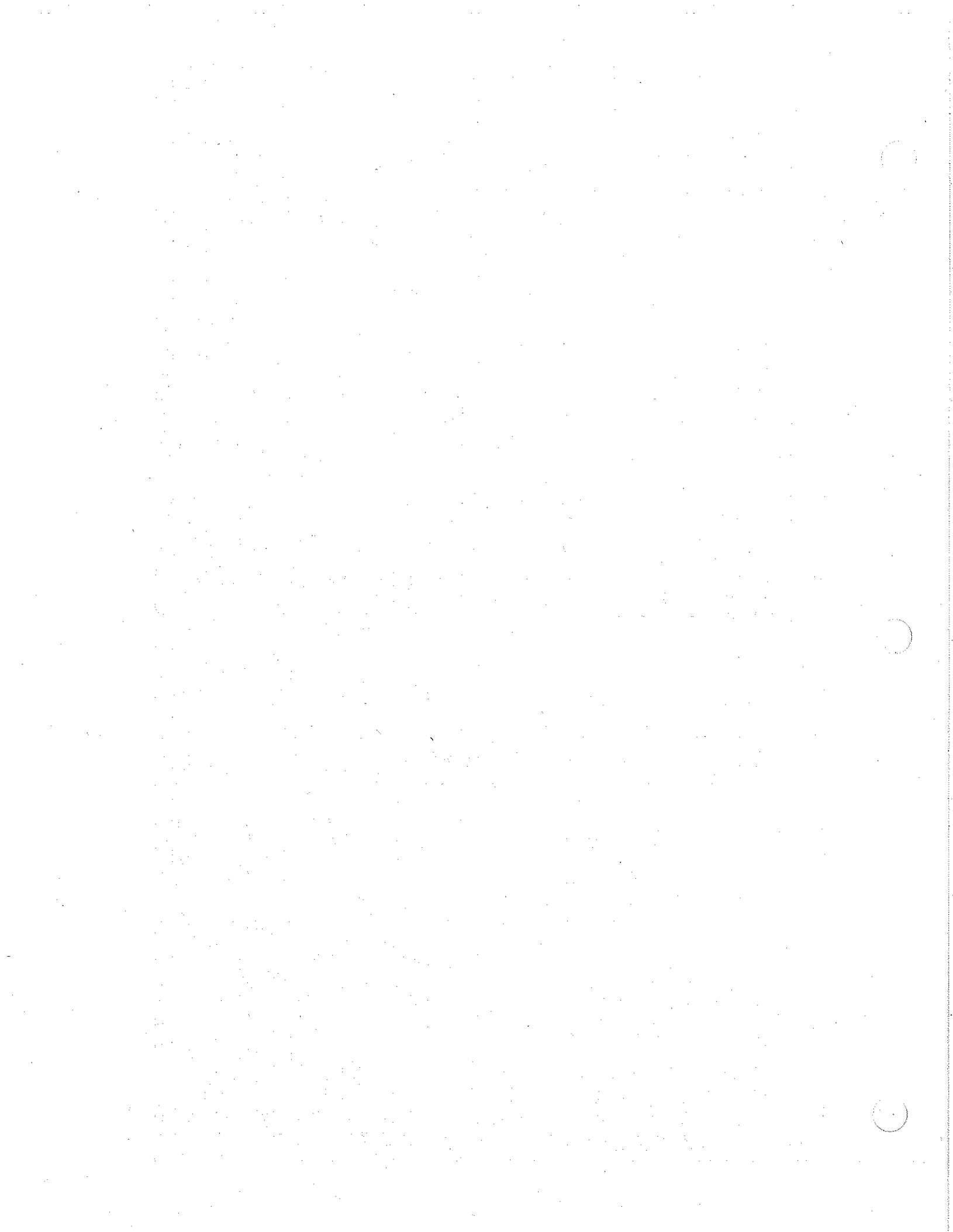
**MEDRAD INC.**  
271 KAPPA DRIVE  
PITTSBURGH, PA 15238

**MARK V INJECTOR SERVO  
CONTROL CARD ASSEMBLY**

DWG. NUMBER: **78101-10-AP-03-1** REV. **JJ**  
DATE COPIED: \_\_\_\_\_ SHEET **3** OF **3**

All technical information equipment design and application data contained within this document is the CONFIDENTIAL Property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.

SIGNATURES	DATE	SEE
WHALEN	9-13-90	SHT. 1
A. VARZEBINSKI	9-25-90	FOR
J. BRAUNSTEIN	9-26-90	REV.
VARZEBINSKI	3-18-92	LEVEL



## SECTION 10 - PPI CARD SUMMARY OF OPERATION

There are two versions of the PPI Card, identified by the card revision level and serial number. Use the chart below to identify card type, then refer to the appropriate description on the following pages.

Revision	Serial Number	PPI Card Type
Revision B	20000 - 25665	Type B
Revision C and above	25666 and Above	Type C

## PPI Card

### PPI Card Type B (serial numbers 20000 to 25665)

For a block diagram of PPI Card type B, refer to page 10 - 3.

#### Bus Buffers

The control bus is buffered by U3, and the address bus is buffered by U1 and U2. The bidirectional data bus is buffered by U4 and U5, with direction defined by one of the lines on the control bus. These data buffers are enabled by the decoder. The data bus for the primary memory chip is controlled by U5; the backup memory chip data bus is controlled by U4.

#### Decoder

Decoder U6 provides the enable lines for the data bus buffers, the memory chips, and the RS-232C circuit.

#### PPI Memory

Primary PPI memory is provided by 8K RAM chip U9. Backup PPI memory is provided by U8. The chips are enabled by signals from the decoder, while outputs are enabled by lines from the control bus.

#### Memory Comparator

The data buses from the primary and backup memory are continuously compared by U7. If the values do not agree, latch U11 sends a signal (MEMERRINT), enabling the CTC on the I/O Card, to interrupt the processor.

#### Analog Switches

The enable lines for the memory chips are switched through U16. When power is applied, the switches are closed, connecting the enable lines. When power is removed or low, the switches are open, disabling the enable lines to preserve power.

#### Power-Down and Battery Circuit

When power is removed, or is low, Q1- Q3 disables U16 and opens the analog switches to disable the memory chips to conserve power as they switch over to battery power. The battery provides backup power for the memory chips to preserve the contents when power is removed.

#### RS-232C Circuit

This circuit consists of Universal Asynchronous Receiver/Transmitter (UART) U13, baud rate clock U12, line receivers U15, and line driver U14. With this circuit, the Main Unit can communicate with the Remote Panel using an RS-232C serial interface.

# PPI Card

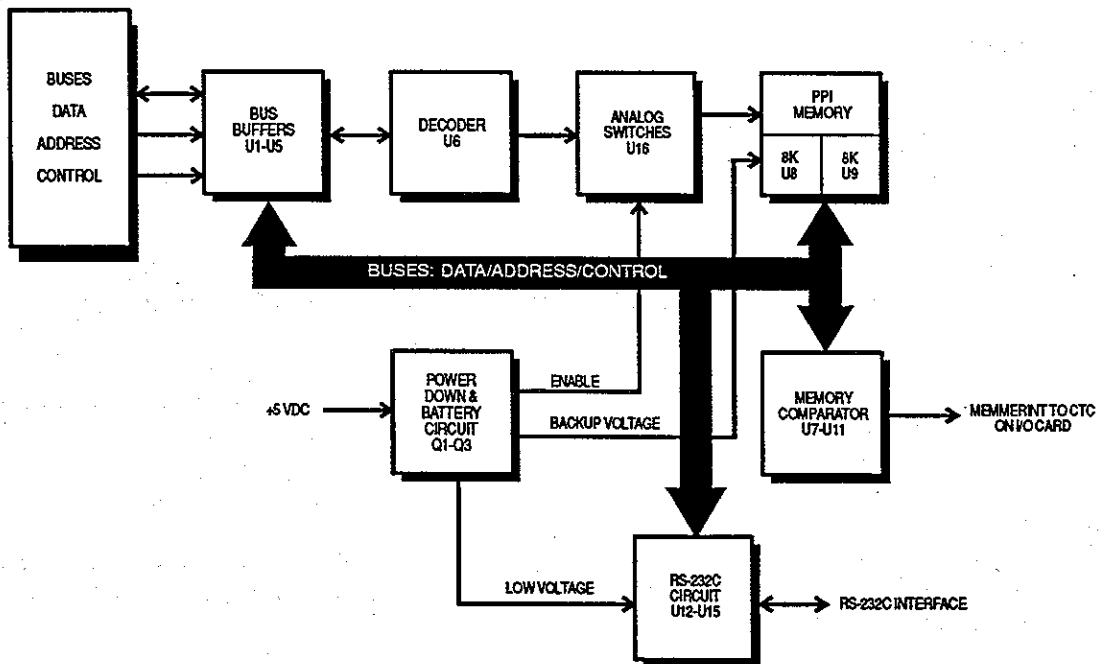


Figure 10.1: Block Diagram: PPI Card Type B

# PPI Card

## PPI Card Type C (serial number 25666 and above)

For a block diagram of PPI Card type C, refer to page 10 - 5.

### Bus Buffers

The control bus is buffered by U3, and the address bus is buffered by U1 and U2. The bidirectional data bus is buffered by U4 and U5, with direction defined by one of the lines on the control bus. These data buffers are enabled by the decoder. The data bus for the primary memory chip is controlled by U5; the backup memory chip data bus is controlled by U4.

### Decoder

Decoder U6 provides the enable lines for the data bus buffers, the memory chips, and the RS-232C circuit.

### PPI Memory

Primary PPI memory is provided by 8K RAM chip U9. Backup PPI memory is provided by U8. The chips are enabled by signals from the decoder, while outputs are enabled by lines from the control bus.

### Memory Comparator

The data buses from the primary and backup memory are continuously compared by U7. If the values do not agree, latch U11 sends a signal (MEMERRINT), enabling the CTC on the I/O Card, to interrupt the processor.

### Analog Switches

The enable lines for the memory chips are switched through U16. When power is applied, the switches are closed, connecting the enable lines. When power is removed or low, the switches are open, disabling the enable lines to preserve power.

### Power-Down and Battery Circuit

When power is removed, or is low, Q1-Q3 disables U16 and opens the analog switches to disable the memory chips to conserve power as they switch over to battery power. The battery provides backup power for the memory chips to preserve the contents when power is removed.

### RS-232C Circuit

This circuit consists of Universal Asynchronous Receiver/Transmitter (UART) U13, baud rate clock U12, U19 RS-232 line driver and receiver, and ferrite beads FB1-FB5 added for reduction of conducted noise. With this circuit, the Main Unit can communicate with the Remote Panel using an RS-232C serial interface.

### Ferrite Beads

FB1-FB5 are used to suppress noise from the power source.

### NOTE: *Clearing of PPI Memory:*

Type C PPI Cards contain a new feature for the simplification of clearing PPI memory. This can be accomplished by removing power, then shorting pin 2 to pin 1 or 3 of JU1. The shorted state should be maintained for several seconds to ensure complete and proper clearing.



# PPI Card

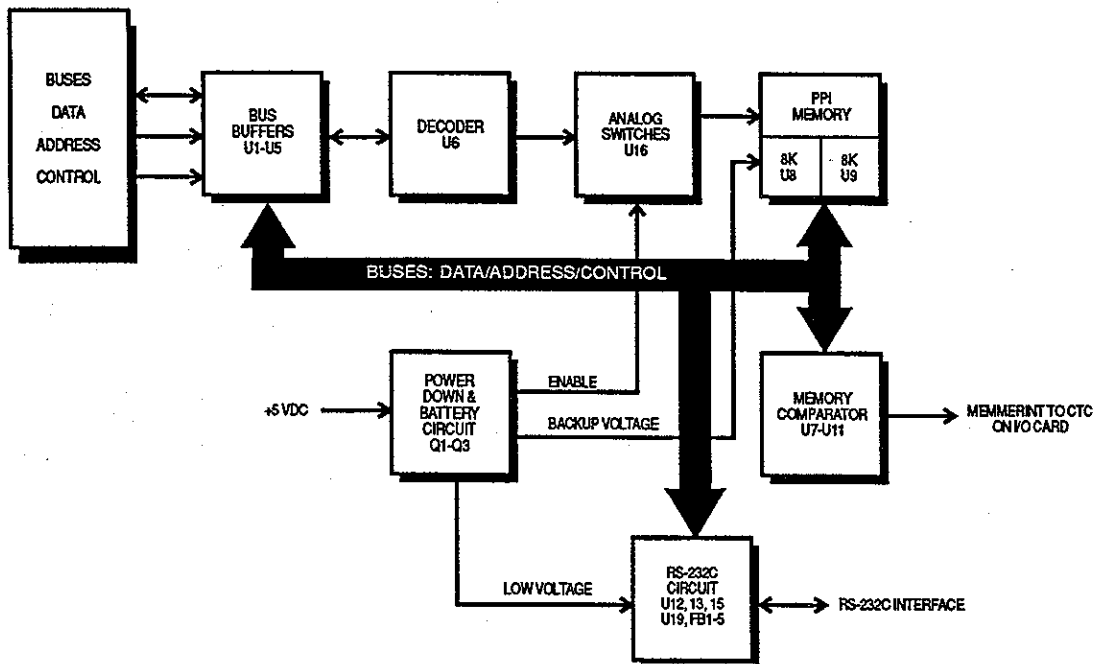


Figure 10.2: Block Diagram: PPI Card Type C

# PPI Card

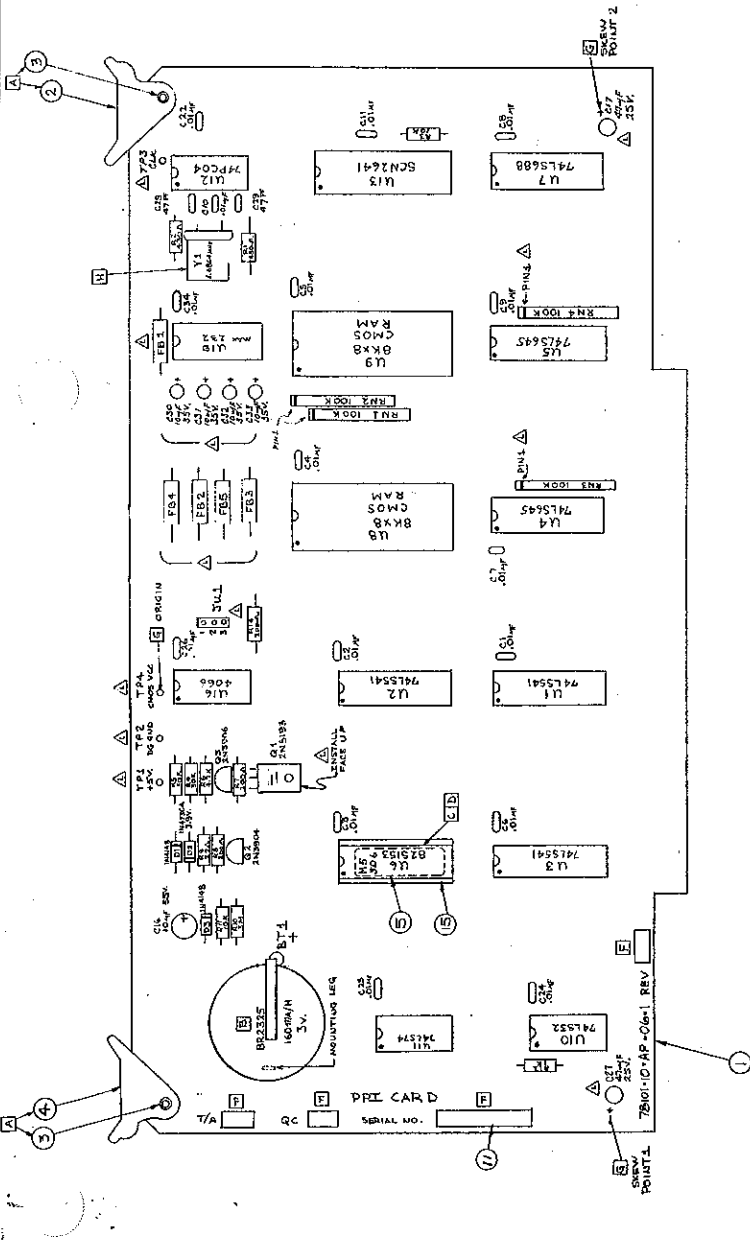
**Notes:**

NOTES:

- 1. INSTALL ITEMS 2314 AFTER WAVE SOLDER
- 2. BATTERY (BAT) NOT TO BE INSTALLED UNTIL ALL ELECTRICAL CHECK-OUT PROCEDURES FOR THIS CARD COMPLETED. BATTERY IS TO BE STORED IN A SAFE MANNER TO PREVENT BATTERY SHORTS OR ACCIDENTAL PREMATURE DISCHARGE. DO NOT WAVE SOLDER BATTERY. HAND SOLDER 5 SECONDS MAX.
- 3. U6 MUST BE PROGRAMMED PRIOR TO INSTALLATION. REFER TO DOCUMENT NO. 89064-00-01 FOR PROGRAMMING INSTRUCTIONS. UNPROGRAMMED U6 IS ON-BIS-300. PAL LABELS ARE TO BE ATTACHED AFTER PROGRAMMING.
- 4. REFER TO DOCUMENT NO. 7810-10-00-04 FOR MASTER DEVICE PROGRAMMING INSTRUCTIONS FOR U6. THIS DOCUMENT FOR ENGINEERING USE ONLY.


PROCESSES	COMPATIBILITY	MARKET
U6	U6	U6
U7	U7	U7
U8	U8	U8
U9	U9	U9
U10	U10	U10
U11	U11	U11
U12	U12	U12
U13	U13	U13
U14	U14	U14
U15	U15	U15
U16	U16	U16
U17	U17	U17
U18	U18	U18
U19	U19	U19
U20	U20	U20
U21	U21	U21
U22	U22	U22
U23	U23	U23
U24	U24	U24
U25	U25	U25
U26	U26	U26
U27	U27	U27
U28	U28	U28
U29	U29	U29
U30	U30	U30
U31	U31	U31
U32	U32	U32
U33	U33	U33
U34	U34	U34
U35	U35	U35
U36	U36	U36
U37	U37	U37
U38	U38	U38
U39	U39	U39
U40	U40	U40
U41	U41	U41
U42	U42	U42
U43	U43	U43
U44	U44	U44
U45	U45	U45
U46	U46	U46
U47	U47	U47
U48	U48	U48
U49	U49	U49
U50	U50	U50
U51	U51	U51
U52	U52	U52
U53	U53	U53
U54	U54	U54
U55	U55	U55
U56	U56	U56
U57	U57	U57
U58	U58	U58
U59	U59	U59
U60	U60	U60
U61	U61	U61
U62	U62	U62
U63	U63	U63
U64	U64	U64
U65	U65	U65
U66	U66	U66
U67	U67	U67
U68	U68	U68
U69	U69	U69
U70	U70	U70
U71	U71	U71
U72	U72	U72
U73	U73	U73
U74	U74	U74
U75	U75	U75
U76	U76	U76
U77	U77	U77
U78	U78	U78
U79	U79	U79
U80	U80	U80
U81	U81	U81
U82	U82	U82
U83	U83	U83
U84	U84	U84
U85	U85	U85
U86	U86	U86
U87	U87	U87
U88	U88	U88
U89	U89	U89
U90	U90	U90
U91	U91	U91
U92	U92	U92
U93	U93	U93
U94	U94	U94
U95	U95	U95
U96	U96	U96
U97	U97	U97
U98	U98	U98
U99	U99	U99
U100	U100	U100

- 5. ALL COMPONENTS MASKED. DO NOT WAVE SOLDER TO HAVE THEM LOCATIONS MASKED PRIOR TO WAVE SOLDER AND BE INSTALLED AFTER WAVE SOLDER. ALL SOCKETED COMPONENTS TO BE INSTALLED AFTER WAVE SOLDER.
- 6. PROMOTION TO COMPLETE SERIAL NO. 4 REV BLOCKS USING PRE-PRINTED LABEL (FROM 11) OR PERMANENT INK T/A 48 SMD. LIGHT GUARD INSERTION ORIENTATION IS T/A. LIGHT GUARDED SMD COMPONENTS TO BE T/A SIDE-SHEW POINT 1, JCT 4 SIDE-SHEW POINT 2.
- 7. USE ADHESIVE TO FASTEN VI TO THE BOARD. USE RTV OR EQUIVALENT REMOVABLE SILICON BASED ADHESIVE.



ITEM	DESCRIPTION	QUANTITY	REVISION RECORD
N	REVISION RECORD		
M	REVISION RECORD		
L	REVISION RECORD		
K	REVISION RECORD		
J	REVISION RECORD		
I	REVISION RECORD		
H	REVISION RECORD		
G	REVISION RECORD		
F	REVISION RECORD		
E	REVISION RECORD		
D	REVISION RECORD		
C	REVISION RECORD		
B	REVISION RECORD		
A	REVISION RECORD		
001	REVISION RECORD		

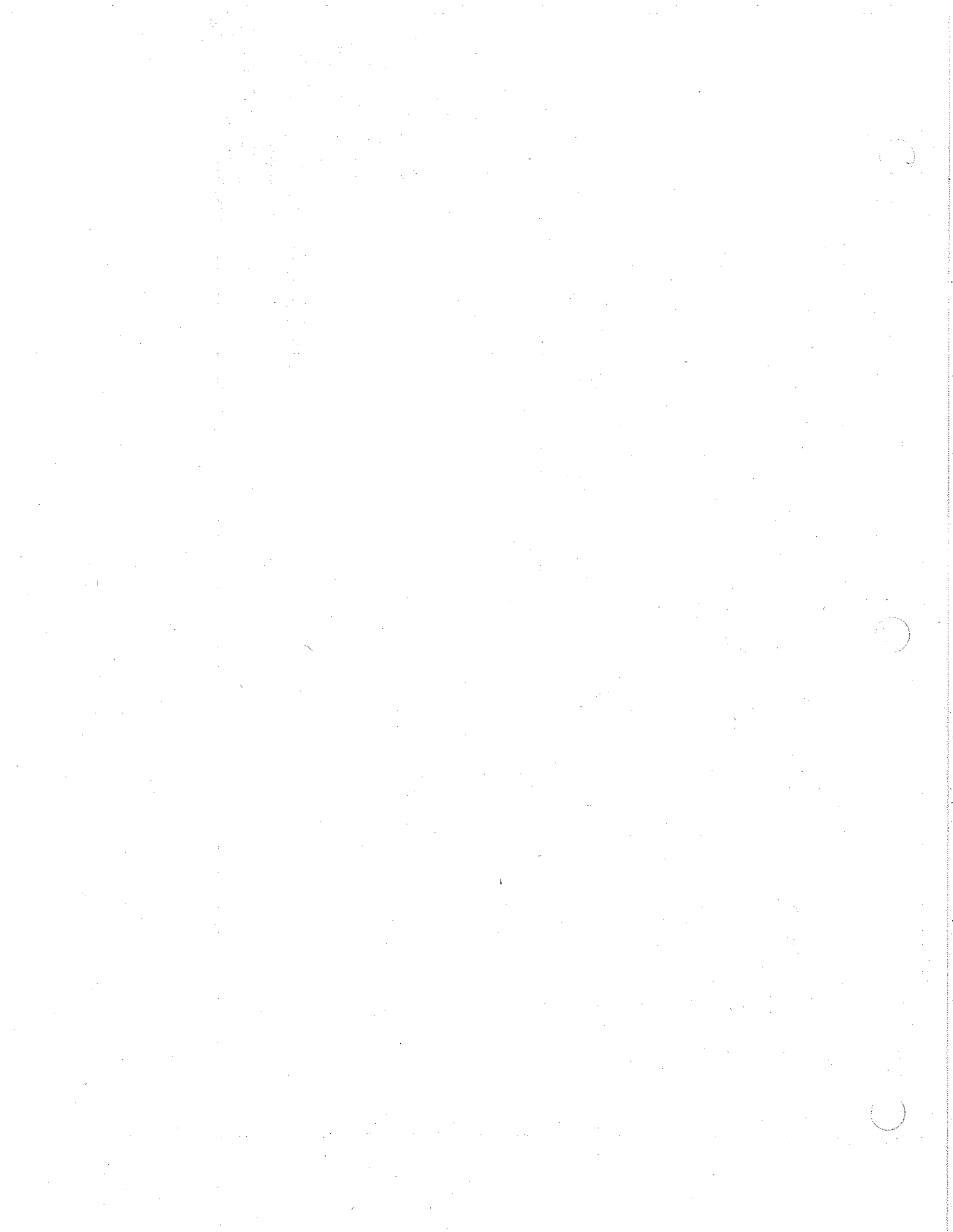All hardware materials used in the assembly of this card are to be controlled under the provisions of MIL-STD-883C, Method 2000, Class B, Level 1, and the provisions of MIL-STD-883C, Method 2000, Class B, Level 1, and the provisions of MIL-STD-883C, Method 2000, Class B, Level 1.



# BILL OF MATERIAL

ITEM	DESCRIPTION	PART NUMBER	QTY.	PROCESS CODE	REF. DESIG ITEM	DESCRIPTION	PART NUMBER	QTY.	PROCESS CODE	REF. DESIG
1	PC CARD, (REV E AND UP)	78107-10-PC-06	1	LG	66	IC, RS232, MAX232	051-0232-000	1	DIP	U19
2	CARD EJECTOR (WHITE)	510-7000-000	1	AW	67	IC, HEX INVERTER, 74PC04	080-74PC-004	1	DIP	U12
3	CARD EJECTOR PIN	510-7001-000	2	AW	68	IC, UART, SCN2641	051-2641-100	1	LG	U13
4	CARD EJECTOR, (YELLOW)	510-7002-000	1	AW	69	IC, COMPARTOR, 74LS688	074-74LS-688	1	DIP	U7
5	LABEL, BLANK PAL	550-0017-619	1	AW	70	IC, PROG. LOGIC ARRAY	071-8215-300	1	AW	NOTE C U6
6	BATTERY, BR2325-1HB	497-2325-100	1	AW	BT1					
10	PIN HEADER, 3 POSITION	518-8001-003	1	BW	JU1					
11	LABEL, BLANK SERIAL NUMBER	550-0093-652	1	BW						
12	TUBULAR TERMINAL	642-1238-000	4	BW	TP1-TP4					
15	IC SOCKET, 20 PIN	545-8201-020	1	LG	U6					
20	CAP. 10UF, 35V	204-0106-135	5	LG	.1-LS C16, C30-C33					
21	CAP. 47UF, 25V	204-0476-025	2	LG	.1-LS C17, C27					
22	CAP. .01UF, 50V	201-0103-050	16	LG	.1-LS C1-01, C21-C26, C28					
23	CAP. CERAMIC, 47PF	201-0470-050	2	LG	.1-LS C28, C29					
25	CRYSTAL, 3.6864 MHZ	180-0037-000	1	LG	Y1					
26	FERRITE BEAD	280-1625-000	5	LG	6 BEND FB1-FB5					
30	TRANSISTOR, 2N5193	125-5193-000	1	BW	O1					
31	TRANSISTOR, 2N3904	120-3904-000	1	BW	O2					
32	TRANSISTOR, 2N3906	121-3906-000	1	BW	O3					
35	DIODE, 1N4148	105-4148-000	2	AXI	5 BEND DI, D3					
36	DIODE, 1N4730A, 3.9V	102-4730-050	1	AXI	5 BEND D2					
40	RESISTOR, 10 KOHM, 1/4 W 5%	303-0103-000	2	AXI	5 BEND R3, R11					
41	RESISTOR, 430 OHM, 1/4 W 5%	303-0431-000	2	AXI	5 BEND R1, R2					
42	RESISTOR, 30 KOHM, 1/4 W 5%	303-0303-000	2	AXI	5 BEND R4, R5					
43	RESISTOR, 4.7 KOHM, 1/4 W 5%	303-0472-000	1	AXI	5 BEND R6					
44	RESISTOR, 200 OHM, 1/4 W 5%	303-0201-000	3	AXI	5 BEND R7, R8, R14					
45	RESISTOR, 27 OHM, 1/4 W 5%	303-0270-000	1	AXI	5 BEND R9					
46	RESISTOR, 3 MOHM, 1/4 W 5%	303-0305-000	1	AXI	5 BEND R10					
47	RESISTOR, 1 KOHM, 1/4 W 5%	303-0102-000	1	AXI	5 BEND R12					
50	RESISTOR ARRAY, 100 KOHM SEP	340-4310-104	4	LG	RN1-RN4					
60	IC, D FLIP-FLOP, 74LS74	074-74LS-074	1	DIP	U11					
61	IC, QUAD 2-INPUT OR 74LS32	074-74LS-032	1	DIP	U10					
62	IC, BUFFER, 74LS541	074-74LS-541	3	DIP	U1, U2, U3					
63	IC, QUAD ANALOG SWITCH, 4066	004-4066-000	1	DIP	U16					
64	IC, CMOS RAM, TC5564-PL15	071-5564-015	2	LG	U8, U9					
65	IC, BUS TRANS, 74LS645	074-74LS-645	2	DIP	U4, U5					

**MEDRAD INC.**  
271 KAPPA DRIVE  
PITTSBURGH, PA 15238

MKV INJECTOR - PPI  
CARD ASSEMBLY

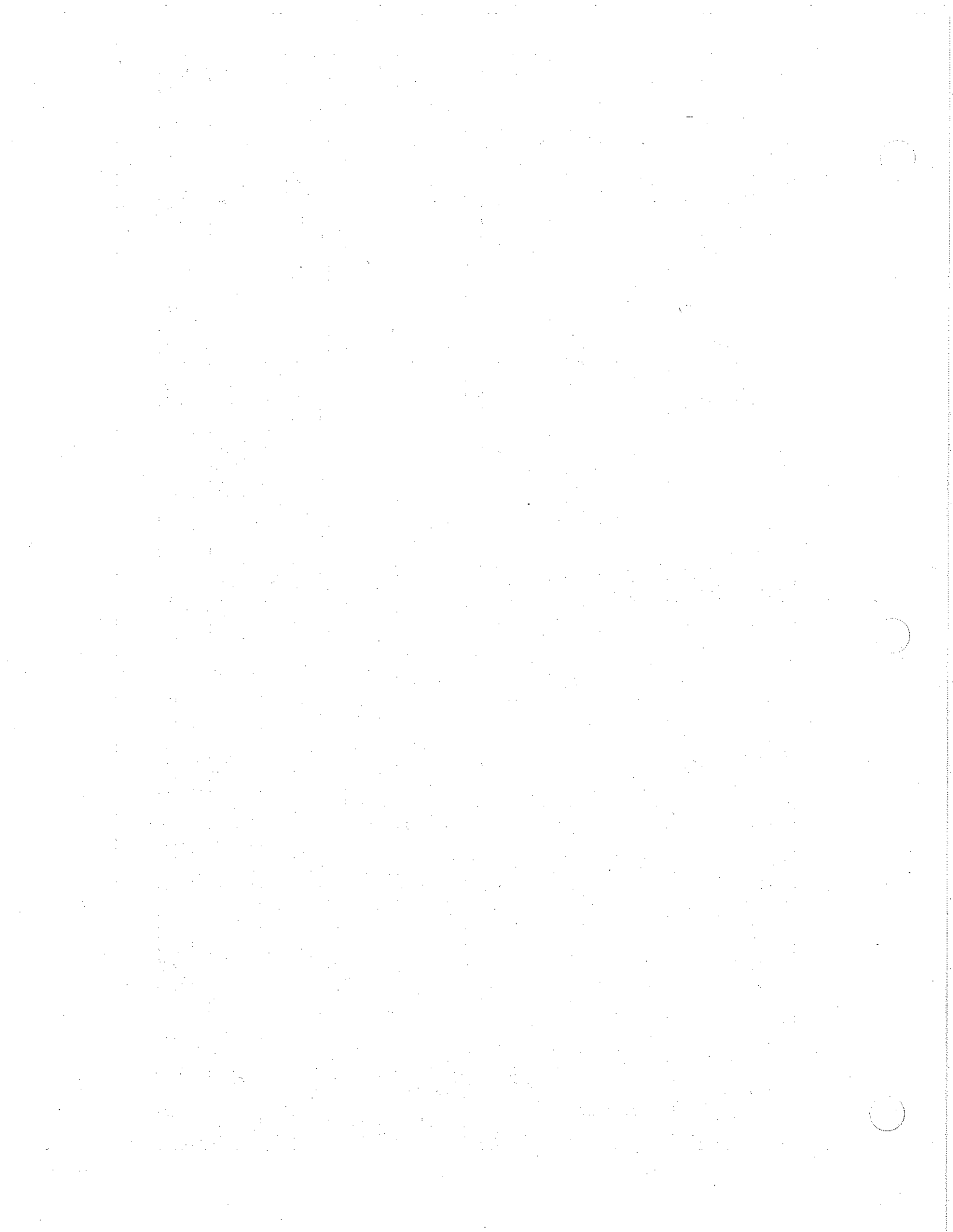
DWG. NUMBER  
**78101-10-AP-006-1**

DATE COPIED  
SHEET 2 OF 2


All technical information equipment design and application data contained within this document is the CONFIDENTIAL Property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.

SEE DATE SHT. 1  
FOR REV. LEVEL

SIGNATURES DATE  
DRAWN MAIESE 3/11/91  
CHECKED KISTE 3/12/91  
ENGR.  
PROJ. MGR.




REVISONS		DATE	APPRD.
REV	ECO NO.	DATE	APPRD.
DESCRIPTION			

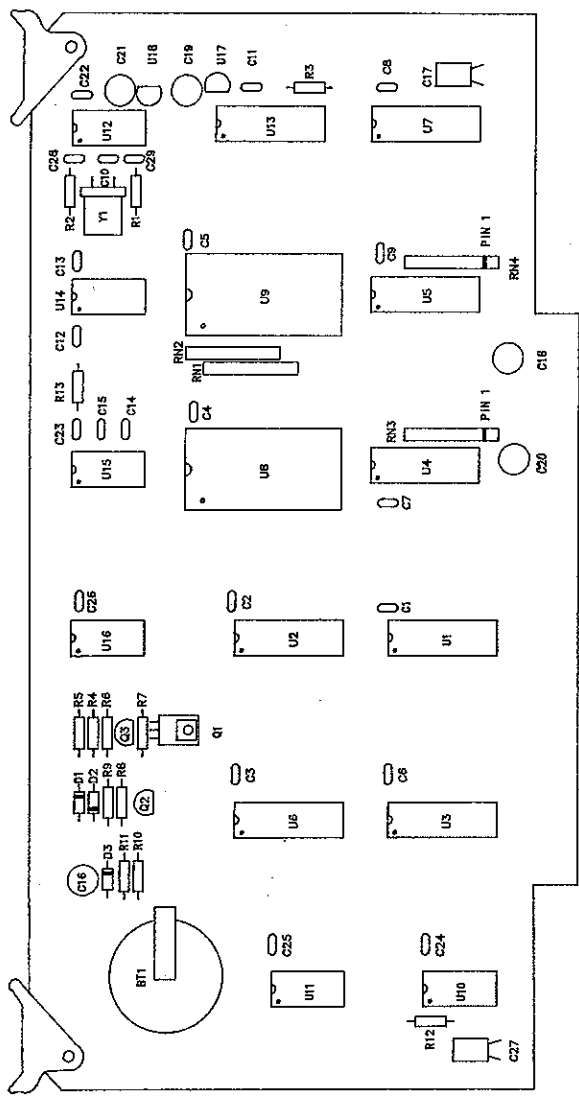

PARTS/MATERIAL LIST			
REF DESIG	DESCRIPTION	PART NUMBER	QTY.
R13	10K 5% 1/4 W	303-0103-000	1
U14	MC1488P	051-1488-000	1
U15	MC1489AN	051-1489-000	1
C18, C19			
C20, C21	10 ufd 35 V	206-0106-350	4
U17	78L12A or PC	052-7812-200	1
U18	79L12A or ACP or CP	052-7912-200	1
C12, C13			
C23	0.01 ufd 100 V	201-0103-100	3
C14, U15	200 pf	201-0201-000	2

Unless otherwise noted, this assembly applicable to boards rev D and below.

		<b>REVISION HISTORY</b> PPI CARD, 20K SERIES	
TITLE REVISION HISTORY PPI CARD, 20K SERIES		SIZE 78101-10-RH-01	
SCALE NTS DATE COPIED		SHEET 1 OF 2	
STANDARD SYMBOLS DIMETER (XX) REF. DIM. DIM. ORIGIN ANGULAR TAPER SLOPE C-BORE C-SINK DEPTH/DEEP DIM. NOT TO SCALE	UNLESS OTHERWISE SPECIFIED DIMENSIONS ARE IN INCHES TOLERANCES ON: 2 PL DECIMALS ± .01 3 PL DECIMALS ± .005 FRACTIONS 0"-8" ± 1/32 FRACTIONS 8"-UP ± 1/16 ANGULAR ± 0.5° MAX. SURFACE TEXTURE 125	SIGNATURES DRAWN CHECKED APPROVED DATE	DO NOT SCALE DRAWING All technical information, equipment design and application data contained within this document is the CONFIDENTIAL property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.

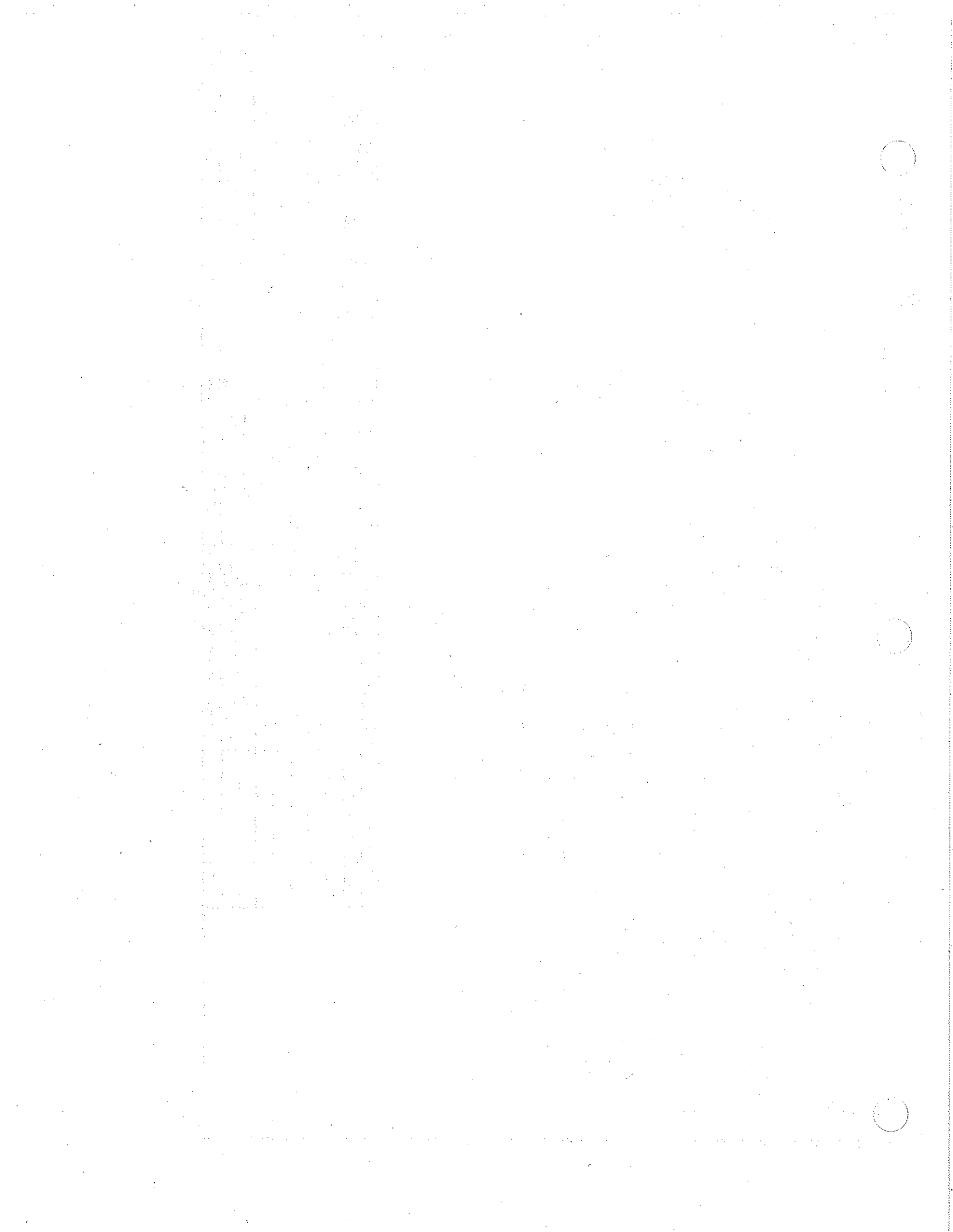



REV. ECO. NO.		DESCRIPTION		DATE	APP'D.

C.R. 78101-10-29. Update artwork negative and assembly drawings to reflect modifications made during the ECG/Scope pilot run. DDN-112. 78101-10-AP-06-1 no rev. to A. 78101-10-PA-06 rev C to D. 78101-10-PD-06 rev E to F.

C.R. 78101-10-32. Added R13, a pull up resistor that allows the PPI card to function properly with a Remote Panel. All cards rev C and up. DDN-115. 78101-10-SE-06 rev D to E.

C.R. 78101-10-42. Drawing corrections only: C14 and C15 were changed from 220 pf to 200 pf and C17 and C27 were changed from 35 V to 20 V as shown in BOM. DDN-147. 78101-10-AP-06-1 rev B to C.

C.R. 78101-10-55. There were 2 part numbers for MC1488. Change 015-1488-000 to 051-1488-000. DDN-234. 78101-10-AP-06-1 rev C to D.

C.R. 78101-10-73. Production related changes. Revise assembly net to install card ejectors and retaining pins after wave soldering. Add PAL labels (item 7) to BOM. Revise note C to reference document 78101-10-PX-02 for programming procedures. DDN-384. 78101-10-AP-06-1 rev D to E.

C.R. 78101-10-90. Revise note C document reference from 78101-10-PX-02 to 78101-00-PX-02. DDN-459. 78101-10-AP-06-1 rev E to F.

C.R. 78101-10-102. Change U1 pin 19 to pin 15. DDN-524. 78101-10-SE-06 rev E to F.

C.R. 10000-NS-358. Revise note C document reference from 78101-00-PX-02 to 89064-00-PX-01. Add note D. C1 through C13, C22, C23 through C26 were changed from 201-0103-050 (0.1" LS) to 201-0103-100 (0.2" LS). DDN-532. 78101-10-AP-06-1 rev F to G.

C.R. 78101-10-105. Change C10 from 201-0103-100 to 201-0103-050 because hole spacing is incorrect. DDN-556. 78101-10-AP-06-1 rev G to H.

C.R. 10000-NS-420. Change C18 through C21 from 206-0106-200 to 206-0106-350. This increased the voltage rating from 20 V to 35 V. DDN-605. 78101-10-AP-06-1 rev H to I.

C.R. 78101-10-115. Added blank label (item 8) and note E. Change voltage ratings of C16, C17, and C27 from 20 V to 35 V, 25 V, and 25 V respectively. Change C28 and C29 from 50 pf to 47 pf. DDN-676. 78101-10-AP-06-1 rev I to J. 78101-10-SE-06 rev G to H.

C.R. 78101-10-123. C17 and C27 were mounted upright. Change orientation so capacitors lay down flat. 78101-10-AP-06-1 rev J to K. DDN-723.

EC-24. U2 had an extra pin 17 that was changed to pin 12. 78101-10-SE-06 rev H to I.

EC-114. PPI card has been relayed out with significant changes. The following components have been deleted on rev E board: R13, U14, U15, C18, C19, C20, C21, U17, U18, C12, C13, C23, C14, C15. The following components have been added: J11, U19, TP1 through TP4, FB1 through FB5, C30 through C33, and R14. RN3 and RN4 were previously oriented with pin 1 toward the card edge. RN3 and RN4 are rotated 180 degrees on boards rev E and up. Q1 was installed face down. On rev E boards and up Q1 will be installed face up. C17 and C27 were mounted flat. Change orientation so capacitors are upright. Rotate C10, C28, and C29 90 degrees. Change C16 from a 100 ufd to 10 ufd. Added notes F and G. 78101-10-AP-06-1 rev K to L. 78101-10-SE-06 rev I to J. Add note 3 to show schematic of boards rev D and below.

10-10

		<b>REVISION HISTORY</b> <b>PPI CARD, 20K SERIES</b>	
TITLE		SIZE DWG. NUMBER <b>B</b> 78101-10-RH-01	
SCALE NTS		DATE COPIED SHEET 2 OF 2	
GEOMETRIC TOLERANCES SYMBOL CHARACTERISTIC	STANDARD SYMBOLS Ø DIAMETER (XX) REF. DIM. DIM. ORIGIN CONICAL TAPER SLOPE C-SHORE DEPTH/DEEP DIM. NOT TO SCALE	UNLESS OTHERWISE SPECIFIED DIMENSIONS IN INCHES TOLERANCES ON: 2 PL DECIMALS ± .01 3 PL DECIMALS ± .005 FRACTIONS 0"-8" ± 1/32 FRACTIONS 8"-UP ± 1/16 ANGULAR ± 0.5° MAX. SPREAD TEXTURE 125/	SIGNATURES DRAWN: [Signature] CHECKED: [Signature] ENGR: [Signature] PROJ. [Signature] DATE: [Date]
<b>DO NOT SCALE DRAWING</b>			
All technical information, equipment design and application data contained within this document is the CONFIDENTIAL property of MEDRAD, INC. No use or disclosure thereof may be made without written permission from MEDRAD, INC.			



## SECTION 11 - CONTROL PANEL SUMMARY OF OPERATION

### Control Panel Summary of Operation

The Control panel is comprised of the keyboard overlay attached to the Display Card. Refer to page 11 - 2 for a block diagram of the Control Panel cards.

### Keyboard Overlay

The keyboard overlay contains rubber membrane keys, interconnected in a matrix of four columns and five rows. This matrix is connected to the Display Card, then to the I/O Card, where PIO2 senses key closures for the system.

### Display Card

The Display Card contains the numeric displays and LEDs that show through the keyboard overlay. These devices convey messages and operating status to the operator.

### Numeric Display Drivers and Numeric Displays

Three display groups are used to display injection parameters on the Control Panel. Each display group consists of a display driver (U1, U2, or U3), and numeric display units. Display drivers also receive the data and control line information from the I/O Card. Group designations are as follows:

- Display Driver U1, displays Volume Remaining (DS1 - DS3), First Phase Flow Rate (DS6 - DS7), and Scan Delay (DS4, DS5).
- Display Driver U2, displays First Phase Volume (DS8 - DS10), Second Phase Flow Rate (DS11, DS12), and Second Phase Volume (DS13 - DS15).
- Display Driver U3, displays Volume Delivered (DS16 - DS18), and Injection Duration (DS19 - DS22).

### LED Driver and LEDs

The remainder of the circuit on the display card controls and drives the LEDs used as indicators and that backlight various words printed on the Control Panel. Display driver U4 and buffer U5 control all LEDs (LP1-LP15), receiving data and control line information from the I/O Card. U4 and U5 control the displays: First Phase, Second Phase, Insufficient Volume, Phase 1, Phase 2, Single Arm, Multi Arm, Armed, and Program Levels 1-6.

# Control Panel

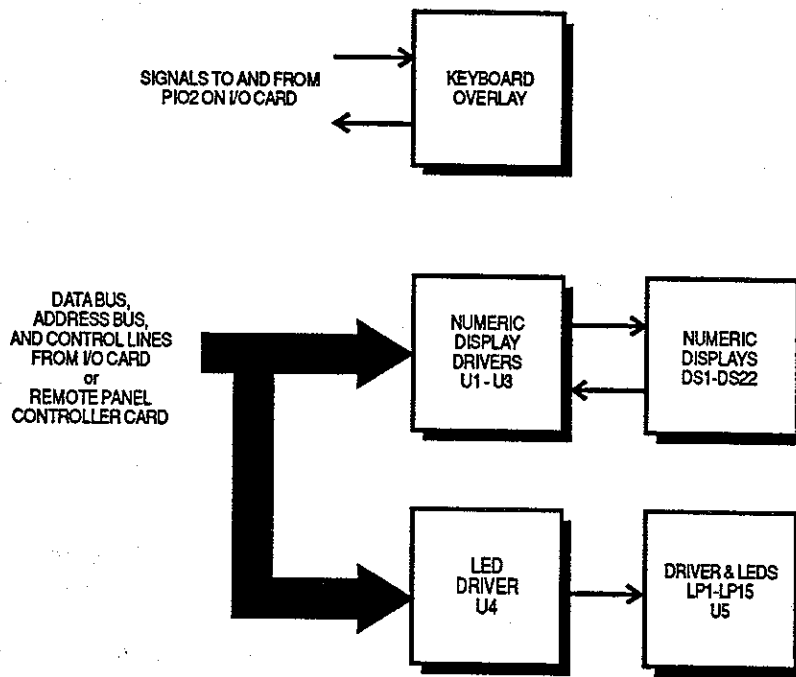


Figure 11.1: Block Diagram: Control Panel